

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			

Site Information

Name: _____ Installation Address: _____ _____ _____ Contact Name: _____ Title: _____ Phone #: _____ Fax: _____ E-mail: _____
--

Training Information (If available)

Operator: _____ Title: _____ Phone #: _____ Fax: _____ E-mail: _____ Experience with automated systems: _____ _____ Operator Training dates: _____ Operator proficiency testing dates on-site: _____
--

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			

Installation Logistics (if available)

<i>Pre-site visit date:</i> _____	<i>Final site inspection date:</i> _____
<i>Installation date:</i> _____	
<i>LIS Supplier:</i> _____	<i>Phone #:</i> _____
<i>E-mail:</i> _____	
<i>IT contact name:</i> _____	<i>Phone #:</i> _____
<i>E-mail:</i> _____	

List of materials needed for the site evaluation:

1. Tape measure (5m or greater)
2. Digital camera (if available)

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			

Installation Requirements

The following information is needed to assess the readiness of a facility before installing the Erytra Eflexis® System. This checklist is to be used by a qualified service representative who should carefully review each item and provide any comments as needed during the inspection.

Receiving/Unloading	Yes/No (☐ one)	Comments
Number of instruments to be delivered		
Delivery hours		
Receiving area access by truck (note any size restrictions)	☐ Yes / ☐ No	
Receiving area for unloading and unpacking the system (1,5 m x 5 m)	☐ Yes / ☐ No	
Is there a forklift or pallet jack available at the site for use with proper notice?	☐ Yes / ☐ No	
Room to store the crate (120 cm x 80 cm x 137 cm) temporarily if necessary	☐ Yes / ☐ No	
Can the site dispose of the crate or will Grifols be required to manage the return, disposal, recycling of the materials?	☐ Yes / ☐ No	
Is there a route to the laboratory with minimum clearance (including doorways) of 91 cm x 71 cm (H x W)? No weight restrictions <200 Kg	☐ Yes / ☐ No	
Does the instrument require to be moved to an elevated level?	☐ Yes / ☐ No	
Can we use an elevator?	☐ Yes / ☐ No	
Is there a suitable elevator (capacity > 300 Kg)	☐ Yes / ☐ No	
Lifting equipment requested? Provide details and indicate if carrier visit is requested	☐ Yes / ☐ No	
Space requirements	Yes/No (☐ one)	Comments
Has an appropriate space out of direct sunlight been reserved for the analyzer(s)? FSE should have clear access to both sides of the instrument.	☐ Yes / ☐ No	

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			

Space requirements	Yes/No (☐ one)	Comments
If UPS required, has space near the instrument been allocated? <i>Note: UPS is standard part of instrument configuration in US</i>	☐ Yes / ☐ No	
Is there a separate designated bench for racking samples?	☐ Yes / ☐ No	
Is a Grifols bench requested (132 cm x 80 cm x 80 cm)? <i>See Erytra properties for external card waste</i>	☐ Yes / ☐ No	
Is the instrument being installed in a region with seismic requirement considerations? <i>For example OSHPD in California.</i>	☐ Yes / ☐ No	
Has a 15°C to 25°C storage area been identified for Wash Solution Fluids, DG Gel cards, consumables, racks and spare containers?	☐ Yes / ☐ No	
Has an adequate 2°C to 8°C storage area been planned for assay reagents? Are additional fridges required?	☐ Yes / ☐ No	
Can the site provide a solution $\geq 0.5\%$ sodium hypochlorite and deionized (DI) water?	☐ Yes / ☐ No	
Power requirements	Yes/No (☐ one)	Comments
Power plug type		
Is there a suitable socket available? (please state distance in m)	☐ Yes / ☐ No	
Has lab confirmed supply voltage either 120 or 240 VAC +/- 10%?	☐ Yes / ☐ No	

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			

Will an UPS be required?	☐ Yes / ☐ No	
Environmental requirements	Yes/No (☐ one)	Comments
Temperature: 15°C to 30°C? Humidity: 15% to 85% including heat dissipation of approximately 0.26 kWh /936 kJ /223.79 kcal per analyzer	☐ Yes / ☐ No	
Elevation ≤ 3,000 m?	☐ Yes / ☐ No	
Direct exposure beneath an air conditioning vent?	☐ Yes / ☐ No	
Sample barcodes	Yes/No (☐ one)	Comments
Samples barcode configuration(s) Code39, Code 128 A, B & C (ISBT 128), Interleaved 2 of 5, Codabar, EAN8/EAN13/JAN8/JAN13. Any apparent quality issues? If possible, take some sample tubes labels to test prior to installation.		
Are all used barcodes types including a checksum?	☐ Yes / ☐ No	
Are there any potential issues with barcode placement? (Multiple labels, conflicting positions for other analyzers, label heights, suspect quality)?	☐ Yes / ☐ No	

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			


Erytra Eflexis[®] System footprint

Customer's Erytra Eflexis[®] placement and laboratory layout.

Sketch the proposed Erytra Eflexis[®] instrument placement and layout on the grid below or submit a digital representation of the proposed layout (for example, CAD, Excel).

Include the following information in the sketch:

- Room location
- Location of the Erytra Eflexis[®] instrument

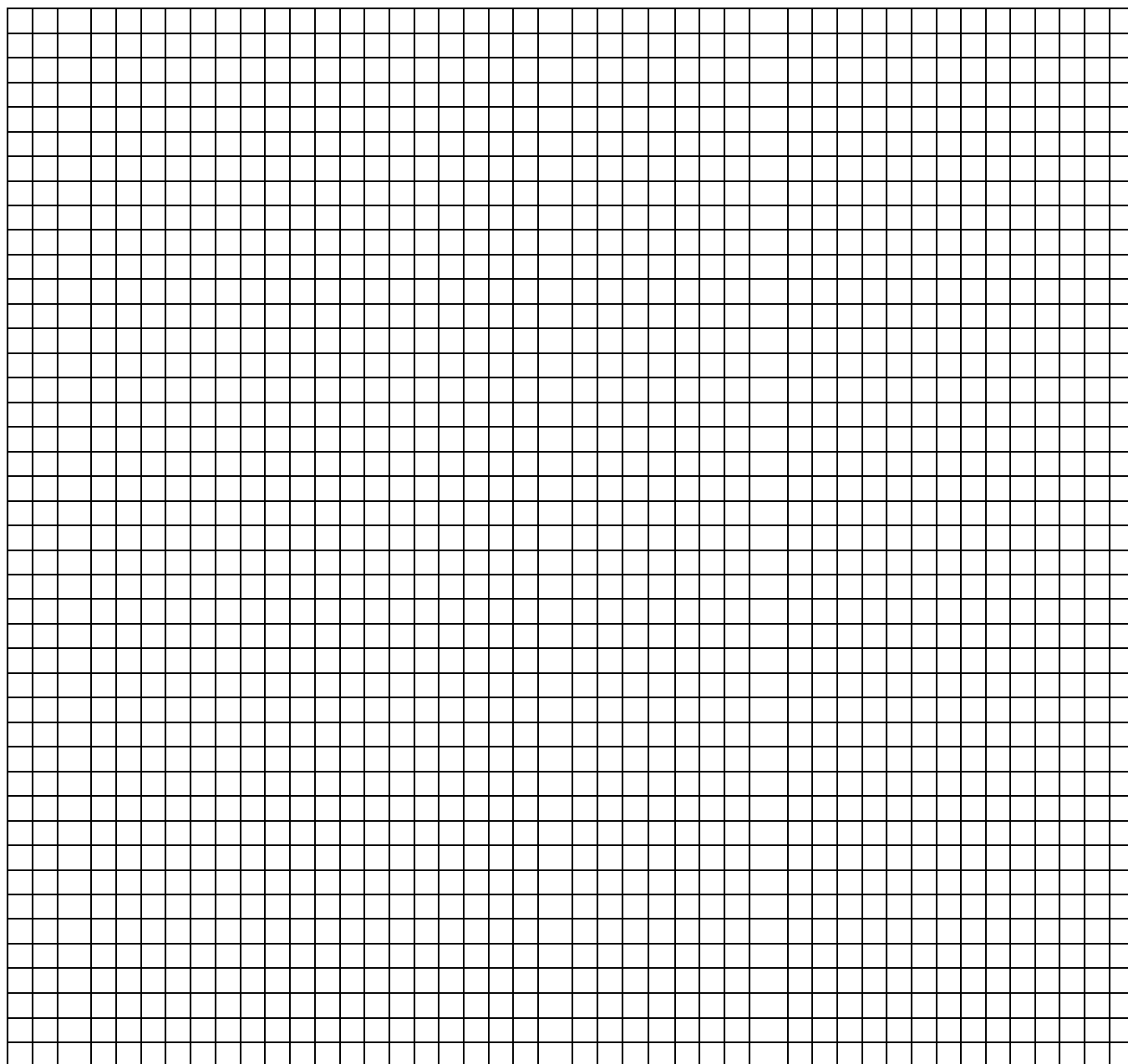

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			

Sample centrifugation	Yes/No (☐ one)	Comments
Laboratory centrifuge can be configured according to Grifols specifications?	☐ Yes / ☐ No	
Network configuration	Yes/No (☐ one)	Comments
Will the Erytra Eflexis® be connected to a network?	☐ Yes / ☐ No	
Is there an Ethernet port (or hub) present near the instrument? If yes, what is the required length of the network cable from the Ethernet port to the instrument?	☐ Yes / ☐ No	
Remote access available?	☐ Yes / ☐ No	
Back Up/Archive space available?	☐ Yes / ☐ No	
LIS configuration	Yes/No (☐ one)	Comments
Select the LIS option (select one)	☐ No LIS ☐ Serial ☐ File ☐ TCP/IP	
Diana batch system active?	☐ Yes / ☐ No	
LIS Vendor name		
Lis Vendor Contact		
LIS Software name		
Has ASTM protocol been provided to vendor?	☐ Yes / ☐ No	
Is LIS connection shared?	☐ Yes / ☐ No	
Indicate if direct connection or shared server PC are used		

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			

Printer	Select option	Comments
Printer connection <i>Note: Printer (USB) is standard part of instrument configuration in US</i>	☐ None ☐ USB ☐ Network	
Erytra properties configuration	Select option	Comments
Service Rack mode	☐ None ☐ 1 ☐ 2 ☐ 4	
Result configuration	☐ SRW ☐ SRDL ☐ SRDP ☐ SROD ☐ RI ☐ PI ☐ CSI ☐ SRI ☐ ALL	
All sample tubes type in use meets requirements (Type: plastic or glass siliconized, Diameter: 9-16 mm, Height: ≤ 100 mm)	☐ Yes / ☐ No	
Sample tubes type (specify characteristics, e.g. diameters and bottom thickness)	☐ NORMAL ☐ PLUNGER ☐ PEDIATRIC	
Samples rack labels to use (specify special needs per each type)	☐ GREY ☐ RED ☐ YELLOW ☐ BLUE (pediatric)	

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			

Erytra properties configuration	Select option	Comments
LIS Mode	☐ Q ALL/ SM True ☐ Q ALL/ SM False ☐ Query by Sample	
QC Warning new lot	☐ Yes / ☐ No	
Results validation button disabled	☐ Yes / ☐ No	
Results Export button disabled	☐ Yes / ☐ No	
Password re-authentication	☐ Yes / ☐ No	
Password expiration days		
Confirmation message on results modification	☐ Yes / ☐ No	
Archive (if yes specify destination and time interval)	☐ Yes / ☐ No	
Delta check (if yes specify criteria and interval)	☐ Yes / ☐ No	
Demographic data management	☐ Yes / ☐ No	
Double tube management	☐ Yes / ☐ No	
QC policy activated	☐ Yes / ☐ No	
QC automatic management	☐ No ☐ Only validation ☐ Validation & export	
Samples auto-validation and export	☐ Yes / ☐ No	
ISBT barcodes	☐ Yes / ☐ No ☐ 13 digits ☐ 16 digits ☐ Default manual	

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			

Erytra properties configuration	Select option	Comments
External drain (specify)	☐ Yes / ☐ No	
External card waste (specify)	☐ Yes / ☐ No	
Fluidic system always ready	☐ Yes / ☐ No	
Tests	Yes/No (☐ one)	Comments
List the tests to enable on Erytra Eflexis®		
Do any of the tests to enable require a DG Clean adaptor?	☐ Yes / ☐ No	

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			

Tests	Yes/No (☐ one)	Comments
Not Unitary Profile management (if yes, specify the composition of profiles)	☐ Yes / ☐ No	
Special results to configure		
Operator training		Comments
Number of basic operators		
Number of advanced operators		
Starting time (morning session) Starting time (afternoon session)		

GRIFOLS	Document #:	405749	Effective Date:	15-Feb-2019
	Version:	4.0	CONFIDENTIAL	
	Class/Type:	GST Service Document - Form		
Title:	Erytra Eflexis Pre-Installation Checklist Form			

INSTALLATION PLANNING

ACTION	Delivery, IQ & OQ		Dates:	
GRIFOLS representative		Customer representative		
ACTION	LIS connection		Dates:	
GRIFOLS representative		Customer representative		
ACTION	Advanced operator training		Dates:	
GRIFOLS representative		Customer representative		
ACTION	Basic operator training		Dates:	
GRIFOLS representative		Customer representative		

Grifols Representative	
Printed name	
Signature	Date

Customer Acceptance Signature	
Printed name	
Signature	Date