

Luminex

Field Service Advisory

Title: LX200 CON1 Failing With Zero Events in xPONENT 4.2.1509

Page 1 of 5

Effective Date: 09/27/2018 Rev. A

Distribution:

World-wide

Reason for Advisory:

CON1 is failing in xPONENT 4.2.1509 with zero events observed.

Reference/Source of FSA:

NA

Applicable Hardware/Software/Reagent:

Luminex 200

xPONENT 4.2.1509

LX200-CAL-K25

LX200-CON-K25

Specialty Equipment Required:

N/A

Part Number(s):

N/A

Details:

It has been discovered that when running Calibration and Verification in xPONENT 4.2.1509 that CAL1, MagCAL, CAL2, MagCON, CON2, and Fluidics will all pass; however, CON1 will fail with zero events observed. This is caused by an unknown issue in the software that corrupts the NVRAM on the Digital Board, and affects new systems with xPONENT 4.2.1509 installed. The root cause investigation is ongoing. If this issue occurs, perform the procedure below:

NOTE: This procedure must be conducted in xPONENT using the Cmds & Routines page, and each command must be run manually.

Luminex

Field Service Advisory


Title: LX200 CON1 Failing With Zero Events in xPONENT 4.2.1509

Page 2 of 5

Effective Date: 09/27/2018 Rev. A

1. Click Eject and place the following reagents in the specified well using the figure below:

NOTE: Ensure the correct lot information has been loaded into xPONENT.

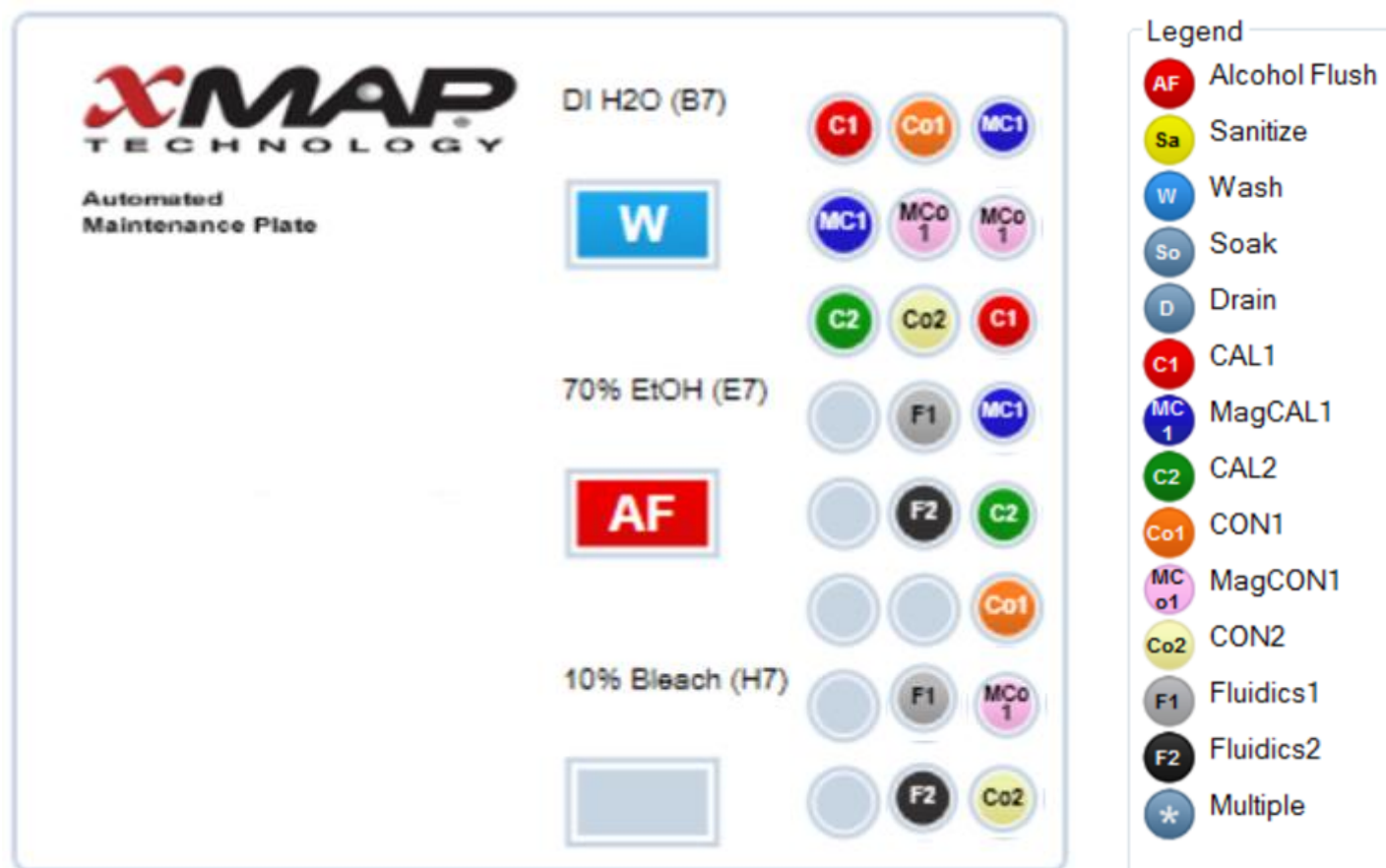

2. Open xPONENT and navigate to the **Maintenance Page > Cmds & Routines tab**.
3. Select well **A12** and click the **MagCAL1** button.
4. Select well **B12** and click the **MagCON1** button.
5. Click **Run** and wait for the commands to be completed.

Maintenance Commands and Routines

Instructions Run single or multiple commands, or run a routine. Create a routine for the selected commands by clicking "Save As".

Routine Name: <None>

Commands

- Wash (W)
- Soak (So)
- CAL1 (C1)
- CON1 (Co1)
- Alcohol Flush (AF)
- Warmup
- MagCAL1 (MC1)
- MagCON1 (MCo1)
- Sanitize (Sa)
- Prime
- CAL2 (C2)
- CON2 (Co2)
- Drain (D)
- Backflush
- Fluidics

Command Sequence

Command	Lot Number	Location	Reagent	Status	Information
MagCAL1	B63932	1,A12		Pending	
MagCON1	B64889	1,B12		Pending	

Run

- Power off the instrument then shut down the PC and wait ~ 1 minute.
- Power on the PC then power on the instrument (allow the instrument to initialize).
- Open xPONENT and navigate to the **Maintenance Page > Ccmds & Routines** tab.
- Select well **C12** and click the **CAL1** button.
- Select well **D12** and click the **MagCAL1** button.

Luminex

Field Service Advisory

Title: LX200 CON1 Failing With Zero Events in xPONENT 4.2.1509

Page 4 of 5

Effective Date: 09/27/2018 Rev. A

11. Select well **E12** and click the **CAL2** button.
12. Select well **F12** and click the **CON1** button.
13. Select well **G12** and click the **MagCON1** button.
14. Select well **H12** and click the **CON2** button.
15. Select well **G11** and click the **FLUIDICS** button. (FLUIDICS2 will automatically populate in well H11).
16. Click **Run** and wait for the commands to be completed.

Maintenance Commands and Routines

Instructions Run single or multiple commands, or run a routine. Create a routine for the selected commands by clicking "Save As".

Routine Name: <None>

Commands

- Wash (W)
- Soak (So)
- CAL1 (C1)
- CON1 (Co1)
- Alcohol Flush (AF)
- Warmup (Warmup)
- MagCAL1 (MC1)
- MagCON1 (MCo1)
- Sanitize (Sa)
- Prime (Prime)
- CAL2 (C2)
- CON2 (Co2)
- Drain (D)
- Backflush (Backflush)
- Fluidics (F)

Well Plate Diagram

	1	2	3	4	5	6	7	8	9	10	11	12	R	1
A														
B														
C														
D														
E														
F														
G														
H														

Command Sequence

Command	Lot Num...	Location	Reagent	Status	Information
CON1	B65915	1,F12		Pending	
MagCON1	B64889	1,G12		Pending	
CON2	B64890	1,H12		Pending	
Fluidics1		1,G11		Pending	
Fluidics2		1,H11		Pending	

Buttons: Self Test, Report, Import, Save As, Cancel, Run

Document Number: 01815 Rev. C

Printed copies of this document are for reference purposes only.

The only controlled copy of this document is available in the Luminex electronic document control system.

17. In xPONENT, navigate to the **Maintenance** page > **Auto Maint** tab.
18. Ensure **Calibration Verification** is selected, and the correct lots are selected on both dropdown menus.
19. Click **Run** and wait for the commands to be completed.
20. Verify passing results. If not passing, contact Hardware Product Services for further troubleshooting assistance.

Luminex User: [] Log Off Help xPONENT SOFTWARE SOLUTIONS

Home Samples Batches Results Protocols **Maintenance** Admin

Auto Maint

Lot Management

Cmds & Routines

Probe & Heater

System Info

System Status

Schedule

Support Utility

Automated Maintenance - Calibration Verification

Instructions Select the Automated Maintenance Option or Maintenance Command. Then, select the appropriate kit if applicable.

Automated Maintenance Options

☒ Calibration Verification ☐ Performance Verification ☐ Fluidics Prep ☐ System Shutdown

Command	Location	Reagent	Status	Information
Warmup		None	Pending	
Prime		None	Pending	
Alcohol Flush	1,E7	70% Alcohol	Pending	
Wash	1,B7	DI H2O	Pending	
Wash	1,B7	DI H2O	Pending	
CAL1	1,A10	CAL1	Pending	
MagCAL1	1,B10	MagCAL1	Pending	
CAL2	1,C10	CAL2	Pending	
Wash	1,B7	DI H2O	Pending	
CON1	1,A11	CON1	Pending	
MagCON1	1,B11	MagCON1	Pending	
CON2	1,C11	CON2	Pending	
Wash	1,B7	DI H2O	Pending	
Fluidics1	1,D11	Fluidics1	Pending	

Reagents

Calibration Kit: B64462 Performance Verification Kit: B66893

Legend

- Alcohol Flush
- Sanitize
- Wash
- Soak
- Drain
- CAL1
- MagCAL1
- CAL2
- CON1
- MagCON1
- CON2
- Fluidics1
- Fluidics2
- Multiple

Automated Maintenance Plate

DI H2O (B7) C1 Co1

70% EtOH (E7) MC1 MC2

AF F1 F2

10% Bleach (H7)

Report Import Kit **Run**

System Status Command: [] System State: [] Thursday 9/27/2018 10:28 AM Stop Pause Eject

Delta Cal Temp: [] Laser: []

Sheath Pressure: [] Region Events: []

XY Status: [] Total Events: [] Warm Up

If you have any questions or need further assistance contact your immediate supervisor, or hardwareprodsvc@luminexcorp.com

Document Number: 01815 Rev. C

Printed copies of this document are for reference purposes only.

The only controlled copy of this document is available in the Luminex electronic document control system.