

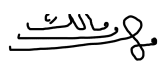


Compo Guard
Yearly MAINTENANCE
CHECKLIST



Work Order No	WO-00161763	PPM Date	6/NOV/2025
Region	Western	Hospital	Andalusia Hospital - Hai Aljamea Jeddah
Serial No	20ACA0013	Department	Blood Bank

Action		Checked / Done	Remarks
General checks			
1	Visual inspection: instrument general status.	☒	
2	Mains Inlet Status	☒	
Touch buttons and LCD functionality Check			
1	Check and verify the functionality of touch buttons	☒	
2	Check if the Display is Working properly	☒	
Decontamination & Cleaning			
1	Turn off the instrument, disconnect the power cord.	☒	
2	If there is spillage, wipe it with disposable absorbent material (for example, paper towels, gauze, or cloth).	☒	
3	Wipe the entire instrument with disposable towels or cotton balls dipped in a 0.5% (v/v) aqueous solution of sodium hypochlorite	☒	
4	Rinse the surface with a soft detergent and water to eliminate odors and any residual solution.	☒	
5	Dry the surface.	☒	
6	Dispose of contaminated materials used during decontamination in a container for biological waste.	☒	

Service Engineer	Malek Aldweik	
-------------------------	---------------	---