

OQ Injectors Cleaning Procedure


Background:

The bleach and buffer tend to buildup crystals around the OQ injector, often leading to an occlusion of the OQ MTU presence sensor. When this occurs, the instrument will crash with a - “The MTU presence sensor in the Output Queue detects a MTU after the loading process” error (081.000.00069). To prevent this, the OQ injector and OLV housing should be cleaned of all crystal buildup.



Procedure:

1. Warning- the OQ is considered a “dirty” area. Be very conscientious of contamination risks.
2. Wear proper PPE.
3. Power down the instrument.
4. Open the Right Side instrument panel.
5. Unplug the OQ OLV Board.

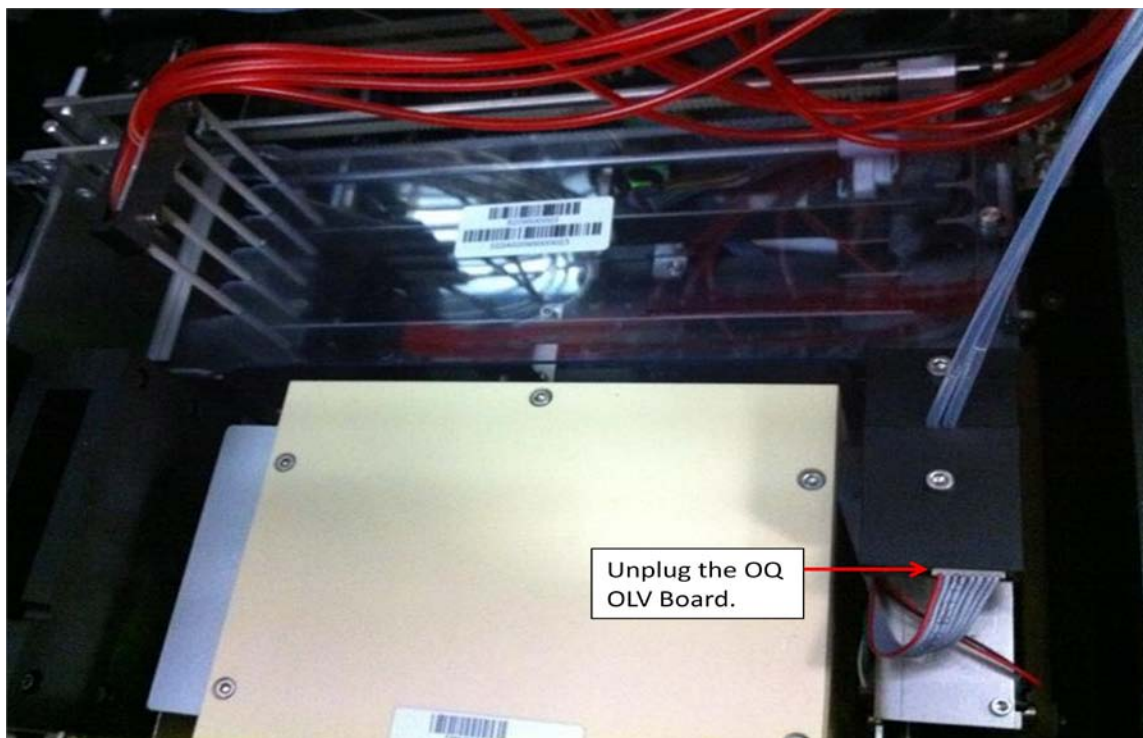


6. Loosen the OQ OLV Housing thumbscrew.




7. Place an absorbent pad over the OQ for protection.
8. Remove OQ OLV Housing/Injector (OQ Injector will still be connected to the pump) and place on the absorbent pad.



9. Clean the Housing and Injector with 50% bleach to prevent any contamination.
10. Change gloves.
11. Rinse with DI Water, and remove any bleach/buffer buildup.
12. Dry the Housing and Injector.
13. Remove absorbent pad.
14. Clean any debris from the OQ OLV mount location with bleach and DI water.

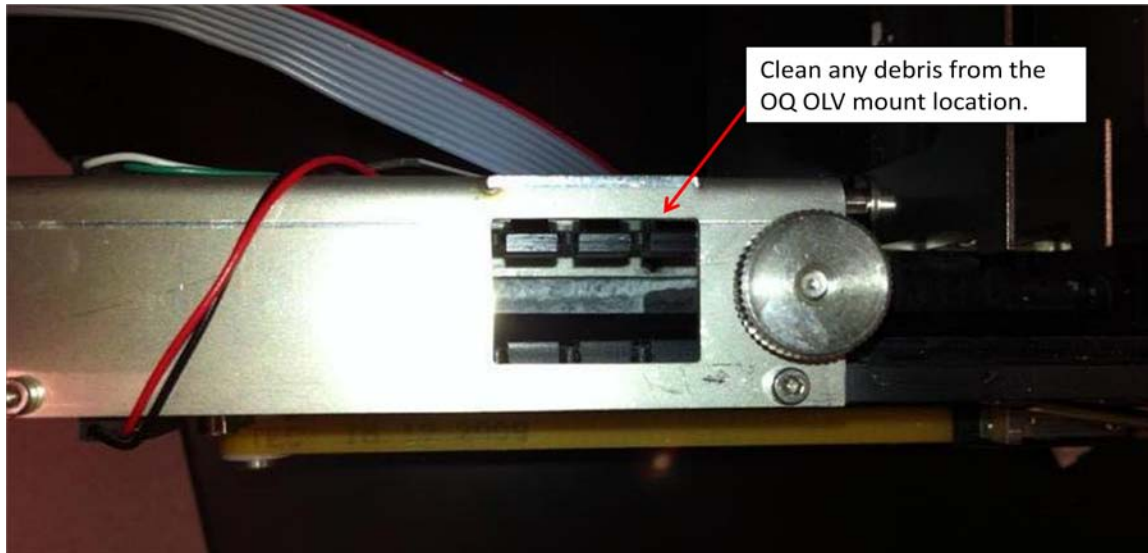

15. Replace OLV Housing/Injector.
16. Tighten the OQ OLV Housing thumbscrew.
17. Plug in the OQ OLV Board.
18. Using a cotton swab and ethanol, wipe clean the OQ MTU presence sensor.

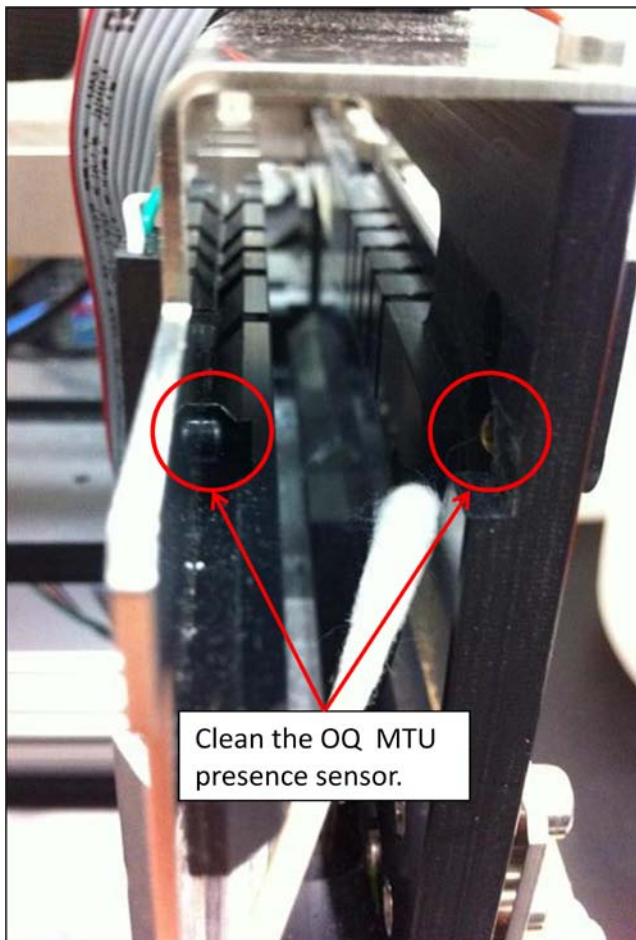

19. Close the Right Side instrument panel.
20. Start the Panther GUI, and perform a full Prime for verification.