

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

DG-57 Service Manual

SOFTWARE INSTALLATION

Procedure description
This procedure describes the software installation steps to follow for both the Matrox 4Sight GP computer and a standard computer. <i>Note: Document REGD-0000777_v6 replaces document REGD-0000777_v5.</i>

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

EQUIPMENT

1.	DG-57 instrument.
2.	Matrox 4Sight GP Computer / Standard computer fulfilling the specified requirements.
3.	Matrox Meteor II / Morphis board
4.	Software CDROM
5.	USB drive

Installation procedures

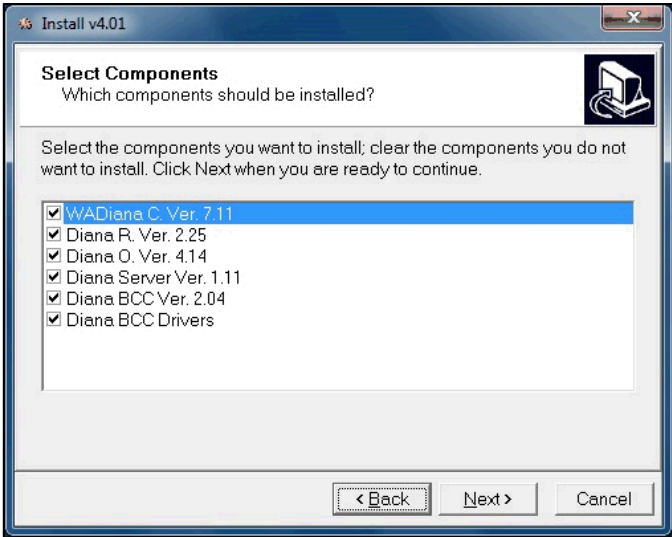
- **PROCEDURE 1: SOFTWARE VERSION $\geq$ 4.2.0**
 - **Case 1: PC Matrox 4Sight GP**
 - **Case 2: Desktop PC (Windows 7)**
 - **Case 3: Desktop PC (Windows XP)**
- **PROCEDURE 2: SOFTWARE VERSION 3.2**

PROCEDURE 1: SOFTWARE VERSION $\geq$ 4.2.0

Case 1: PC Matrox 4Sight GP

Process	Expected Results
1. Copy the content of the Diana System SW CD into the root of a USB drive and connect it to the PC. The system automatically detects the installation program, prompting the AutoPlay window. Select Run install.exe (Picture 1). Picture 1 If the system does not detect the installation program, then run Install.exe from the USB drive contents.	The installation wizard is started.

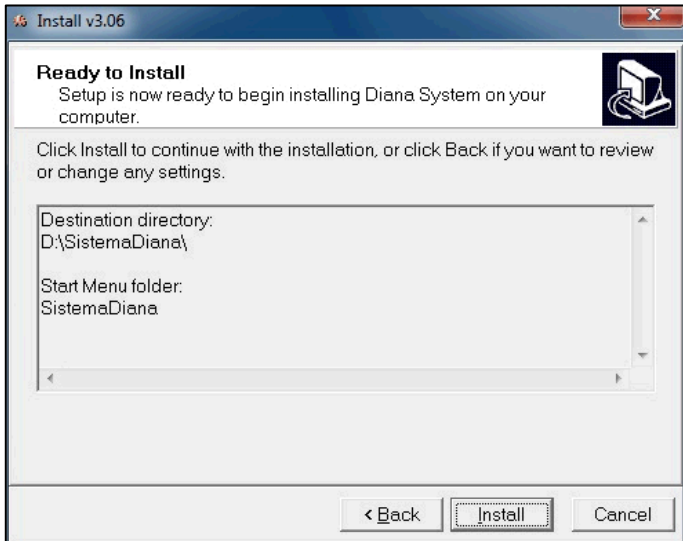
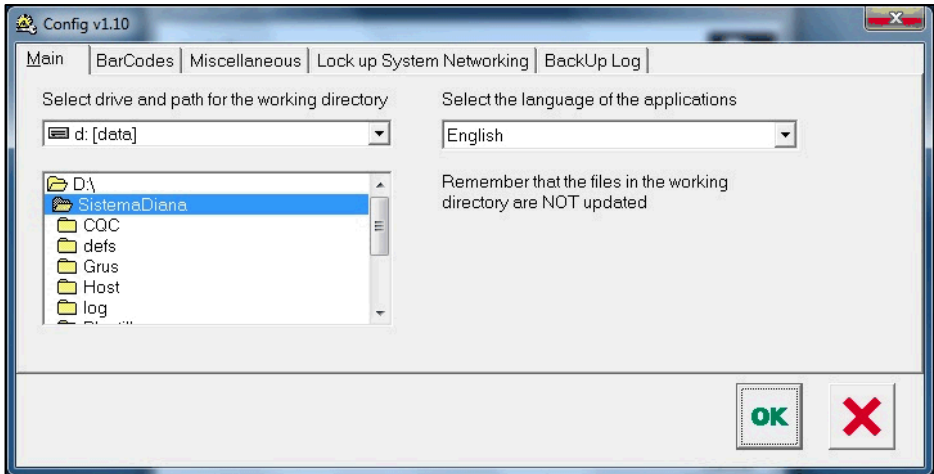
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
2. The "Welcome to the Diana System Setup Wizard" window will be displayed.  Picture 2 Select 'Next' button.	Next installation window appears.
3. "Select Components" window: Select the checkboxes of the components that are going to be installed. Check all of them for a full installation and click 'Next'  Picture 3 Note: Matrox drivers are not installed during the software installation. They are installed along with the O.S. installation.	Next installation window appears.

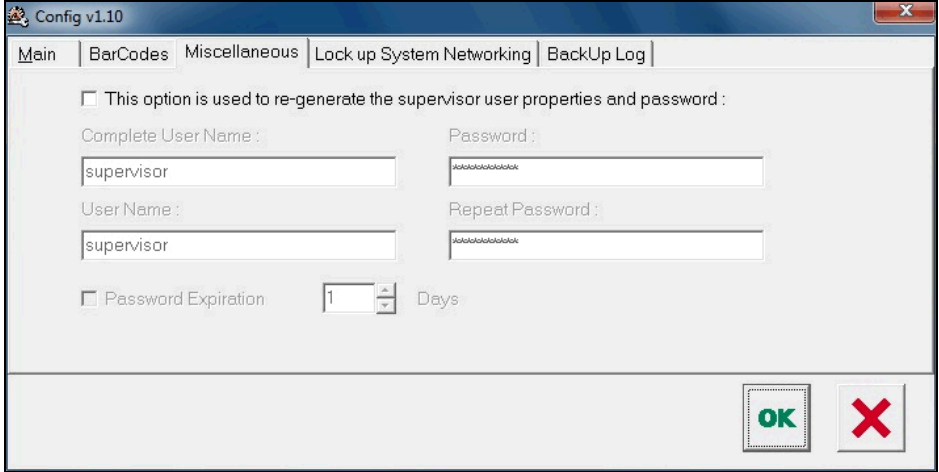
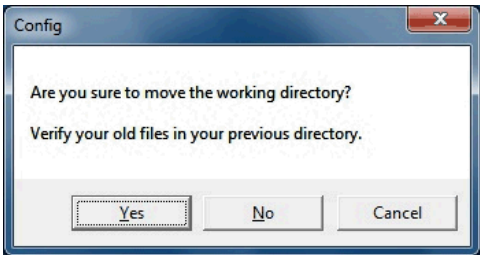
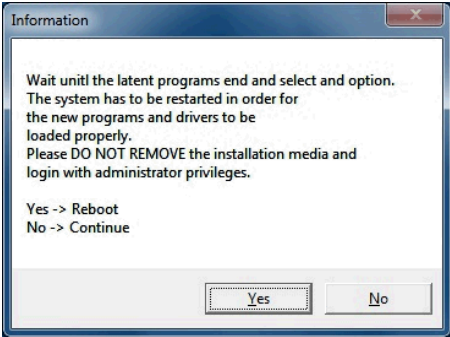
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
4. The next window indicates the destination directory where the software will be installed. The default path should not be changed. Picture 4 Select the 'Next' button to accept the installation folder "D:\SistemaDiana"	Select D:\SistemaDiana directory (default directory) Next installation window appears.
5. Select the menu folder where the program shortcuts will be created: Picture 5 Select the 'Next' button to accept the Start Menu Folder.	Start menu folder will be 'SistemaDiana'. Next installation window appears.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
6. The installation program shows the destination folder for the programs installation and the shortcuts destination too:  Picture 6 Click 'Install'. Note: If the software doesn't find 'password2.def' file in the working directory, an alert message will be prompted indicating that it will be created. Click 'ok'.	Software installation starts and 'config.exe' program will be executed.
7. Main tab shows working directory selection (D:\SistemaDiana) which is the folder where all the processed information will be saved. The language used for the software applications can also be chosen  Picture 7	--

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
8. Select the "Miscellaneous" tab (picture 8) and check the "This option is used ..." checkbox.  Picture 8 Set user name and password for the Supervisor user and press 'OK'.	A new config window appears.
9. If the working directory has to change its location, select 'Yes' in the Config window.  Picture 9	The installation process will end.
10. User will be prompted to reboot the computer without removing the installation media (Picture 10).  Picture 10 Select 'Yes'. Note: DO NOT extract the USB drive from the PC yet, as the installation process will be completed in the next reboot (although no new windows will appear).	--

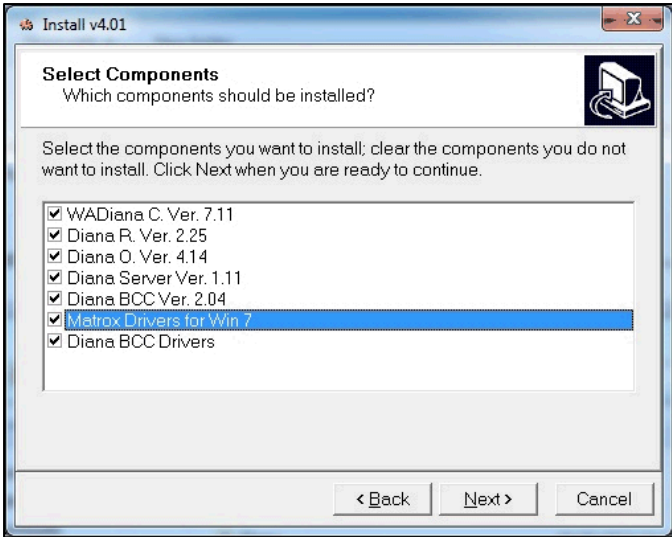
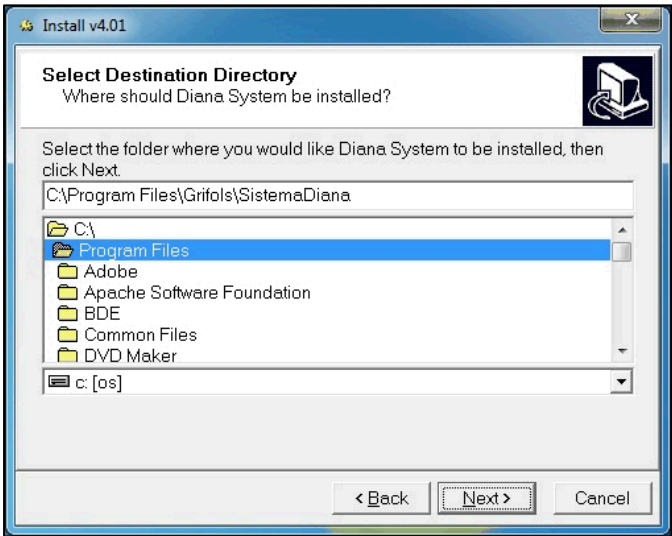
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
11. After PC reboot Diana System software installation is completed	--
12. In order for Windows users other than the Administrator to use any Diana System application, administrator privileges should be given to every application contained in the working directory (i.e. all the '.exe' files except for the 'Install.exe' file). With the O.S. write filter in <i>Disabled</i> status (refer to REGD-0010049): - Access D:\SistemaDiana folder - Right-click on the application and select 'Properties' - Click the 'Change settings for all users' button in the 'Compatibility' tab. - Check the 'Run this program as administrator' checkbox in the 'Privilege level' section. - Apply changes and click 'Ok' to exit the 'Properties' window. Reboot again the PC, making sure that the O.S. write filter is set back to Enabled status.	All Diana System applications can be executed with administrator privileges.
13. Install the desired Reagents Pack version following the instructions given in the corresponding INR/FSB	--
14. Restore the desired backed up data to the working directory (if any)	--

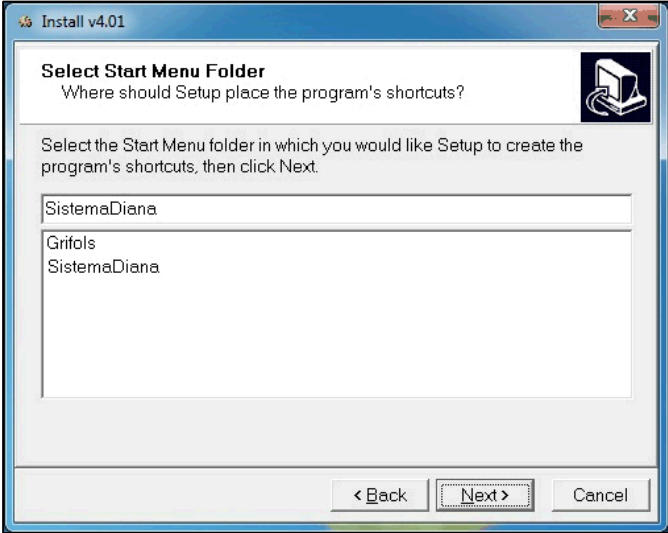
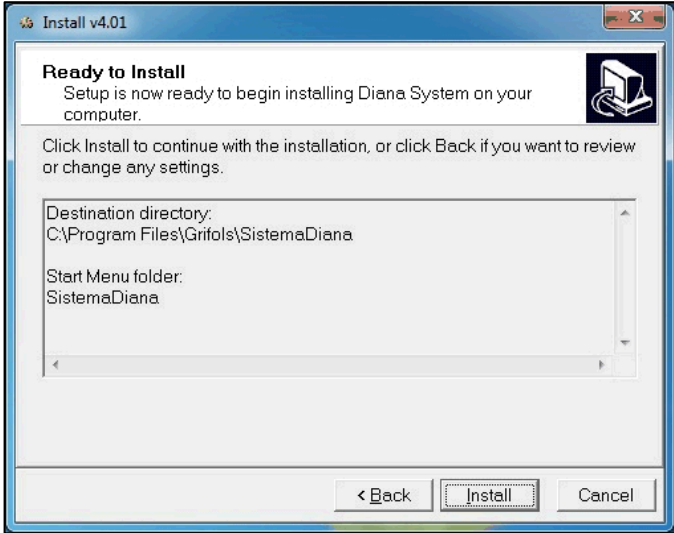
Case 2: Desktop PC (Windows 7)

Process	Expected Results
1. Copy the content of the Diana System software CD into the root of a USB drive and connect it to the PC. If no USB drive will be used during the installation, insert the software CD into a CDROM reader unit. If the system does not detect the installation program, then run Install.exe from the USB drive / CD contents	The installation wizard is started.
2. The "Welcome to the Diana System Setup Wizard" window will be displayed. Picture 11 Select 'Next' button.	Next installation window appears.

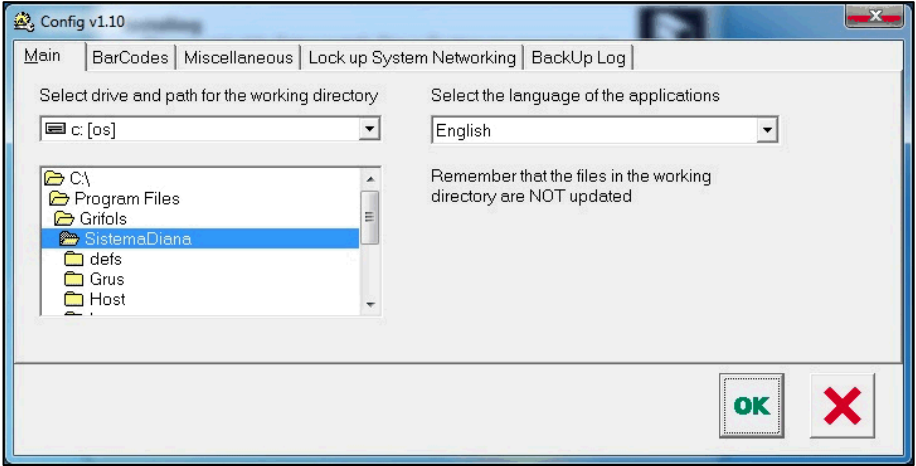
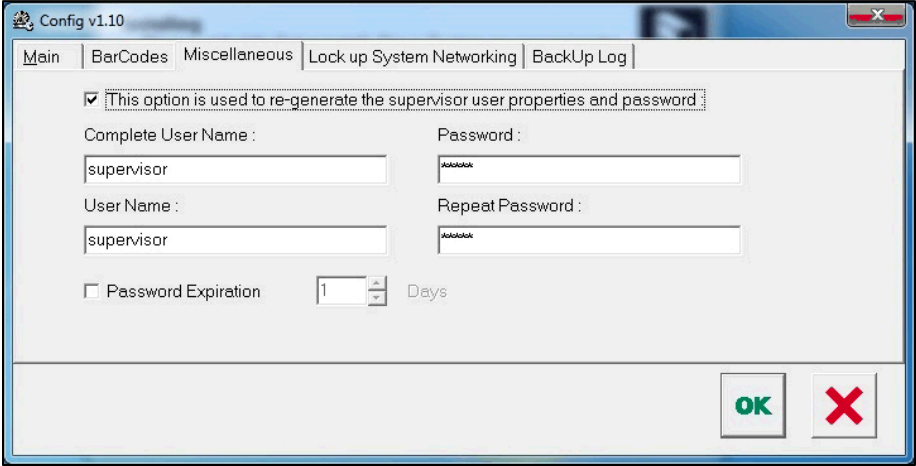
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
3. "Select Components" window: Select the checkboxes of the components that are going to be installed. Check all of them for a full installation and clic 'Next'  Picture 12 Note: Select Matrox drivers option only when no Matrox drivers are already installed in the computer.	Next installation window appears.
4. The next window indicates the destination directory where the software will be installed.  Picture 13 Select the 'Next' button to accept the installation folder	Next installation window appears.

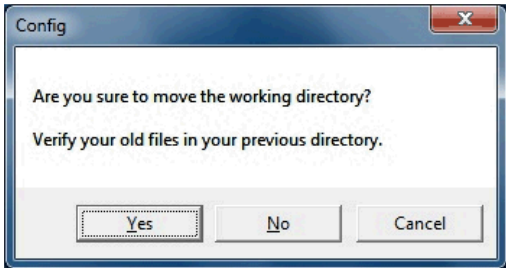
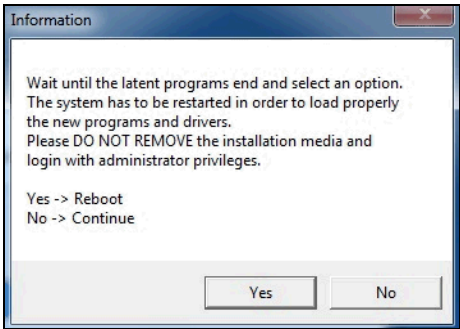
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
5. Select the menu folder where the program shortcuts will be created:  Picture 14 Select the 'Next' button to accept the Start Menu Folder.	Start menu folder will be 'SistemaDiana'. Next installation window appears.
6. The installation program shows the destination folder for the programs installation and the shortcuts destination too:  Picture 15 Click 'Install'. Note: If the software doesn't find 'password2.def' file in the working directory, an alert message will be prompted indicating that it will be created. Click 'ok'.	Software installation starts and 'config.exe' program will be executed.

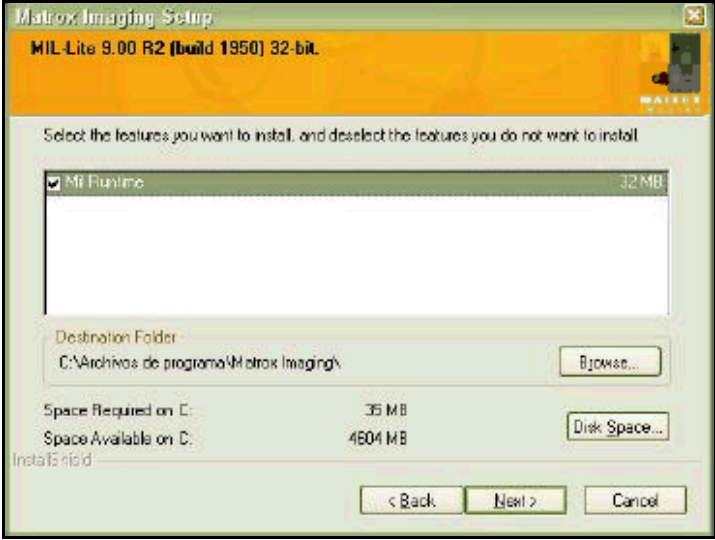
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
7. Main tab shows working directory selection, which is the folder where all the processed information will be saved. The language used for the software applications can also be chosen  Picture 16	--
8. Select the “Miscellaneous” tab (picture 17) and check the “This option is used ...” checkbox.  Picture 17 Set user name and password for the Supervisor user. Then press the 'OK' button.	A new config window appears.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
9. If the working directory has to change its location, select 'Yes' in the Config window.  Picture 18	The installation process will end.
10. User will be prompted to reboot the computer without removing the installation media (Picture 19).  Picture 19 Select 'Yes'. Note: DO NOT extract the installation USB drive / CD from the PC yet, as the installation process will be completed in the next reboot	--
11. After PC reboot Diana System software installation is completed	--
12. In case of using a desktop PC running under Windows 7: In order for Windows users other than the Administrator to use any Diana System application, administrator privileges should be given to every application contained in the working directory (i.e. all the '.exe' files except for the 'Install.exe' file). With the O.S. write filter in <i>Disabled</i> status (refer to REGD-0010049): - Access the working directory folder - Right-click on the application and select 'Properties' - Click the 'Change settings for all users' button in the 'Compatibility' tab. - Check the 'Run this program as administrator' checkbox in the 'Privilege level' section. - Apply changes and click 'Ok' to exit the 'Properties' window. Reboot again the PC, making sure that the O.S. write filter is set back to <i>Enabled</i> status.	All Diana System applications can be executed with administrator privileges.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

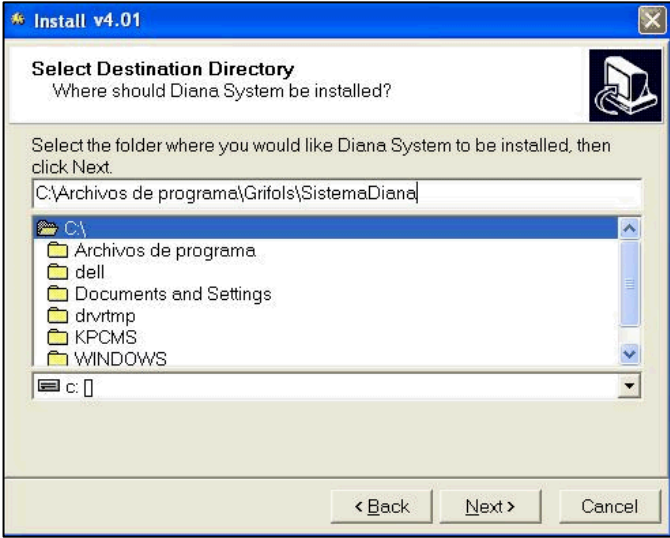
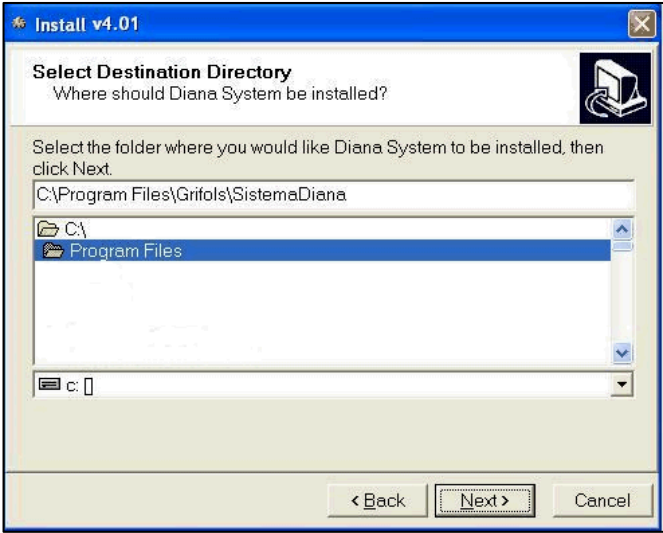
Process	Expected Results
13. If the Matrox drivers have been selected to be installed in step 3, the installation wizard continues with the following window:  Picture 20 Select 'Next'	Next installation window appears.
14. Select the 'Morphis' checkbox and click 'Next'	Next installation window appears.
15. Click on 'I agree' for the licence terms	Next installation window appears.
16. When the process is finished, you will be asked to restart the computer. Select 'Yes, I want to restart computer' option.	Next installation window appears.
17. Install the desired Reagents Pack version following the instructions given in the corresponding INR/FSB	--
18. Restore the desired backed up data to the working directory (if any)	--

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

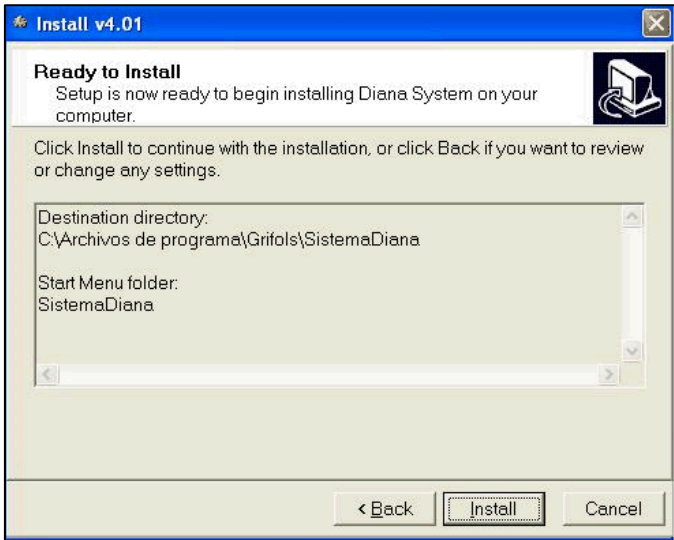
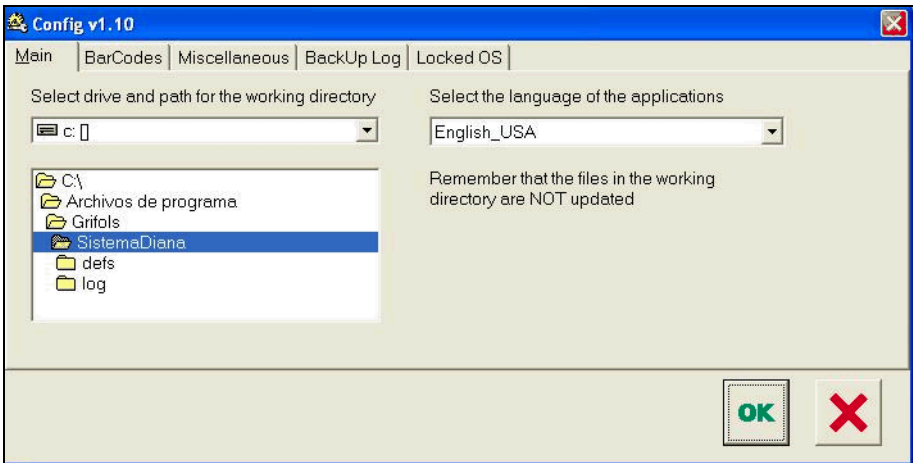
Case 3: Desktop PC (Windows XP)

Process	Expected Results
1. Insert the CD software into a CDROM reader unit If the system does not detect the installation program, then run <i>Install.exe</i> from the CD contents	The installation wizard is started.
2. The “Welcome to the Diana System Setup Wizard” window will be displayed. Picture 21 Select ‘Next’ button.	Next installation window appears.
3. “Select Components” window: Select the checkboxes of the components that are going to be installed. Check all of them for a full installation and clic ‘Next’ Picture 22 Note: Select Matrox drivers option only when no Matrox drivers are already installed in the computer.	Next installation window appears.

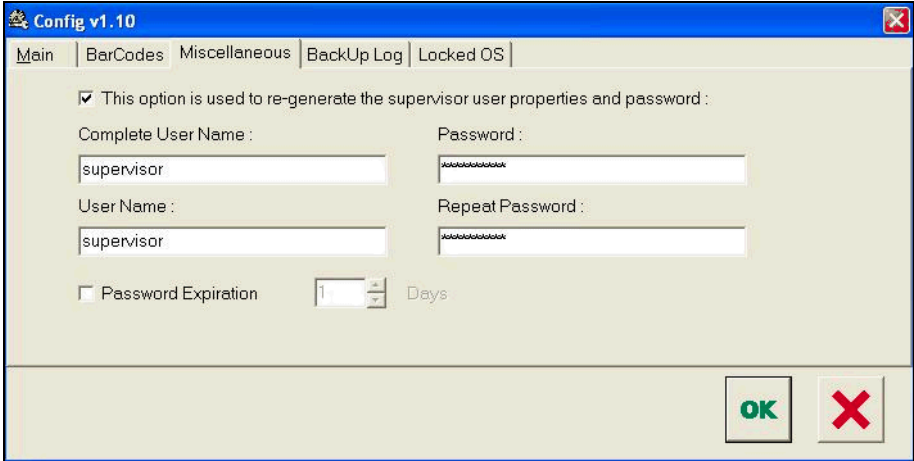
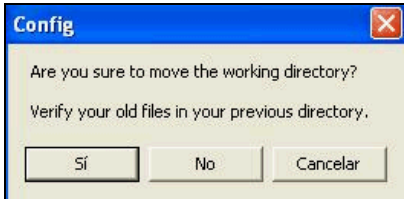
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
4. The next window indicates the destination directory where the software will be installed.  Picture 23 Select the 'Next' button to accept the installation folder	Next installation window appears.
5. Select the menu folder where the program shortcuts will be created:  Picture 24 Select 'Next' to accept the Start Menu Folder	Start menu folder will be 'SistemaDiana'. Next installation window appears.

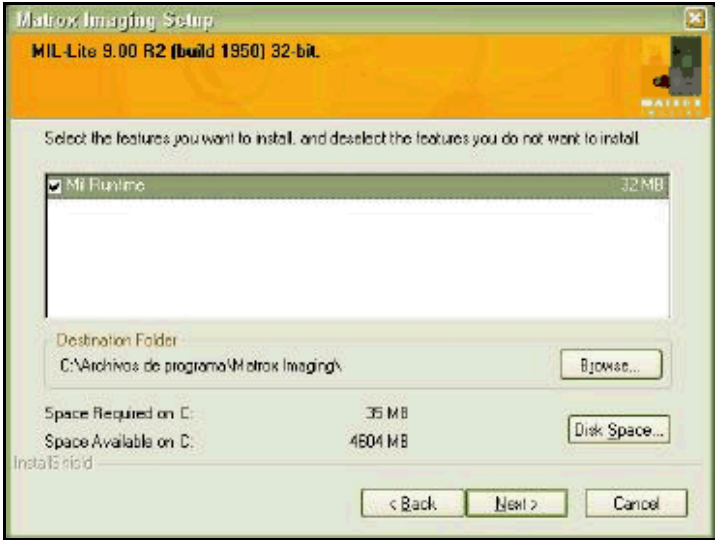
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
6. The installation program shows the destination folder for the programs installation and the shortcuts destination too:  Picture 25 Click 'Install'. Note: If the software doesn't find 'password2.def' file in the working directory, an alert message will be prompted indicating that it will be created. Click 'ok'.	Software installation starts and 'config.exe' program will be executed.
7. Main tab shows working directory selection, which is the folder where all the processed information will be saved. The language used for the software applications can also be chosen  Picture 26	--

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
8. Select the "Miscellaneous" tab (picture 29) and check the "This option is used ..." checkbox.  Picture 27 Set user name and password for the Supervisor user. Then press the 'OK' button.	A new config window appears.
9. If the working directory has to change its location, select 'Yes' in the Config window.  Picture 28	The installation process will end.
10. User will be prompted to reboot the computer without removing the installation media (Picture 31).  Picture 29 Select 'Yes'. Note: DO NOT extract the installation CD from the PC yet, as the installation process will be completed in the next reboot	--

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
11. After PC reboot Diana System software installation is completed	
12. If the Matrox drivers have been selected to be installed in step 3, the installation wizard continues with the following window:  Picture 30 Select 'Next'	Next installation window appears.
13. Depending on the model of the image acquisition board installed on the computer, select the Meteor II or Morphis checkbox and click 'Next'.	Next installation window appears.
14. Click on 'I agree' for the licence terms Note: When using Windows XP with a Matrox Meteor II board, a message saying 'This update is only compatible with Morphis board' will appear at the end of the process. Accept it to continue (no problems should appear once the installation process is finished)	Next installation window appears.
15. When the process is finished, you will be asked to restart the computer. Select 'Yes, I want to restart computer' option.	--
16. Install the desired Reagents Pack version following the instructions given in the corresponding INR/FSB	--
17. Restore the desired backed up data to the working directory (if any)	--

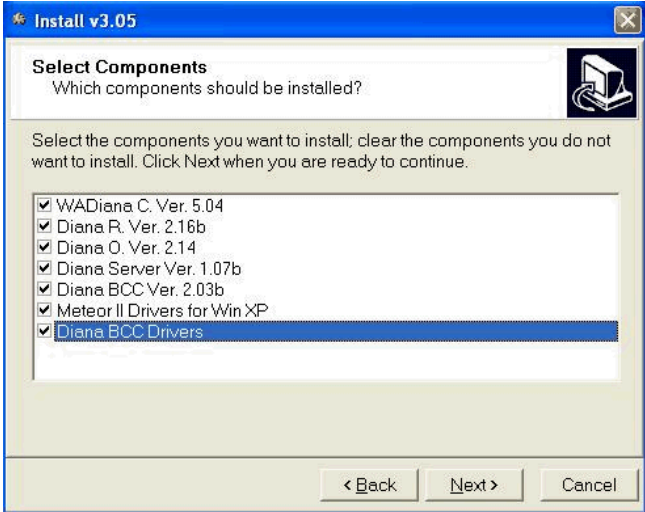
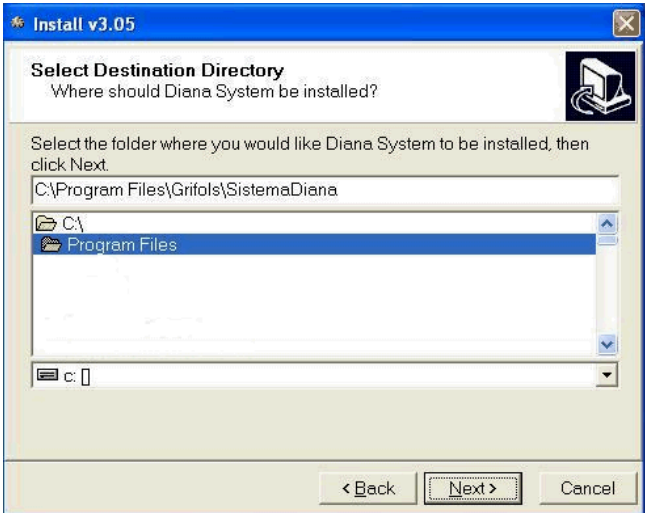
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

PROCEDURE 2: SOFTWARE VERSION 3.2

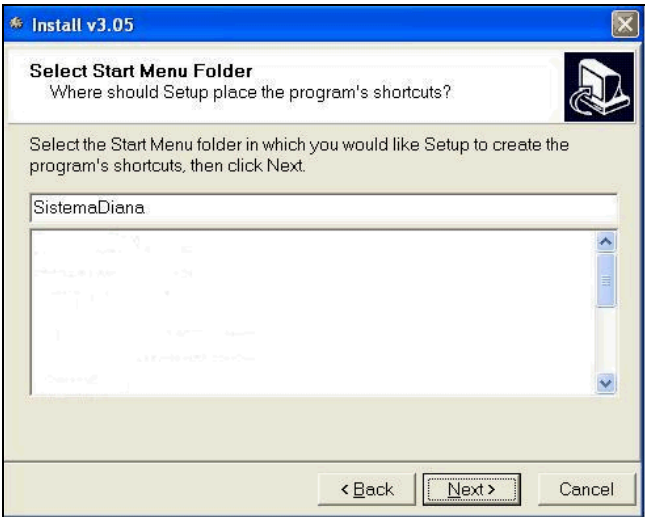
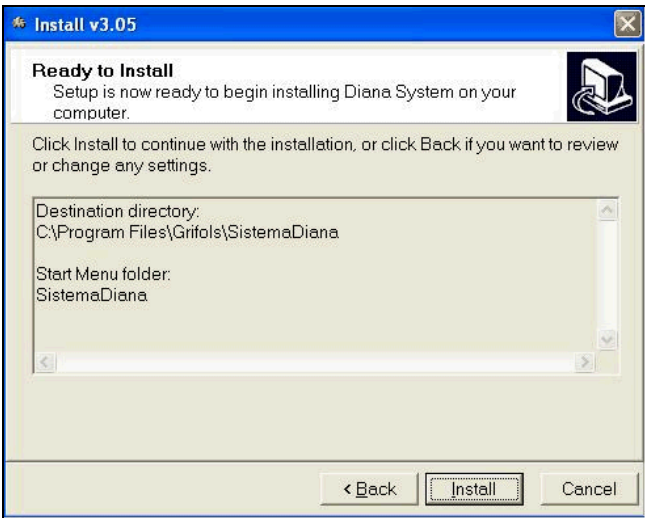
Process		Expected Results
1. Insert the CD software into a CDROM reader unit.		The installation program runs automatically.
2. If the operating system is not configured properly to auto run the installation program, open an explorer (file management) window and find the CDROM reader unit and run manually the 'install.exe' application. Picture 1		The installation program runs.
3. Click on 'Next' button or press Alt+N. Picture 2		Next installation window appears.

Confidential Information. It is not allowed the total or partial reproduction of this document without prior authorization.

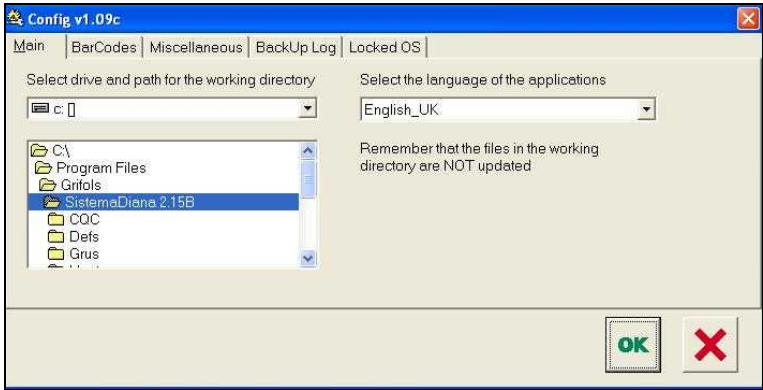
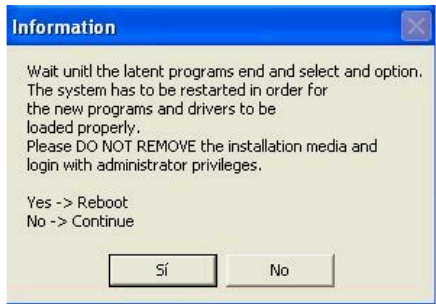
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
4. Check the checkboxes to install the appropriate software to work with the instrument and handle the information. Click 'Next'.  Picture 3 NOTE: Matrox drivers must be installed along with the software. Do not install the drivers separately from the rest of the software, as wrong files may be loaded in the PC.	The next installation window appears.
5. Select the destination directory where the software will be installed. Click 'Next'.  Picture 4	Next installation window appears.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
6. Select the menu folder to create the shortcuts to the programs. Click next.  Picture 5	Next installation window appears.
7. The installation program shows the destination folder for the programs installation and the shortcuts destination too. Click Install.  Picture 6	Software installation starts (Picture 7) and config.exe' program executes.
8. Select from the Main tab the working directory, the place where all the processed information will be saved and the language to be used by the applications. Select from the BarCodes tab the barcodes you'll want to use with the instrument. Select Miscellaneous tab and check the checkbox to create the privileged user of the software. You have to enter a username and a real complete name. The password has to be longer than 4 characters and is case sensitive. If apply select an expiration days for this password. (You'll be prompted to enter the actual password and you'll have to change for another one different).	The first installation process will end and user will be prompted to reboot the computer.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

Process	Expected Results
Select the BackUp Log tab to activate the log backup system to ensure all the activity with the software will be able to be traced. Also select the drive and path where it will be stored. Ensure that you have permissions over this directory and you have free space (at least 5MB. You need more than this to run properly the applications). Select the Locked OS to configure the Closed Operating System. (only W2000) Click the OK button.  Picture 7	
9. If it's the first time you're installing the software click yes. If you have the software installed previously and you're installing Matrox drivers click yes. If you have the software installed previously and you're installing programs you can click no. IMPORTANT: If you click yes, DO NOT REMOVE the software CDROM from the CDROM reader unit.  Picture 8	If click yes, the system reboots. If click no, the installation program ends.
10. After computer reboot (and maybe user successful password entry) you'll see the second installation process. The Borland database engine software installation.	Installation process ends.
11. Restore the desired backed up data to the working directory (if any)	--

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000777		
	Version:	6.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software installation			

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	<i>Document REGD-0000777_v2 replaces document 57004V0M000.</i>	Document reference changed.
3.0	2011/05/04	Added 'previous version' 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
3.0	2011/05/04	Updated the Document Change History: Date format changed to YYYY/MM/DD Dates changed to document version effective date.	Adoption of ISO 8601 date format to avoid ambiguity.
3.0	2011/05/04	ID-GELSTATION and PROVUE erased from the title.	Service manual intended for Grifols market.
3.0	2011/05/04	<i>Document REGD-0000777_v3 replaces document REGD-0000777_v2.</i>	Document version changed.
3.0	2011/05/04	Removed Meteor I from 'equipment' and 'procedure' sections.	Meteor I no longer available.
3.0	2011/05/04	Pictures updated.	Pictures updated with Diana System software version 3.2 installation.
4.0	May 2011	<i>Document REGD-0000777_v4 replaces document REGD-0000777_v3.</i>	Document version changed.
4.0	May 2011	<i>Replaced commercial name by project name.</i>	Delete commercial name.
4.0	May 2011	<i>Header title has changed.</i>	Previous header title has been substituted by project and document title.
5.0	December 2014	Updated installation procedures	New Matrox 4Sight GP computer and software Sistema Diana v4.1
6.0	April 2016	Updated installation procedures	Updated procedures related to latest software version (Sistema Diana v4.2.0). Added step to run software with administrator privileges.