

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000778		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software uninstallation			

DG-57 Service Manual

SOFTWARE UNINSTALLATION

Procedure description
This procedure describes how to uninstall the software for both the Matrox 4Sight GP computer and a standard desktop computer. The uninstallation process does not uninstall the image acquisition board drivers. Refer to Image acquisition board uninstallation document (REGD-0000780) <i>Note: Document REGD-0000778_v5 replaces document REGD-0000778_v4.</i>

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000778		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software uninstallation			

EQUIPMENT

1. Matrox 4Sight GP computer / Standard computer fulfilling the specified requirements
--

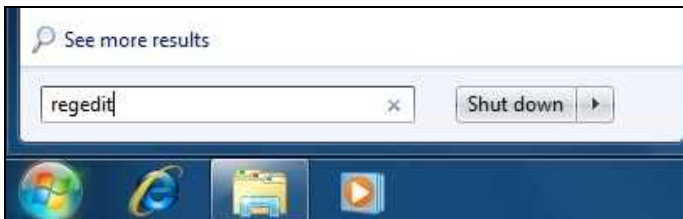
SET UP

1. PC turned on.

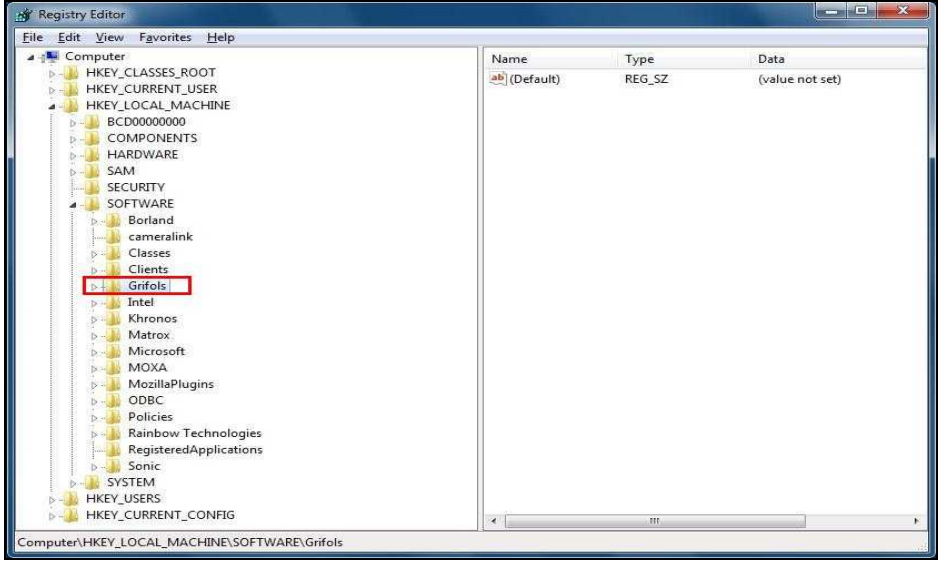
Uninstallation Procedures

- **PROCEDURE 1: PC MATROX 4SIGHT GP**
- **PROCEDURE 2: DESKTOP PC (WINDOWS 7)**
- **PROCEDURE 3: DESKTOP PC (WINDOWS XP)**

PROCEDURE 1: PC MATROX 4SIGHT GP

Process	Expected Results
1. Make a manual copy of the following listed folders (They will be located in 'D:/SistemaDiana' directory): a. Tandas b. Usades c. Folder selected to save the traceability. d. Reference image	---
2. Make sure that 'C:' drive is protected (refer to REGD-0010049)	Write Filter Status is Enabled
3. Select 'D:/SistemaDiana' folder and delete it.	---
4. Erase the specific configurations stored in the Windows registry system: Open the Registry Editor by clicking the Start button, typing 'regedit' into the search box, and then pressing Enter.  Picture 1	---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000778		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software uninstallation			

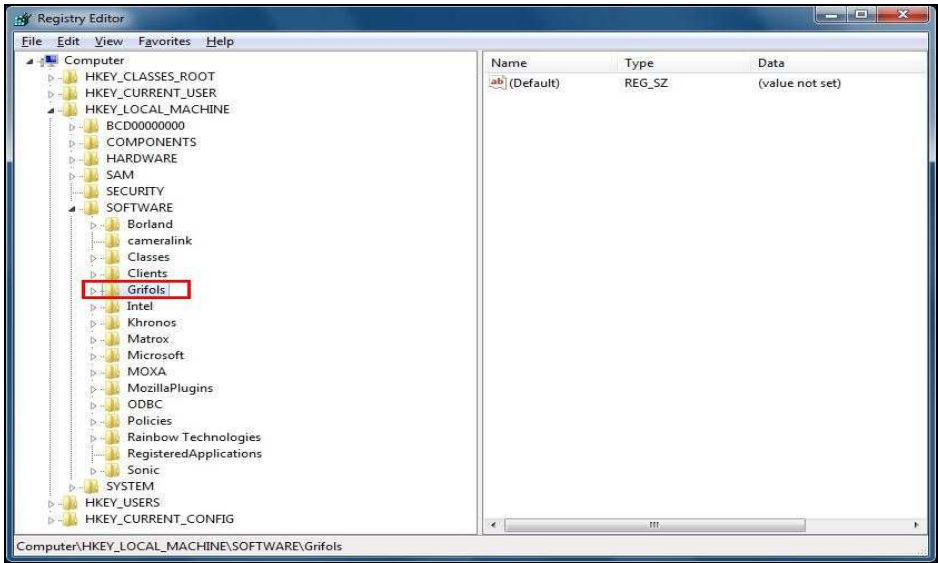
Process	Expected Results
5. The following window will appear:  Picture 2 Search on the tree: HKEY_LOCAL_MACHINE → Software → Grifols Delete Grifols registry.	---
6. Restart the system	---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000778		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software uninstallation			

PROCEDURE 2: DESKTOP PC (WINDOWS 7)

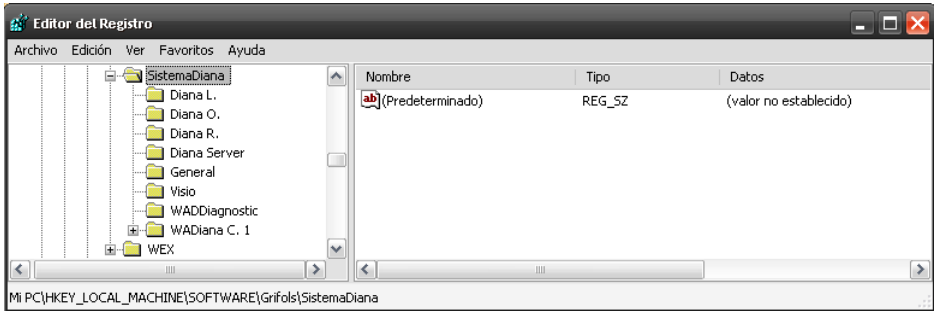
Process	Expected Results
1. Make a manual copy of the following listed folders (They will be located in 'D:/SistemaDiana' directory): - Tandas - Usades - Folder selected to save the traceability. - Reference image 2. Access Start → Control Panel → Uninstall a Program in 'Programs' section. In the programs list, Diana System software can be found. Select it and click on '<i>Uninstall / Change</i>' button. Note: The following message could appear: Picture 3 Once accepted, the uninstallation process will continue.	--- The system asks for confirmation.
3. Click 'Yes' to accept the confirmation message.	The system uninstalls the software.
4. Select your working directory folder and delete it. For instance (default directory) delete C:\Program Files\Grifols\SistemaDiana .	The system will erase all the information related to software.
5. Erase the specific configurations stored in the Windows registry system: Open the Registry Editor by clicking the Start button, typing 'regedit' into the search box, and then pressing Enter. Picture 4	The regedit main window will appear (Picture 5).

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000778		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software uninstallation			

Process	Expected Results
6. The following window will appear:  Picture 5 Search on the tree: HKEY_LOCAL_MACHINE → Software → Grifols Delete Grifols registry.	The key and subkeys will be erased.
7. Restart the system	---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000778		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software uninstallation			

PROCEDURE 3: DESKTOP PC (WINDOWS XP)

Process	Expected Results
1. Access Start → Control Panel → Add / Remove programs.	---
2. In the programs list, instrument software can be found. Select it and click on 'Modify / Remove' button.	The system asks for confirmation.
3. Click 'Yes' button.	The system uninstalls the software.
4. During the software uninstallation, the system asks for confirmation. Accept it to proceed.	The system uninstalls the BDE software using its own uninstaller program.
5. All the software is now erased but not the data created by it. So in order to delete or copy to another machine, open an explorer and find your work directory and there you'll find all this data. In order to delete everything select the working directory and press <i>Delete</i> keyboard button.	The system will erase all the information related to software.
6. The last part of any installation is to erase the specific configurations. These configurations are stored in the 'Windows registry system'. So click on Start, Run. Type 'regedit' and press Enter.	The regedit main window will appear.
 Picture 6	
7. On the left part of the main window, search on the tree HKEY_LOCAL_MACHINE, SOFTWARE, Grifols.	---
8. Delete all 'Grifols entry by pressing <i>Delete</i> keyboard button.	The key and subkeys will be erased.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000778		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	28-Apr-2016
Title:	SM DG57. Software uninstallation			

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	<i>Document REGD-0000778_v2 replaces document 57005V0M000.</i>	Document reference changed.
3.0	May 2011	<i>Document REGD-0000778_v3 replaces document REGD-0000778_v2.</i>	Document version changed.
3.0	May 2011	Document title.	Commercial reference substituted by project reference.
3.0	May 2011	Header title.	Increase information to easy document search.
3.0	May 2011	Picture updated.	Updated with latest software configuration.
4.0	December 2014	Updated uninstallation procedures	New Matrox 4Sight GP computer and software Diana System v4.1
5.0	April 2016	Updated uninstallation procedure	Added checking step

Confidential Information. It is not allowed the total or partial reproduction of this document without prior authorization.

Print date:	28-Apr-2016	Page	8 of 8
--------------------	-------------	-------------	--------