

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000779		
	Version:	4.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board installation			

DG-57 Service Manual

IMAGE ACQUISITION BOARD INSTALLATION

Procedure description
This procedure describes how to install the image acquisition board software drivers, along with the physical installation of the Matrox Morphis board in the Matrox 4Sight GP computer Note: Document REGD-0000779_v4 replaces document REGD-0000779_v3

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000779		
	Version:	4.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board installation			

EQUIPMENT

1. Matrox Meteor II / Matrox Morphis Image acquisition board
--

SET UP

1. PC turned off

INSTALLATION PROCEDURES

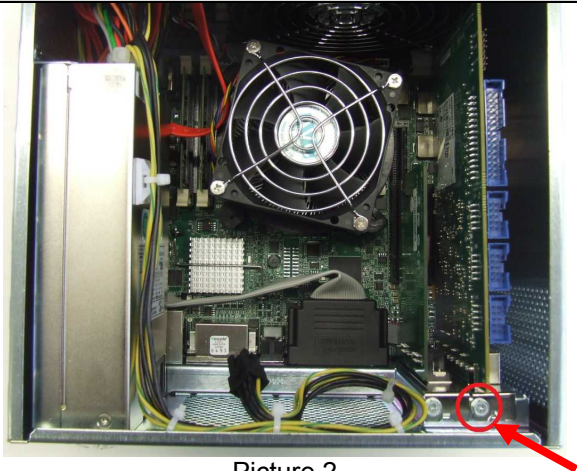
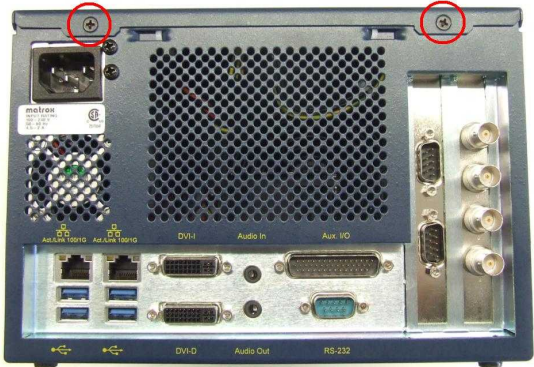
- **PROCEDURE 1: MATROX 4SIGHT GP COMPUTER**
 - Physical installation
 - Drivers installation
- **PROCEDURE 2: STANDARD COMPUTER**
 - CASE 1 MATROX MORPHIS (Diana System versions ≥ 3.3)
 - CASE 2 MATROX METEOR II (only Windows XP)

PROCEDURE 1: MATROX 4SIGHT GP COMPUTER

Physical installation

Process	Expected Results
1. Once an image acquisition board has been uninstalled from the computer (see REGD-0000780), insert the new image board in the free PCI slot of the computer as shown in the following picture: A close-up photograph showing a person's hand, wearing a white sleeve, carefully inserting a green PCI image acquisition board into a slot within a computer chassis. The board has blue connectors on its edge. Other components like a fan and other circuit boards are visible in the background. Picture 1	--
2. Insert the screw that holds the image board to the chassis of the PC	--

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000779		
	Version:	4.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board installation			

Process	Expected Results
 Picture 2	
3. Slide the top cover forward (from rear to front part of the computer) and screw the two black screws located at the back of the PC  Picture 3	--

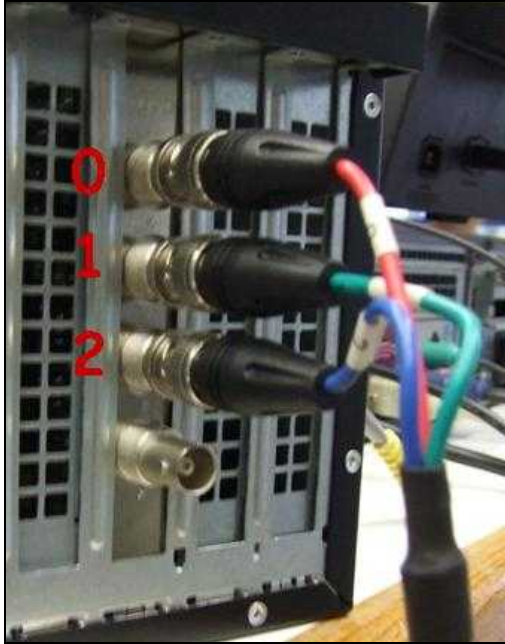
Drivers installation

Process	Expected Results
1. Matrox Morphis drivers are installed when installing the OS (Windows Embedded Standard 7). Refer to REGD-0016476 <i>Windows Embedded Standard 7 Multi-language Installation</i> for more information.	--
2. Once the new board is installed, perform the fine adjustment of the reader camera optics (points 21 to 24 in <i>REGD-0000812 Reader Camera Adjustment</i>).	--
3. Create a new reference image following the procedure described in <i>REGD-0000814 Reference Image creation</i>	--

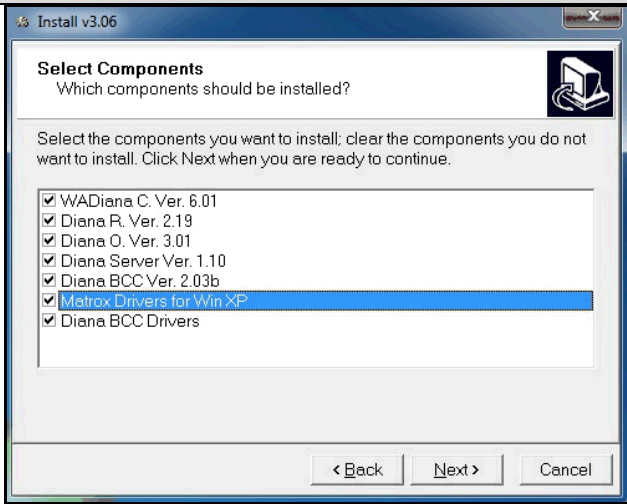
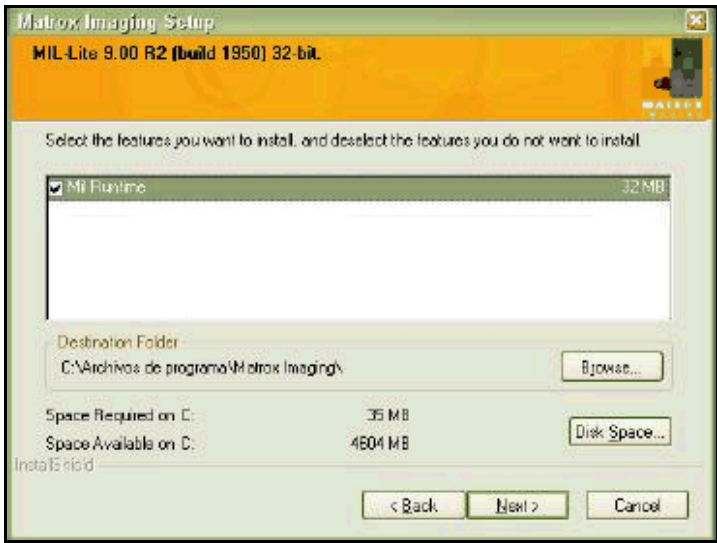
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000779		
	Version:	4.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board installation			

PROCEDURE 2: STANDARD COMPUTER

CASE 1: MATROX MORPHIS (Diana System versions ≥ 3.3)

Process	Expected Results
1. With the computer switched off, insert the Matrox Morphis board into one of the PC's free PCI slots and connect the video cable checking its labels as shown in the following picture:  Picture 1 Note: This procedure will vary depending on the PC brand and model	--
2. Turn on the computer and reject any automatic hardware installation for Matrox board.	
3. Insert the Diana System SW CD and run the Install program	The 'Select components' window is shown
4. Select the 'Matrox Drivers for WinXP' and press 'Next'. (Picture 2). Do not check the rest of options if Diana System software is not required to be installed	The installation will start.

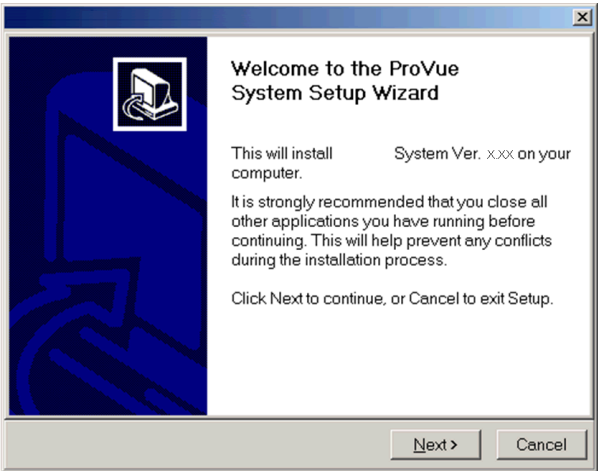
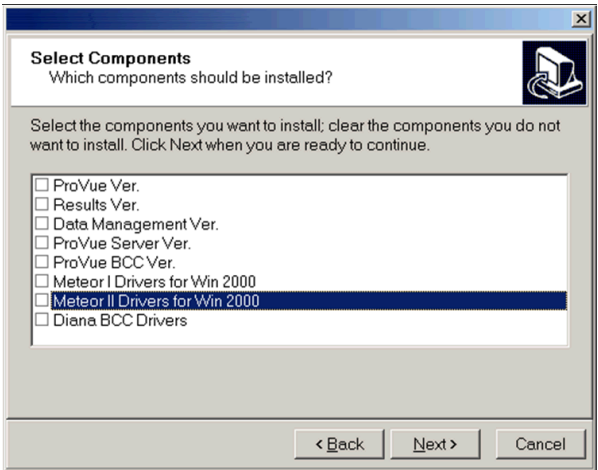
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000779		
	Version:	4.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board installation			

Process		Expected Results
 Picture 2		
5.	When the installation ends the system will ask to reboot the system. Click on 'yes'.	The system will reboot.
6.	The installation wizard continues with the following window:	Next installation window appears.
 Picture 3 Select 'Next'		
7.	Select the 'Morphis' checkbox and click 'Next'	Next installation window appears.
8.	Click on 'I agree' for the licence terms	Next installation window appears.
9.	When the installation ends the system will ask to reboot the system. Click on 'yes'.	The system will reboot.
10.	Once the new board is installed, perform the fine adjustment of the reader camera optics (points 21 to 24 in <i>REGD-0000812 Reader Camera Adjustment</i>).	--
11.	Create a new reference image following the procedure described in <i>REGD-0000814 Reference Image creation</i>	--

Confidential Information. It is not allowed the total or partial reproduction of this document without prior authorization.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000779		
	Version:	4.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board installation			

CASE 2: MATROX METEOR II (only Windows XP)

Process	Expected Results
1. With the computer switched off, insert the Matrox Meteor II board into one of the PC's free PCI slots. Note: This procedure will vary depending on the PC brand and model	--
2. Just after the system reboots it will detect the new hardware installed. Press the 'cancel' button.	The system will continue without the drivers.
3. Insert the Instrument software CDROM and it will automatically run the 'install.exe' file. Should this not be the case, execute it manually from software CD.	The installation program will appear.
4. Press 'next'. (Picture 1)  Picture 1	The next installation step appears.
5. Select the meteor drivers to install depending on the installed hardware and press 'next'. (Picture 2)  Picture 2 Note: For Windows XP, only "Meteor II Drivers for Win 2000" drivers are available and compatible.	The installation will start.
6. When the installation ends the system will ask to reboot the system. Click on 'yes'.	The system will reboot.

Confidential Information. It is not allowed the total or partial reproduction of this document without prior authorization.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000779		
	Version:	4.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board installation			

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
1.0	2005/02/17	NA	1st version of the document
2.0	2011/04/28	Updated pictures 1 and 2. Removed information about versions.	The pictures showed versions that are not current.
2.0	2011/04/28	Step 4: Added information about compatibility between Windows XP and Meteor II.	It was not clear the compatibility between Windows XP and "Meteor II drivers for Windows 2000".
2.0	2011/04/28	Updated the Document Change History: - Date format changed to YYYY/MM/DD - Dates changed to document version effective date. - Added the first version reference.	Adoption of ISO 8601 date format to avoid ambiguity.
2.0	2011/04/28	ID-GELSTATION and PROVUE erased from the title.	Service manual intended for Grifols market.
3.0	May 2011	Document REGD-0000779_v3 replaces document REGD-0000779_v3	Document version changed.
3.0	May 2011	Document title.	Commercial reference substituted by project reference.
3.0	May 2011	Header title.	Increase information to easy document search.
4.0	December 2014	Updated installation procedures	New Matrox 4Sight GP computer