

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000780		
	Version:	7.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board uninstallation			

DG-57 Service Manual

IMAGE ACQUISITION BOARD UNINSTALLATION

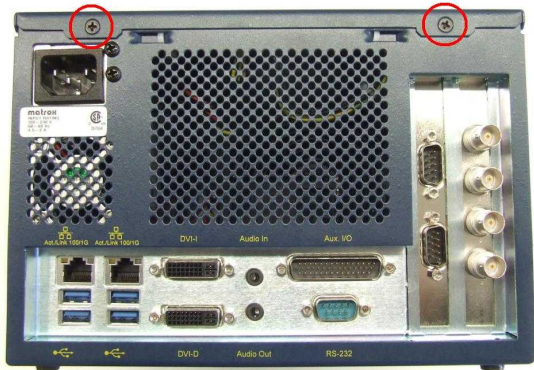
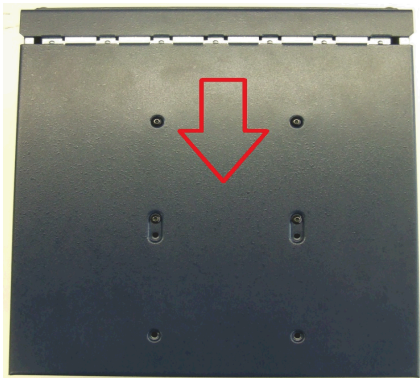
Procedure description
This procedure describes how to uninstall the image acquisition board software drivers for both Morphis and Meteor II boards. The physical uninstallation of the Matrox Morphis board from the Matrox 4Sight GP computer is also detailed. Note: Document REGD-0000780_v7 replaces document REGD-0000780_v6

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000780		
	Version:	7.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board uninstallation			

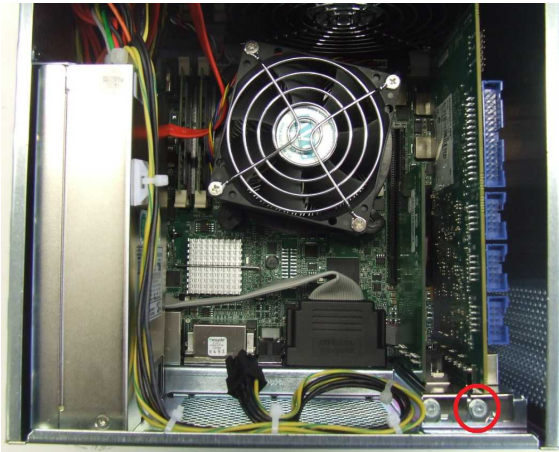
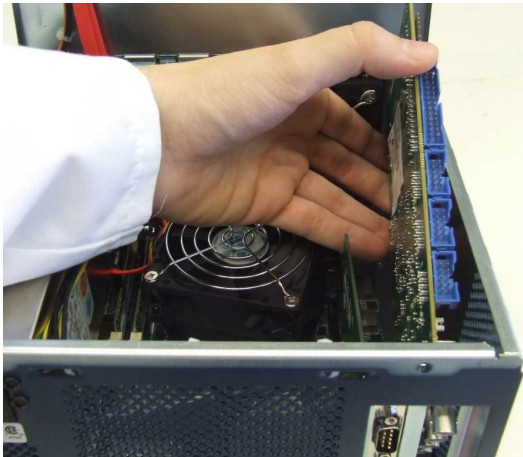
UNINSTALLATION PROCEDURES

- **PROCEDURE 1: PHYSICAL UNINSTALLATION (MATROX 4SIGHT GP COMPUTER)**
- **PROCEDURE 2: SOFTWARE UNINSTALLATION**
 - **CASE 1: MATROX MORPHIS / METEOR II (MIL 9.0 drivers)**
 - **CASE 2: MATROX METEOR II (MIL 6.1 and 7.5 drivers)**
 - MS Windows 2000 (driver version 6.1)
 - MS Windows 2000 (driver version 7.5)
 - MS Windows XP (driver version 7.5)

PROCEDURE 1: PHYSICAL UNINSTALLATION (MATROX 4SIGHT GP COMPUTER)

Process	Expected Results
1. Unscrew the two black screws located at the back of the PC.  Picture 1	--
2. Slide the top cover backwards (picture 2) and lift it so that it remains in vertical position (picture 3) allowing the access to the internal components.  Picture 2	--

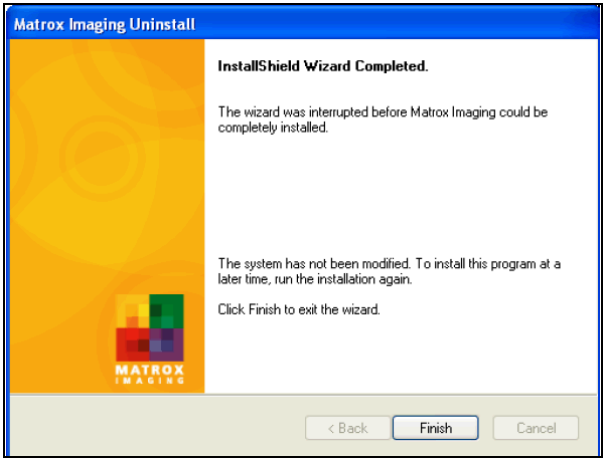
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000780		
	Version:	7.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board uninstallation			

Process	Expected Results
 Picture 3	
3. Remove the screw that holds the image board to the chassis of the PC.  Picture 4	--
4. Unplug the image board from the Main board and remove it by using the free space in the right part of the PC, as shown in Picture 5:  Picture 5	--

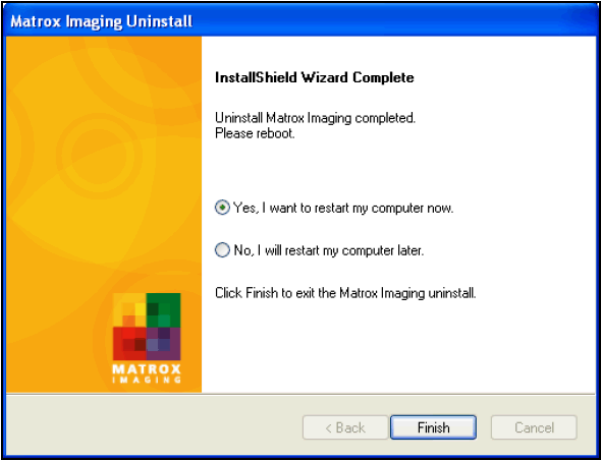
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000780		
	Version:	7.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board uninstallation			

PROCEDURE 2: SOFTWARE UNINSTALLATION

CASE 1: MATROX MORPHIS / METEOR II (MIL 9.0 drivers)

Process	Expected Results
1. Access Control Panel and then click on 'Add or Remove Programs' (Windows XP) or 'Uninstall a Program' (Windows 7). Select 'Matrox Imaging' to remove:  Picture 1	The uninstallation wizard starts
2. Uninstallation process will start. Select 'Yes' for the confirmation message and the process will end:  Picture 2 Select 'Finish'	The uninstallation of the software components will begin

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000780		
	Version:	7.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board uninstallation			

Process	Expected Results
3. The components will be automatically uninstalled and, after accepting another confirmation message, the following window will appear:  Picture 3 Select 'Yes, I want to restart my computer now' and then 'Finish'	The system will reboot
4. Access Windows registry and delete the 'MatroxDrvVersion' register from <i>HKEY_LOCAL_MACHINE\SOFTWARE\Grifols\Sistema Diana\General</i>  Picture 4	--

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000780		
	Version:	7.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board uninstallation			

CASE 2: MATROX METEOR II (MIL 6.1 and 7.5 drivers)

MS Windows 2000 (driver version 6.1)

Process	Expected Results
1. Once the Matrox Meteor board has been removed from the computer, open the Start Menu and go to Settings, Control Panel, Add/Remove Hardware.	A selection window will appear.
2. Select Uninstall/Remove component and click 'Next'.	A selection window will appear.
3. Select Uninstall component and click 'Next'.	A list of installed components will appear.
4. Select the Matrox Meteor PCI item and click 'Next'	A question will appear.
5. Select 'yes I wish to remove this device' and click 'Next'	A final window will appear.
6. Click on 'Finish'.	The drivers will be removed.
7. If there more than one Matrox Meteor board installed, repeat the process to remove them.	
8. Remove the following listed files if they are present, located in the Windows main directory (usually 'winnt') and inside of 'inf' directory: - " Meteor_II.inf " - " Meteor_II.pnf " - " Meteor.inf " - " Meteor.pnf " (example 'C:\winnt\inf\Meteor_II.inf')	--
9. Remove the following listed files if they are present, located in the Windows main directory (usually 'winnt') and inside of 'System32' and 'drivers' directory: - " corona.sys " - " meteor.sys " - " Mtxdma0.sys " - " Mtxdma1.sys " (example 'C:\winnt\System32\drivers\corona.sys')	--
10. Restart the computer. Cancel the possible messages from the plug and play system. Open a MS-DOS emulation session (Start->Run... "command.com").	A MS-DOS like window will appear.
11. Insert the software CDROM. Type the letter of the CDROM unit followed by a ':' sign and press enter. (example "D:")	Unit change.
12. Type "cd meteor\ver6_1\drivers\win2000".	Directory change.
13. Type "mildrv meteor_ii -remove" and click to the message box.	Remove registry links to meteor.
14. Type "mildrv mtxdma -remove" and click to the message box.	Remove registry links to meteor.
15. Type "mildrv meteor -remove" and click ok to the message box.	Remove registry links to meteor.
16. Open Windows registry and access HKEY_LOCAL_MACHINE\SOFTWARE\Grifols\SistemaDiana\General folder	--
17. Remove parameter MatroxDrvVersion if present	MatroxDrvVersion parameter is removed.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000780		
	Version:	7.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board uninstallation			

MS Windows 2000 (driver version 7.5)

Process	Expected Results
18. Once the Matrox Meteor board has been removed from the computer, right click on My Computer icon. Select "Properties" option.	My Computer properties window appears
19. Access to "Hardware" tab. Click on "Device Manager" button.	Device Manger window appears
20. Double click on "Matrox Imaging Adapter".	Meteor_II PCI Frame Grabber option appears
21. Right click on "Meteor_II PCI Frame Grabber". Select "Uninstall" option.	User is prompted to confirm driver deletion
22. Click on "Yes" to remove the drivers.	Driver is removed.
23. If there more than one Matrox Meteor board installed, repeat the process to remove them.	--
24. Open Windows registry and access <i>HKEY_LOCAL_MACHINE\SOFTWARE\Grifols\SistemaDiana\General</i> folder	--
25. Remove parameter MatroxDrvVersion if present	MatroxDrvVersion parameter is removed.

MS Windows XP (driver version 7.5)

Process	Expected Results
26. Once the Matrox Meteor board has been removed from the computer, right click on My Computer icon. Select "Properties" option.	My Computer properties window appears
27. Access to "Hardware" tab. Click on "Device Manager" button.	Device Manger window appears
28. Double click on "Matrox Imaging Adapter".	Meteor_II PCI Frame Grabber option appears
29. Right click on "Meteor_II PCI Frame Grabber". Select "Uninstall" option.	User is prompted to confirm driver deletion
30. Click on "Yes" to remove the drivers.	Driver is removed.
31. If there more than one Matrox Meteor board installed, repeat the process to remove them.	--
32. Open Windows registry and access <i>HKEY_LOCAL_MACHINE\SOFTWARE\Grifols\SistemaDiana\General</i> folder	--
33. Remove parameter MatroxDrvVersion if present	MatroxDrvVersion parameter is removed.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000780		
	Version:	7.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Image acquisition board uninstallation			

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
4.0	2005/05/25	Updated document 57007V0M001	<i>Document REGD-0000780_v4 replaces document 57007V0M001</i>
4.0	2005/05/25	Changed procedures step numbering and CD-ROM folder path for Windows NT. Add procedure for Windows 2000 and version drivers 7.5.	Windows NT requires changes of the path. New driver version 7.5 and Windows 2000 needs new uninstall procedures.
5.0	2011/04/28	Removed Procedures for MS Windows 98 and Windows NT.	Service manual intended for Windows 2000 and XP, not for Windows 98 and NT.
5.0	2011/04/28	Added Procedure for MS Windows XP.	The uninstall procedure for Windows XP was not present in the document.
5.0	2011/04/28	Updated the Document Change History: - Date format changed to YYYY/MM/DD - Dates changed to document version effective date.	Adoption of ISO 8601 date format to avoid ambiguity.
5.0	2011/04/28	ID-GELSTATION and PROVUE erased from the title.	Service manual intended for Grifols market.
6.0	May 2011	Document REGD-0000780_v6 replaces document REGD-0000780_v5.	Document version changed.
6.0	May 2011	Document title.	Commercial reference substituted by project reference.
6.0	May 2011	Header title.	Increase information to ease document search.
7.0	December 2014	Updated uninstallation procedures	New Matrox 4Sight GP computer