

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000781	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	23-May-2011
Title:	SM DG57-DG517. Networking.			Page 2 of 6

DG-57 / DG-517 Service Manual

NETWORKING

Procedure description
This procedure describes how to share the working directory data with other computers over a network. <i>Note: Document REGD-0000781_v2 replaces document 57008V0M000.</i>

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000781	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	23-May-2011
Title:	SM DG57-DG517. Networking.			Page 3 of 6

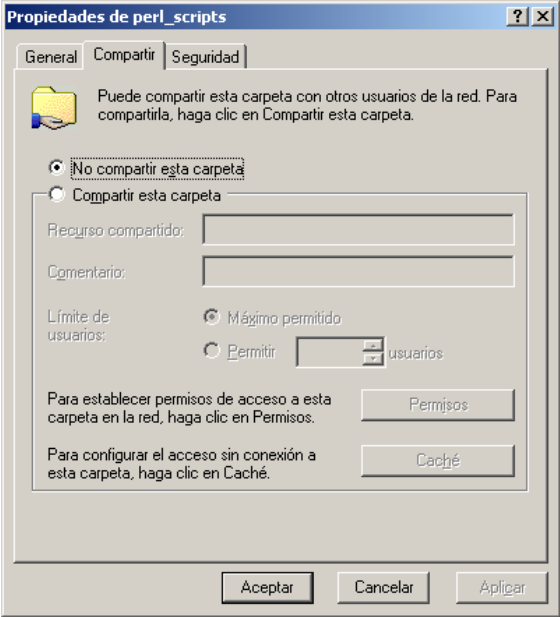
EQUIPMENT

1. Network card installed, configured and running.

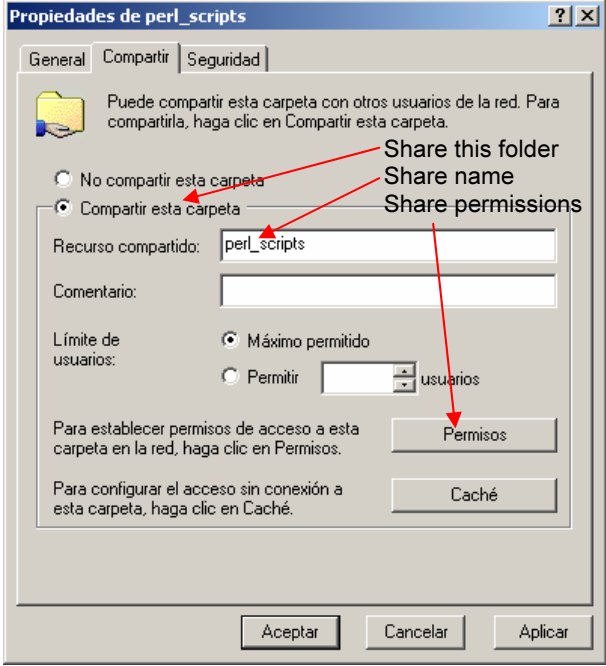
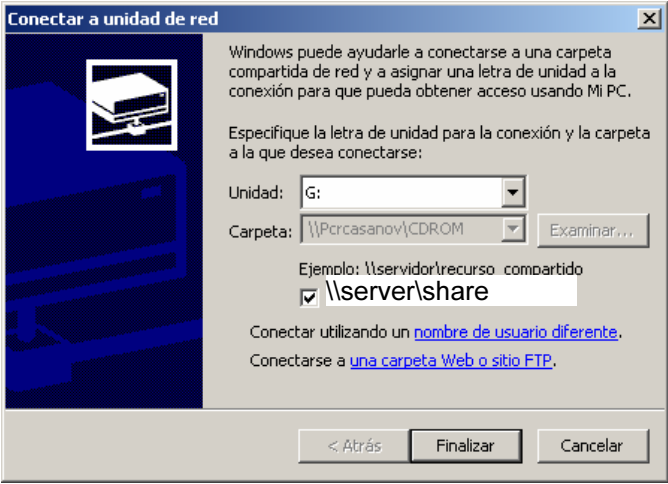
SET UP

1. PC turned on and connected to the network.

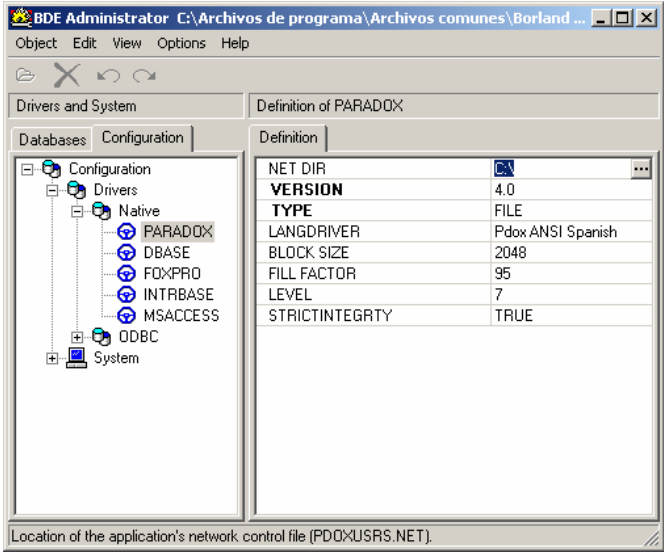
PROCEDURE

Process	Expected Results
1. Install the software in all the computers that will access to the working directory. 2. Share the instrument working directory of the desired instrument to share. This has to be done in MS Windows way. Open an 'explorer' (Start->Run... "explorer") session and browse to find where the working directory is located. 3. Over the shared directory press right button and select the 'Share...' option.  Picture 1	The showing image will appear.
4. Select the 'share this folder' and fill the share name with some use to remember name. Check the permissions if the network is insecure and select the users who will have write permissions and which not. Click OK. (image 01)	Share will be available.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000781	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	23-May-2011
Title:	SM DG57-DG517. Networking.			Page 4 of 6

Process	Expected Results
 Picture 2 5. To configure the local system of the computer whom will connect to that remote computer working directory it's necessary to change the working directory. First it's necessary to map a unit letter to the desired computer destination. Open an explorer browser and at 'My Network' find the shared directory's computer and over the shared directory click right button and select 'Map network unit'.  Picture 3 6. Now to change the working directory to the remote one, we have to run the Config application and select the previous remote mapped unit. Click OK. 7. To configure the computer sharing the working directory, access the Regedit program (Start button → Run → Type "Regedit") 8. Unfold the following path: HKEY_LOCAL_MACHINE → SOFTWARE → GRIFOLS → SistemaDiana → General. 9. Double click on "DirBase" parameter.	Window as Picture 3 will appear. It will change the working directory. Regedit program opens. Options in selected path appear. A new window

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000781	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	23-May-2011
Title:	SM DG57-DG517. Networking.			Page 5 of 6

Process	Expected Results
10. Type the network path to the shared directory (\\server\share) in the example above. Click OK and close Regedit. 11. Run the application named 'bdeadmin.exe' (located at the windows program files, by default in MS Windows in English language will be named 'Program Files', and 'Common files\Borland Shared\bde'). Follow the path to 'Configuration->Drivers->Native->Paradox' and the property named 'NET DIR' should be changed by the name of the unit mapped to the remote share. See Picture 4.  The screenshot shows the 'BDE Administrator' window. On the left, a tree view shows the hierarchy: Configuration > Drivers > Native > PARADOX. The 'PARADOX' driver is selected. On the right, the 'Definition of PARADOX' tab is active, showing a list of properties: NET DIR, VERSION (4.0), TYPE (FILE), LANGDRIVER (Pdox ANSI Spanish), BLOCK SIZE (2048), FILL FACTOR (95), LEVEL (7), and STRICTINTEGRITY (TRUE). The 'NET DIR' property is highlighted with a mouse cursor. Picture 4 In the computer sharing the working directory, network path (\\server\share) must be written instead of a mapped directory.	appears allowing parameter edition Working directory is configured It will change the database searching directory.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000781	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	23-May-2011
Title:	SM DG57-DG517. Networking.			Page 6 of 6

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	Document update.	Differentiate local and remote computer configuration.
3.0	2005/05/09	<i>Document REGD-0000781_v2 replaces document 57008V0M000.</i>	Document reference changed.
3.0	See document header.	Added previous version 'Data', 'change', and 'rational' inside document change history section.	Update 'chance history' section with previous document version modifications.
3.0	See document header.	<i>Document REGD-0000777_v3 replaces document REGD-0000777_v2.</i>	Document version changed.
3.0	See document header.	Document title.	Commercial reference substituted by project reference.
3.0	See document header.	Header title.	Increase information to easy document search.