

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000795	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	07-Jun-2011
Title:	SM DG57-DG517. Home sensor adjustments.			Page 2 of 7

DG-57 / DG-517 Service Manual

HOME SENSORS ADJUSTMENT

Procedure description
This procedure describes how to adjust the mechanical position and electronic values of Samples and Reagents Home detectors (PCA Samples & Reagents Version2). <i>Note: Document REGD-0000795_v4 replaces document REGD-0000795_v3.</i>

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000795	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	07-Jun-2011
Title:	SM DG57-DG517. Home sensor adjustments.			Page 3 of 7

EQUIPMENT

1. Gauge
2. Multimeter.



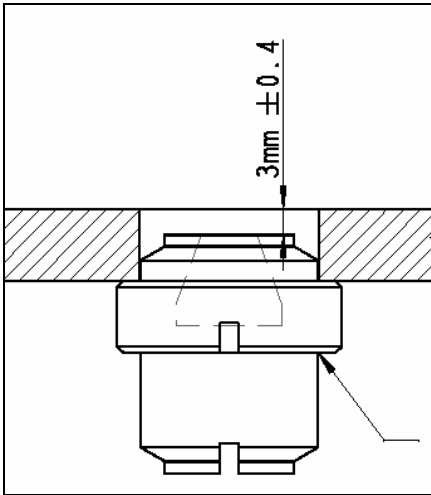
NOTE:

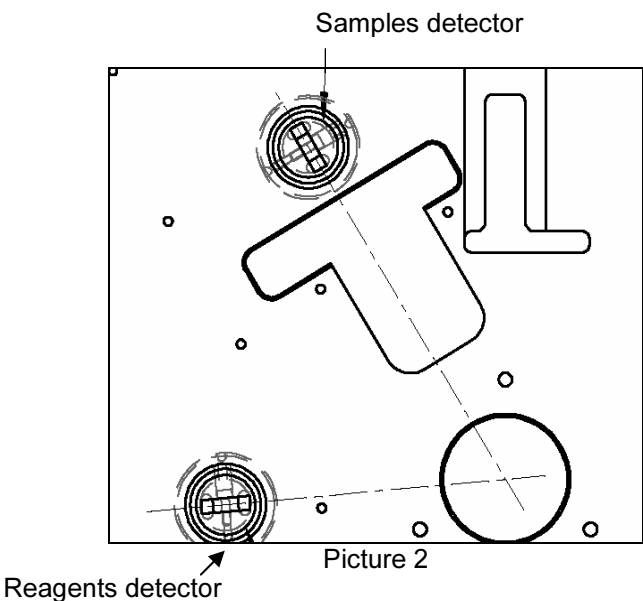
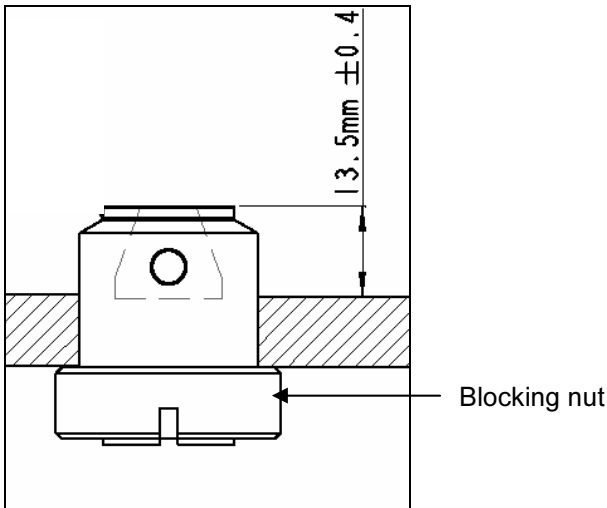
Ground in all measurements is in TP20. (see Annex PCA Layout)

SET UP

1. Equipment turned on.
2. Technical Service software started. See Document REGD-0000797.

PROCEDURE

Process	Expected Results
Adjusting Samples Home threshold. 1. Measure the voltage in TP12 with the multimeter. (see Annex PCA Layout) 2. If required, adjust P3 until the measured value reaches expected result. (see Annex PCA Layout)	$2,100 \pm 0,020 \text{ V}^*$ $2,100 \pm 0,020 \text{ V}^*$
Adjusting Reagents Home threshold. 3. Measure the voltage in TP11 with the multimeter. (see Annex PCA Layout) 4. If required, adjust P6 until the measured value reaches expected result. (see Annex PCA Layout)	$2,500 \pm 0,020 \text{ V}$ $2,500 \pm 0,020 \text{ V}$
Adjusting Samples Home position. 5. Measure the distance between the top of the detector and the top of the instrument basis.	---
 Picture 1	$3 \pm 0.4 \text{ mm}$
6. If required, adjust the height of the detector: • Loose blocking nut. Detector can be moved up/down by screwing/unscrewing to the basis. • Adjust detector height to expected result. • Detector orientation should be facing to carousel center (see Picture 2).	$3 \pm 0.4 \text{ mm}$

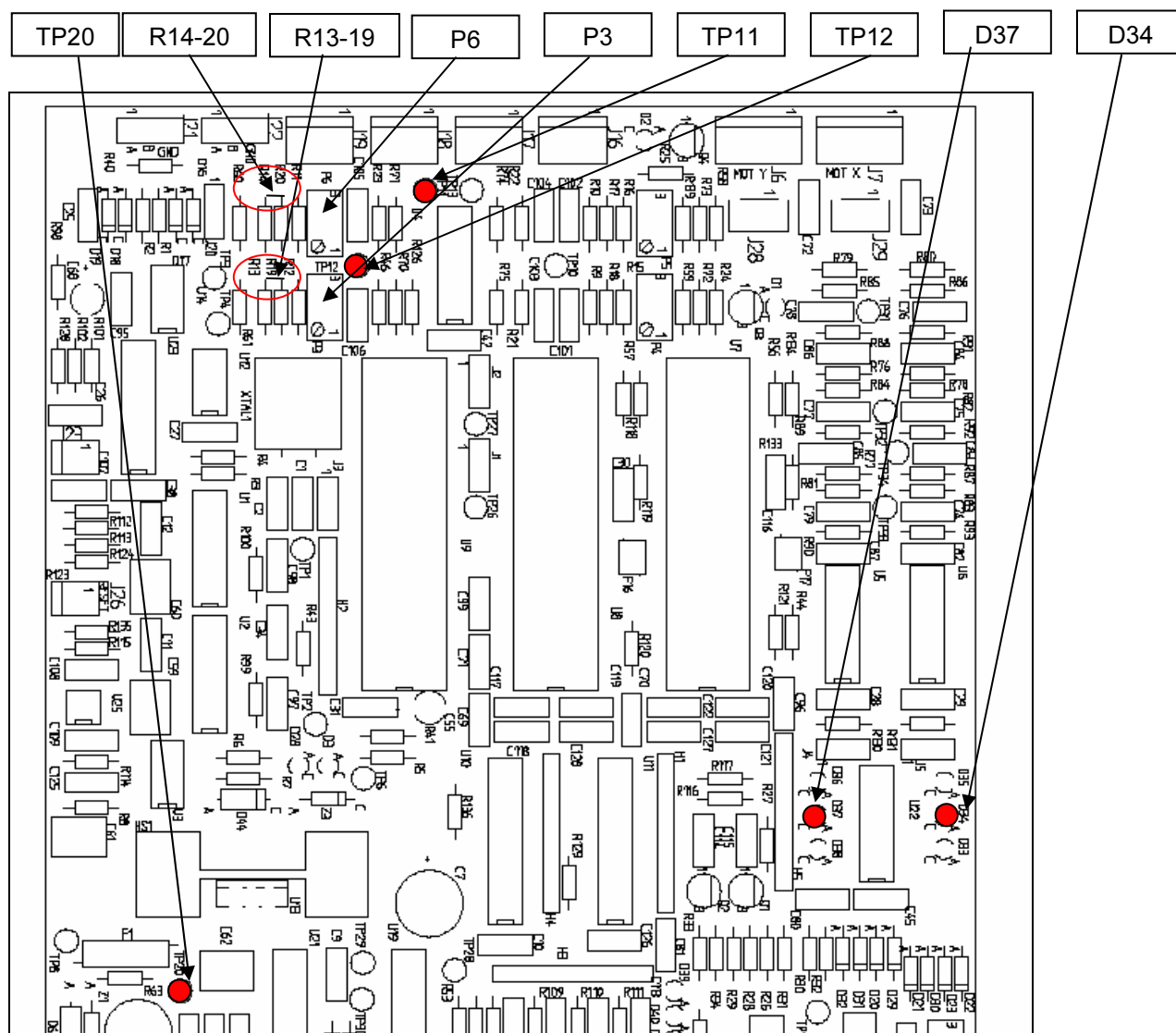
Process	Expected Results
- Tighten blocking nut.  Picture 2 Check the Reagents Home position. 7. Measure the distance between the top of the detector and the top of the instrument basis.  Picture 3 8. If required, adjust the height of the detector: - Loose blocking nut. Detector can be moved up/down by screwing/unscrewing to the basis. - Adjust detector height to expected result. - Detector orientation should be facing to carousel center. (See Picture 2) - Tighten blocking nut.	--- 13,5 ± 0.4 mm 13,5 ± 0.4 mm

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000795	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	07-Jun-2011
Title:	SM DG57-DG517. Home sensor adjustments.			Page 5 of 7

Process	Expected Results
Checking the Samples Home sensitivity. 9. Disable Samples and Reagents motors from Technical Service software: Press <i>Idle</i> button in Carousel tab. 10. Rotate manually the samples carousel until the led D37 is ON. 11. Measure the voltage in R13-19 (see Annex PCA Layout) and compare to the expected value. 12. Rotate manually the samples carousel until the led D37 is OFF. 13. Measure the voltage in R13-19 (see Annex PCA Layout) and compare to the expected value.	Motors are Disabled D37 ON < 1.3V D37 OFF > 3.8V
Checking the Reagents Home sensitivity. 14. Disable Samples and Reagents motors from the Technical Service software: Press <i>Idle</i> button in Carousel tab. 15. Rotate manually the samples carousel until the led D34 is ON 16. Measure the voltage in R14-20 (see Annex PCA Layout) and compare to the expected value. 17. Rotate manually the samples carousel until the led D34 is OFF 18. Measure the voltage in R14-20 (see Annex PCA Layout) and compare to the expected value.	Motors are Disabled D34 ON < 1.3V D34 OFF > 3.8V

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000795	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	07-Jun-2011
Title:	SM DG57-DG517. Home sensor adjustments.			Page 6 of 7

ANNEX PCA LAYOUT EST-CAR-2



GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000795	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	07-Jun-2011
Title:	SM DG57-DG517. Home sensor adjustments.			Page 7 of 7

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
3.0	2005/05/09	Document REGD-0000795_v3 replaces document 57023V0M001.	Document reference changed.
3.0	2005/05/09	- Added reference to PCA S&R Version 2 in Procedure Description. - included note for machine s/n samples home threshold dependence. - changed Samples and Reagents home LED's names.	Document updated.
4.0	See document header.	Added previous version 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
4.0	See document header.	<i>Document REGD-0000795_v4 replaces document REGD-0000795_v3.</i>	Document version changed.
4.0	See document header.	Document title.	Commercial reference substituted by project reference.
4.0	See document header.	Removed notes related to instruments with s/n< 057-114-xxxx to avoid serial numbers confusion.	Document available for DG-57 and DG-517.
4.0	See document header.	Document header.	Increase information to easy document search.