

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000796	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	07-Jun-2011
Title:	SM DG57-DG517. Reagents pipetting position offset adjustment.			Page 2 of 4

DG-57 / DG-517 Service Manual

REAGENT PIPETTING POSITION OFFSET ADJUSTMENT

Procedure description

This procedure describes how to adjust the reagent pipetting position offset, ensure that the probe is always perfectly centered to the vial neck.

Note: Document REGD-0000796_v3 replaces document REGD-0000796_v2.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000796	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	07-Jun-2011
Title:	SM DG57-DG517. Reagents pipetting position offset adjustment.			Page 3 of 4

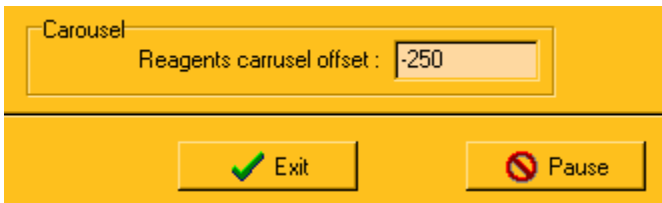
EQUIPMENT

1. An instrument connected with a PC with the software installed.

SET UP

1. Turn on the instrument.
2. Execute the ST software (Technical Service software).
3. Enter the Technical Service password (d+g_diana).

PROCEDURE

Process	Expected Results
Checking probe position in Reagent vial neck 1. Place an empty 10 ml Reagent vial in the reagent position nr.R1 2. Access on "Procedures" tab. 3. Select "Pipetting Test". 4. Press "Reagent" button (aspirate) and select 0 µl to aspirate. 5. Make a visual check for the position of the probe. 6. Press the "close" button. Adjusting probe position in Reagent vial neck 7. If the probe is not centered to the vial neck click on 'Settings' tab. 8. At the bottom of the screen there is a text window with the Reagents carousel offset. (See Picture 1). The default value is -250.  Picture 1 9. Increase or decrease offset value in intervals of ± 50 counts. Negative values mean CCW offset and positive values mean CW offset. 10. Recheck probe position (steps 1 to 6). Repeat the whole procedure until the probe is centered to the vial neck 11. Readjust Reagents camera following procedure REGD-0000811	The probe will go to the bottom of the Reagent vial and it will stop at the bottom giving a level detection error. --- --- --- Probe centered to vial neck ---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000796	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	07-Jun-2011
Title:	SM DG57-DG517. Reagents pipetting position offset adjustment.			Page 4 of 4

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	<i>Document REGD-0000796_v2 replaces document 57024V0M000.</i>	Document reference changed.
3.0	See document header.	Added previous version 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
3.0	See document header.	<i>Document REGD-0000796_v3 replaces document REGD-0000796_v2.</i>	Document version changed.
3.0	See document header.	Document title.	Commercial reference substituted by project reference.
3.0	See document header.	Document header.	Increase information to easy document search.