

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000799	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	16-Jun-2011
Title:	SM DG57-DG517. Temperature adjustment and verification.			Page 2 of 4

DG-57 / DG-517 Service Manual

TEMPERATURE ADJUSTMENT AND VERIFICATION

Procedure description
This procedure describes how adjust the temperature in the block incubator module. <i>Note: Document REGD-0000799_v3 replaces document REGD-0000799_v2.</i>


GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000799	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	16-Jun-2011
Title:	SM DG57-DG517. Temperature adjustment and verification.			Page 3 of 4

EQUIPMENT

1. Digital thermometer Hibok-15 MiniDin 4V (230903).
2. Two Temperature Probe PT100 for cards- MiniDin 4V (230905).

SET UP

1. Instrument turned on.
2. Technical Service software started. See Document REGD-0000821.
3. Place first Card Probe in position 8 of incubator 1 (16 cards) and second Card Probe in position 21 of incubator 2 (8 cards) at least 20 minutes before the adjustment takes place.

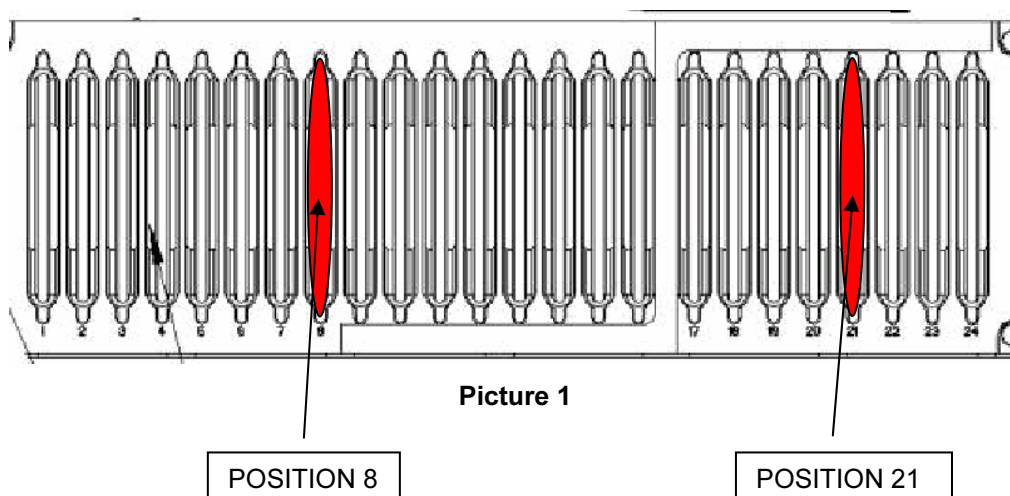

4. Room Temperature between 18°C and 28°C and instrument in user conditions, i.e. main door closed and upper cover mounted.

PROCEDURE

Process	Expected Results
Adjusting the incubator module:	
1. Press Start button in the box called Temperature of the incubator in Settings tab in Technical Service software. Follow the instructions given on the screen.	---
2. Step 1 of 5: Wait until temperature becomes stable (around 24°C).	---
3. Step 2 of 5: Follow the instructions given on the screen. Decimals are entered using a dot.	---
4. Step 3 of 5: Wait until temperature becomes stable (around 37°C).	---
5. Step 4 of 5: Follow the instructions on the screen. Decimals are entered using a dot.	---
6. Step 5 of 5: Accept the adjustment. The computer calculates values within the following limits. If not, carry out the adjustment again.	2200 to 2700 3600 to 3850
Checking the adjust of the incubator module:	
7. Press "On" button in Heating control option of I/O folder.	---
8. Select 24°C radio button in Temperature 1 and 2.	---
9. Wait until the temperature value stabilizes.	---
10. Check that the temperature of the card falls within the expected results.	24±3°C
11. Select 37°C radio button in Temperature 1 and 2.	---
12. Wait until the temperature value stabilizes.	---
13. Check that the temperature of the card falls within the expected results.	37±1°C

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000799	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	16-Jun-2011
Title:	SM DG57-DG517. Temperature adjustment and verification.			Page 4 of 4

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	<i>Document REGD-0000799_v2 replaces document 57027V0M000.</i>	Document reference changed.
3.0	See document header.	Added previous version 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
3.0	See document header.	<i>Document REGD-0000799_v3 replaces document REGD-0000799_v2.</i>	Document version changed.
3.0	See document header.	Required equipment updated.	Thermometer and temperature probe have changed.
3.0	See document header.	Updated procedure description.	Updated radio button name in step 8.
3.0	See document header.	Document title.	Commercial reference substituted by project reference.
3.0	See document header.	Document header.	Increase information to easy document search.