

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000800	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	16-Jun-2011
Title:	SM DG57-DG517. PCA IO57-2 and PCA Incubator fuse status verification.			Page 2 of 7

DG-57 / DG-517 Service Manual

PCA IO57-2 AND PCA INCUBATOR FUSE STATUS VERIFICATION

Procedure description
This procedure describes how check the fuses in block incubator module in order that the equipment is ready to use. <i>Note: Document REGD-0000800_v3 replaces document REGD-0000800_v2.</i>

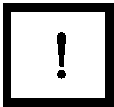
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000800	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	16-Jun-2011
Title:	SM DG57-DG517. PCA IO57-2 and PCA Incubator fuse status verification.			Page 3 of 7

EQUIPMENT

1. Multimeter.

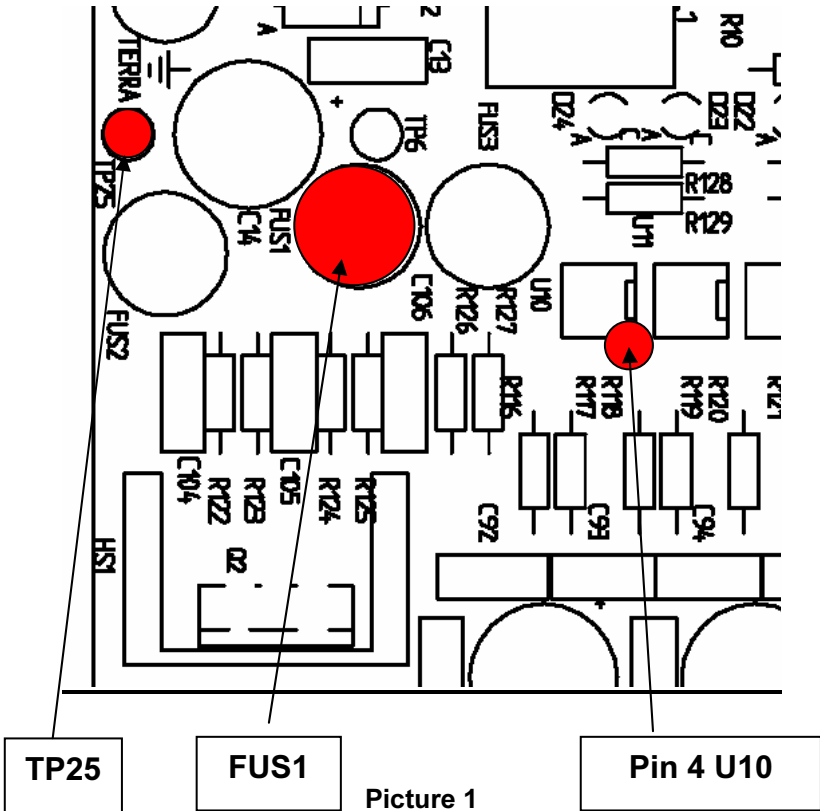
SET UP

1. Turn instrument off.
2. Open the rear panel
3. Turn instrument on.

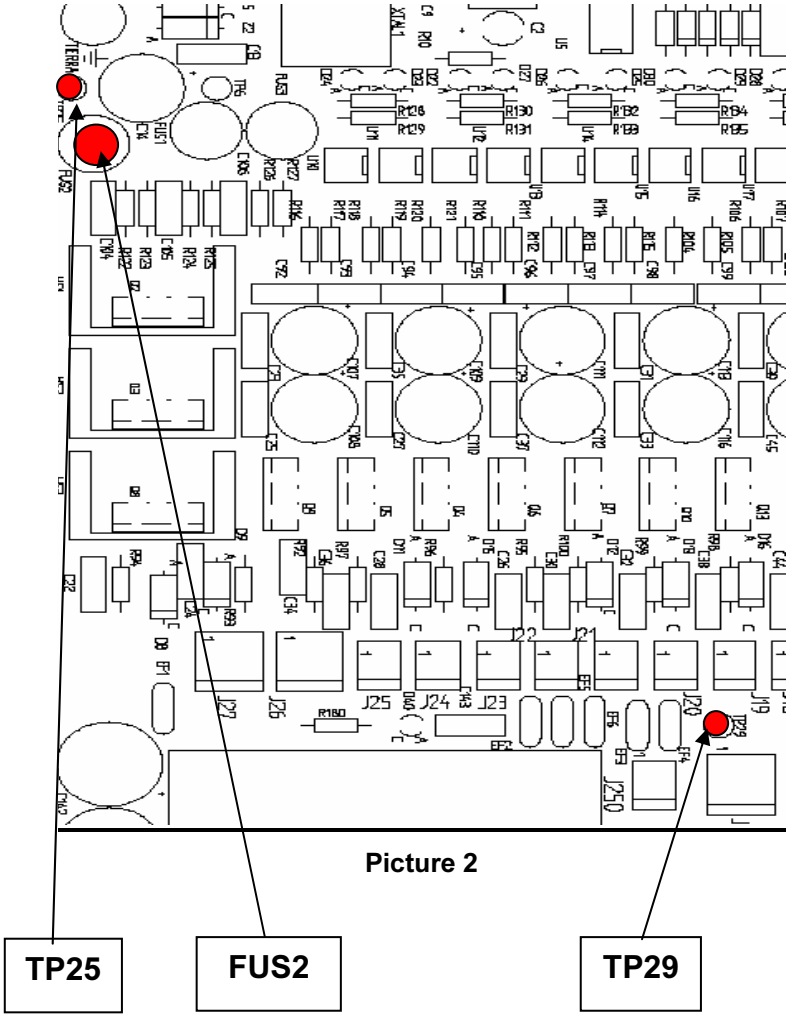


CAUTION!
Short-circuit between any powered part and ground, may damage the instrument.

PROCEDURE

Process	Expected Results
Checking the FUS1 fuse (Power incubator fuse) PCA IO57-2 1. To check the fuse status, measure the voltage between TP25 (negative) and Pin 4 of U10 (positive) with the multimeter (see Picture 1) 	23.95 ± 0,05 V ---
2. If the measured value is not the expected, then the FUS1 is broken.	---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000800	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	16-Jun-2011
Title:	SM DG57-DG517. PCA IO57-2 and PCA Incubator fuse status verification.			Page 4 of 7

Process	Expected Results
Checking the FUS2 fuse (Power Peltier fuse) PCA IO57-2 3. To check the fuse status, measure the voltage between TP25 (negative) and TP29 (positive) with the multimeter (see Picture 2)  Picture 2 TP25 FUS2 TP29	$12 \pm 0,1 \text{ V}$ ---
4. If the measured value is not the expected, then the FUS2 is broken.	

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000800	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	16-Jun-2011
Title:	SM DG57-DG517. PCA IO57-2 and PCA Incubator fuse status verification.			Page 6 of 7

Process	Expected Results
Checking the Thermal FUS fuse in incubator 1 (16 cards) 7. Turn off the instrument. 8. Disconnect the cable number 34. (see Picture 4) 34 35 Picture 4 9. Measure conductivity in cable connector. Compare the value with expected result. 10. If measured value is not the expected one the thermal fuse is blown. Checking the Thermal FUS fuse in incubator 2 (8 cards) 11. Turn off the instrument. 12. Disconnect the cable number 35 (see Picture 4). 13. Measure conductivity in cable connector. Compare the value with expected result. 14. If measured value is not the expected one the thermal fuse is blown.	-- --- SHORT CIRCUIT --- -- --- SHORT CIRCUIT ---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000800	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	16-Jun-2011
Title:	SM DG57-DG517. PCA IO57-2 and PCA Incubator fuse status verification.			Page 7 of 7

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	Document REGD-0000800_v2 replaces document 57028V0M000.	Document reference changed.
3.0	See document header.	Added previous version 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
3.0	See document header.	Document REGD-0000800_v3 replaces document REGD-0000800_v2.	Document version changed.
3.0	See document header.	Modified procedure description.	Points 9, 10, 13, 14 and its expected results are modified.
3.0	See document header.	Document title.	-Commercial reference substituted by project reference. -Title updated to contents.
3.0	See document header.	Document header.	Increase information to easy document search.