

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000803	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	08-Jul-2011
Title:	SM DG57. Unbalance detector adjustment and verification.			Page 2 of 4

DG-57 Service Manual

UNBALANCE DETECTOR ADJUSTMENT AND VERIFICATION

Procedure description
This procedure describes how to adjust the centrifuge unbalance detector in all instruments from serial number 57-98-xxx included and all those that have new silentblocks installed (see INR-0057#0033). <i>Note: Document REGD-0000803_v4 replaces document REGD-0000803_v3.</i>


GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000803	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	08-Jul-2011
Title:	SM DG57. Unbalance detector adjustment and verification.			Page 3 of 4

EQUIPMENT

1. An instrument connected with a PC with the software installed.
2. A Rack of Gel cards.

SET UP

1. Turn on the instrument.
2. Execute the ST software (Technical Service software).
3. Enter the Technical Service password (d+g_diana).
4. Loose the fastening screws, adjust the switch touching without pressure the centrifuge motor.

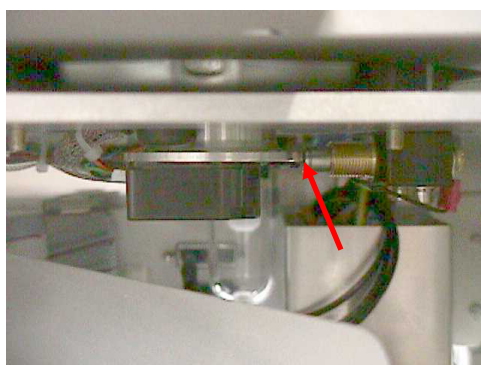


Picture 1

5. Tighten the fastening screws without moving the sensor.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000803	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	08-Jul-2011
Title:	SM DG57. Unbalance detector adjustment and verification.			Page 4 of 4

PROCEDURE.

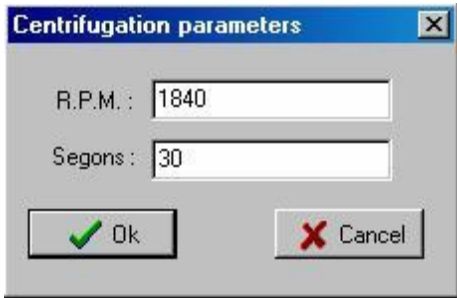Process	Expected Results
- Manually introduce 6 balanced cards on the centrifuge and close the front door. - Click on <i>Centrifuge</i> tab. Press <i>Reset</i> and then <i>Home</i> to initialize the centrifuge. Press <i>Centrifugation</i> button and select 1840 r.p.m and 30 seconds.  Picture 2	--- Switch not activated.
- If the CentErrCod N appears verify the position described in the "SETUP" section. If necessary, loose the fastening screws and relocate the sensor moving it away from the motor. Tighten the screws and repeat the step 2. - Manully take one card from the centrifuge and close the front door. - Repeat step 2.	--- --- CentErrCod N Appear.
If the switch is not activated (CentErrCod N doesn't appear) verify the position described in the "SETUP" section. If necessary, loose the fastening screws and relocate the sensor approaching it to the motor. Tighten the screws and repeat the procedure from the step 2.	---

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
3.0	2005/05/09	<i>Document REGD-0000803_v3 replaces document 57031V0M001.</i>	Document reference changed.
4.0	See document header.	Added previous version 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
4.0	See document header.	<i>Document REGD-0000803_v4 replaces document REGD-0000803_v3.</i>	Document version changed.
4.0	See document header.	Procedure updated.	Procedure updated only with 8 well cards.
4.0	See document header.	Document title.	Commercial reference substituted by project reference.
4.0	See document header.	Document header.	Increase information to easy document search.