

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000804	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	17-Jun-2011
Title:	SM DG57. Centrifuge module verification.			Page 2 of 6

DG-57 Service Manual

CENTRIFUGE MODULE VERIFICATION

Procedure description
This procedure describes how to check the centrifugal module parts: home detector, positioning performance, centrifuging performance and unbalancing detector. <i>Note: Document REGD-0000804_v3 replaces document REGD-0000804_v2.</i>

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000804	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	17-Jun-2011
Title:	SM DG57. Centrifuge module verification.			Page 3 of 6

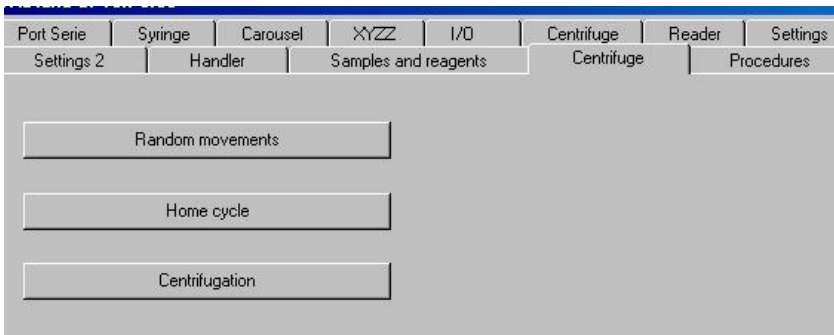
EQUIPMENT

1. An instrument connected with a PC with the software installed

SET UP

1. Turn on the instrument.
2. Execute the ST software (Technical Service software).
3. Enter the Technical Service password (d+g_diana).

PROCEDURE 1 - Checking home

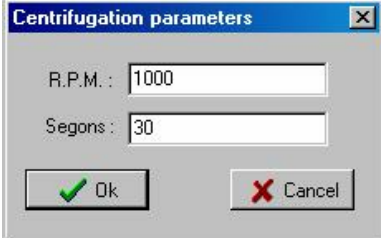
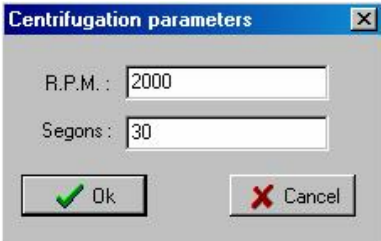
Process	Expected Results
1. Access high level 'Centrifuge' tab in ST software.  Picture 1 2. Press 'Home cycle' button. (see Picture 1)	The software looks like Picture 1. The centrifuge spins counterclockwise up to 1 turn many times until 'stop' button is pressed.

PROCEDURE 2 - Checking positioning

Process	Expected Results
1. Press 'Random movements' button. (see Picture 1)	The centrifuge moves randomly until 'stop' button is pressed.

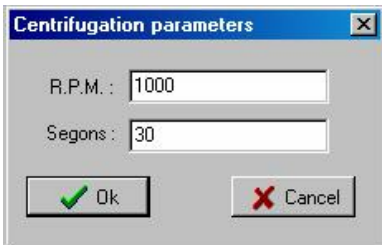
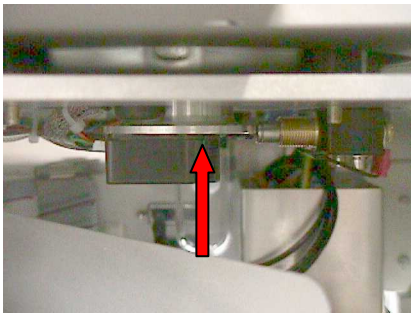
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000804	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	17-Jun-2011
Title:	SM DG57. Centrifuge module verification.			Page 4 of 6

PROCEDURE 3 - Checking centrifuging

Process	Expected Results
1. Press 'Centrifugation' button. (see Picture 1)  Picture 2	The Picture 2 window appears.
2. Press 'Ok' button. (see Picture 2)	The centrifuge accelerates up to 1000 rpms and maintains this speed during 30 seconds. When the centrifugation ends it starts a new centrifugation cycle until the stop button is pressed.
3. Press 'Stop' button.	Centrifugation stops after a centrifugation cycle.
4. Press 'Centrifugation' button. (see Picture 1)  Picture 3	The Picture 3 window appears.
5. Enter '2000' in the R.P.M. edit box and press 'Ok' button. (see Picture 3)	The centrifuge accelerates up to 2000 rpms and maintains this speed during 30 seconds. When the centrifugation ends it starts a new centrifugation cycle until the stop button is pressed.
6. Press 'Stop' button.	Centrifugation stops after a centrifugation cycle.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000804	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	17-Jun-2011
Title:	SM DG57. Centrifuge module verification.			Page 5 of 6

PROCEDURE 4 - Checking unbalancing detector

Process	Expected Results
1. Open the instruments door. 2. Press 'Centrifugation' button. (see Picture 1)  Picture 4	--- The Picture 4 window appears.
3. Locate the unbalancing detector just under the centrifuge.  Picture 5	See Picture 5.
4. Press 'Ok' button. (see Picture 4)	The centrifuge accelerates up to 1000 rpms.
5. Before 30 seconds, manually press unbalancing detector.	The centrifugation stops immediately. 'CentErrCod N' appears in the software,

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000804	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	17-Jun-2011
Title:	SM DG57. Centrifuge module verification.			Page 6 of 6

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	<i>Document REGD-0000804_v2 replaces document 57032V0M000.</i>	Document reference changed.
3.0	See document header.	Added previous version 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
3.0	See document header.	<i>Document REGD-0000804_v3 replaces document REGD-0000804_v2.</i>	Document version changed.
3.0	See document header.	Pictures updated.	Pictures updated with new software version.
3.0	See document header.	Document title.	Commercial reference substituted by project reference.
3.0	See document header.	Document header.	Increase information to easy document search.