

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000807	Version:	5.0, CURRENT
	Status:	Effective	Effective Date:	20-Jun-2011
Title:	SM DG57-DG517. Gripper position and probe height adjustment in incubator.			Page 2 of 7

DG-57 / DG-517 Service Manual

GRIPPER POSITION AND PROBE HEIGHT ADJUSTMENT IN INCUBATOR

Procedure description

This procedure describes how to adjust the gripper working position on the incubator and the probe height. With this procedure absolute coordinates value are stored for:

- gripper in position 1 of incubator rack (manually);
- gripper in positions 8, 16 and 24 of incubator rack (automatically);
- probe height on incubator;

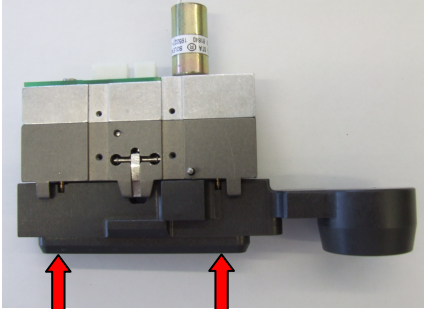
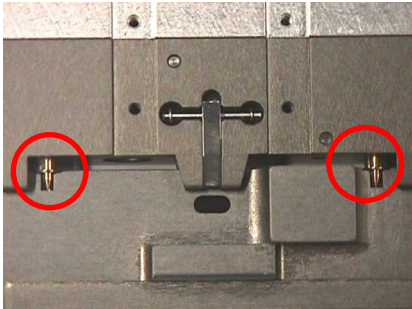
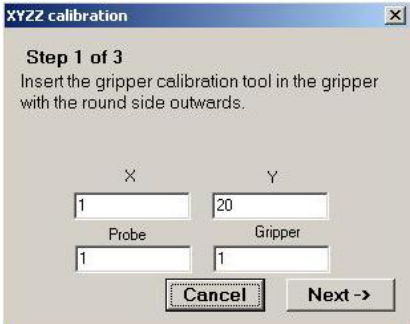
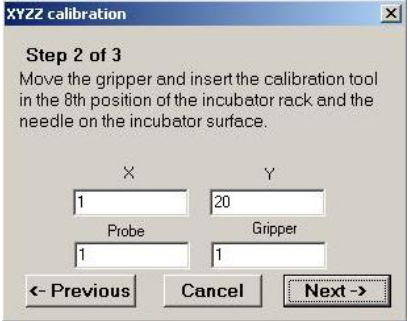
Any other positioning of probe and gripper is calculated relatively to this couple of values.

Disposing of the probe adjusting tool, a probe fine adjustment can be performed to improve probe positioning on main pipetting positions (see document *REGD-0001221*), but present adjustment procedure must be performed in advance.

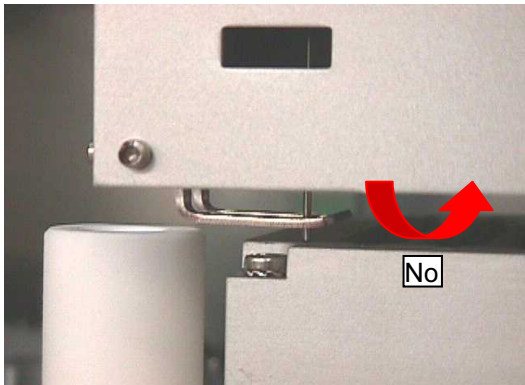
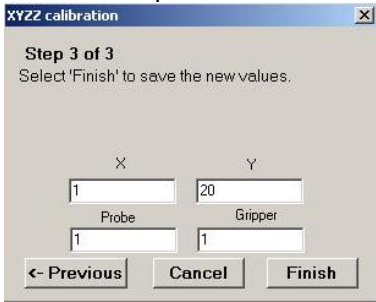
Once a probe X-Y fine adjustment has been performed (*REGD-0001221*), repeating the gripper adjusting procedure (*REGD-0000807*) does not affect the values stored with the former.

Note: Document REGD-0000807_v5 replaces document REGD-0000807_v4.

Reference:	REGD-0000807	Version:	5.0, CURRENT
Status:	Effective	Effective Date:	20-Jun-2011
Title:	SM DG57-DG517. Gripper position and probe height adjustment in incubator.	Page	4 of 7

Process	Expected Results
3. Follow the instructions on the screen. Insert the gripper calibration tool in the gripper with the round side outwards. To insert the gripper calibration tool is necessary to align the pins to it's placements in the tool and push the tool inside the gripper:    Picture 3 Picture 4 Picture 5 Picture 6 4. Click on "Next" button. 5. Move the gripper and insert the calibration tool in the 8th position of the incubator rack. Move the probe down by turning manually the pulley until it touches the incubator surface. Do not the force the probe down, beware to move up the instrument arm when the Probe touches the incubator surface.   Picture 7 Picture 8	The arm moves to the home position. New window appears (Picture 7) New window appears (Picture 10)

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000807	Version:	5.0, CURRENT
	Status:	Effective	Effective Date:	20-Jun-2011
Title:	SM DG57-DG517. Gripper position and probe height adjustment in incubator.			Page 5 of 7

Process	Expected Results
 Picture 9 6. Click on “Next” button. 7. Click the “Finish” button to finish the procedure.  Picture 10 8. The program automatically start an adjust procedure to calibrate the gripper on incubator positions 1 and 16 of block 1 and positions 17 and 24 of block 2. Procedure in correspondence of a single incubator position can be repeated several times, depending on the instrument. 9. At the end of previous procedure, remove the Gripper adjustment tool from the gripper. 10. Press the “Start” button from the “Gripper delay” screen.	New window appears (Picture 10) The arm moves to home position. The gripper is adjusted on incubator blocks edge positions New window appears (Picture 11)

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000807	Version:	5.0, CURRENT
	Status:	Effective	Effective Date:	20-Jun-2011
Title:	SM DG57-DG517. Gripper position and probe height adjustment in incubator.			Page 7 of 7

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
3.0	2005/07/26	Document REGD-0000807_v3 replaces document 57035V0M000.	Document reference changed.
3.0	2011/05/04	Updated procedure with gripper auto adjustment on incubator positions 1, 16, 17 and 24 and reference to new probe X-Y fine adjustment procedure (REGD-0001221).	Updated procedure.
4.0	2011/05/04	Added `previous version 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
4.0	2011/05/04	Updated the Document Change History: Date format changed to YYYY/MM/DD Dates changed to document version effective date.	Adoption of ISO 8601 date format to avoid ambiguity.
4.0	2011/05/04	ID-GELSTATION and PROVUE erased from the title.	Service manual intended for Grifols market.
4.0	2011/05/04	Document REGD-0000807_v4 replaces document REGD-0000807_v3.	Document version changed.
4.0	2011/05/04	Spare part updated.	Modified spare parts codes.
4.0	2011/05/04	Pictures updated.	Pictures updated with Diana System software version 3.2 and with gripper adjustment tool Ver.3..
5.0	See document header.	<i>Document REGD-0000807_v5 replaces document REGD-0000807_v4.</i>	Document version changed.
5.0	See document header.	Document title.	Commercial reference substituted by project reference.
5.0	See document header.	Document header.	Increase information to easy document search.