

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000808	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	20-Jun-2011
Title:	SM DG57-DG517. Gripper adjustment in centrifuge position.			Page 2 of 7

DG-57 / DG-517 Service Manual

GRIPPER ADJUSTMENT IN CENTRIFUGE POSITION

Procedure description
This procedure describes how to adjust the gripper working position in the centrifuge.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000808	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	20-Jun-2011
Title:	SM DG57-DG517. Gripper adjustment in centrifuge position.			Page 3 of 7



EQUIPMENT

1. An instrument connected with a PC with the software installed.
2. A Gripper adjustment tool (ver 3). (P/N: 232245).

SET UP


1. Turn on the instrument.
2. Start the ST software (Technical Service software).
3. Enter the right Technical Service password (d+g_diana).

PROCEDURE

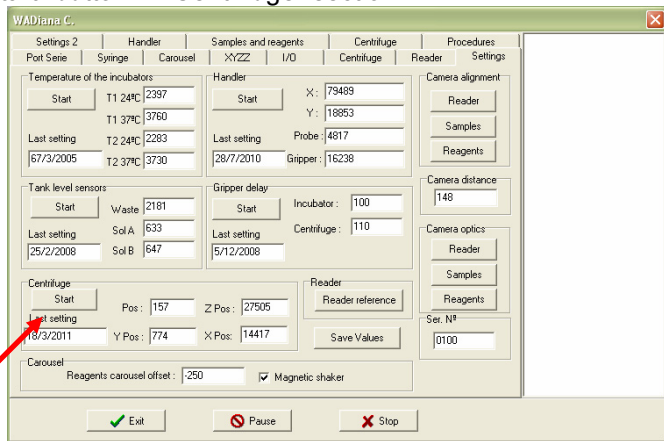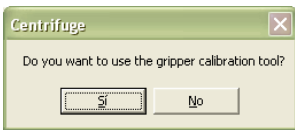
Process	Expected Results
1. Click on "Settings" tab. Click on "Start" button in "Centrifuge" section.  Picture 1	New window appears (Picture 2).
2. Click on "Yes" button.  Picture 2	New window appears (Picture 5).
3. Insert the gripper calibration tool in the gripper with the round side inwards. To insert the gripper calibration tool is necessary to align the pins to it's placements in the tool and push the tool inside the gripper:	The arm and the centrifuge rotor moves to home position.

Process


Expected Results

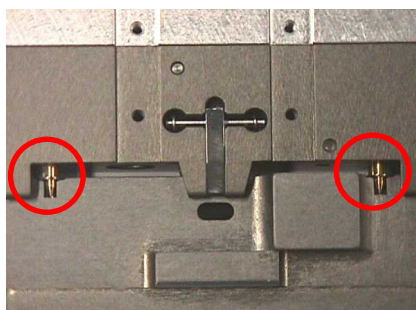

Picture 3

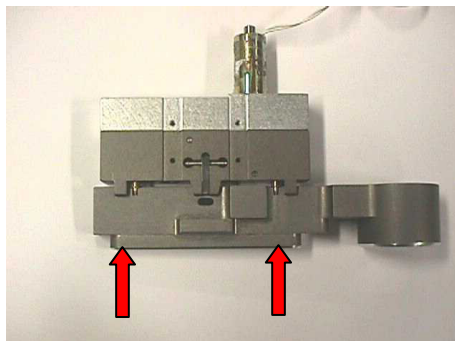

Picture 4

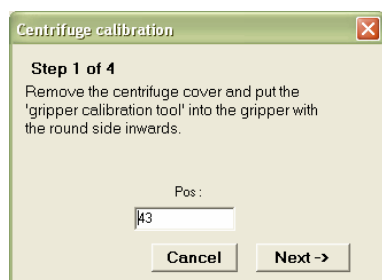


Picture 5



Picture 6

4. Click on "Next" button.
5. Remove all the centrifuge card holders. Place the instrument arm over the centrifuge and move down the gripper to fit the calibration tool in the centrifuge head.

New window appears (Picture 7).

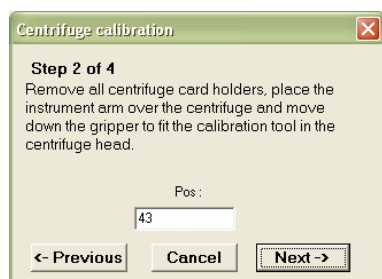


Picture 7



Picture 8


Process**Expected Results**

The tool have to fit in the centrifuge head like in Pictures 9 and 10:

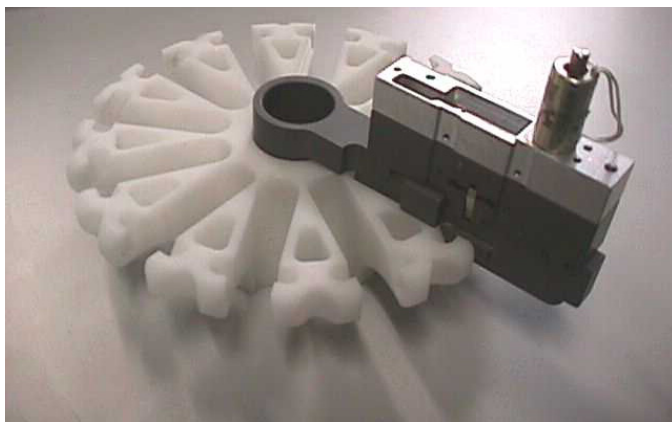

Picture 9


Beware special attention to fit the round side and the centering device:

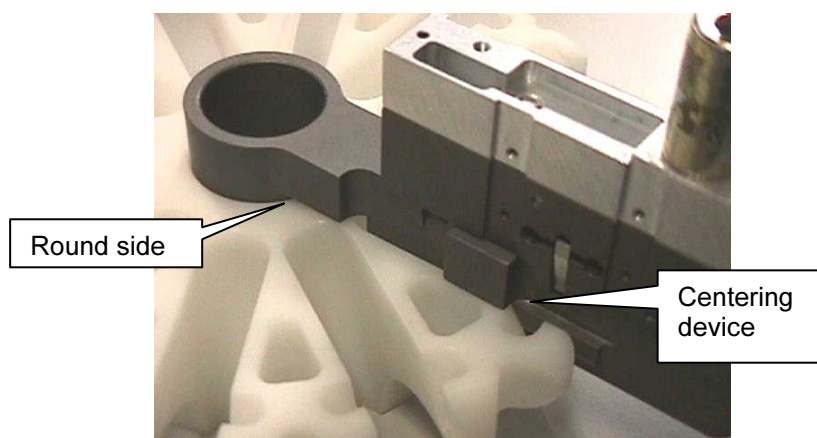

Picture 10

Don't force the instrument arm when fitting the tool inside the rotor.
Move the gripper up and down to fit the tool inside the rotor until the tool doesn't hit the rotor.

6. Click on "Next" button.
7. Wait until the calibration procedure ends.

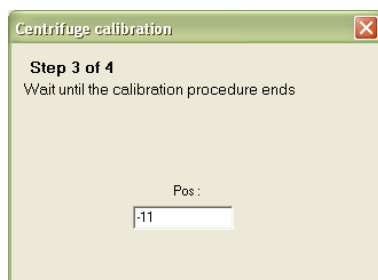

Picture 11

New window appears (Picture 11).

The instrument adjusts automatically all the positions of the centrifuge. New window appears (Picture 12).

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000808	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	20-Jun-2011
Title:	SM DG57-DG517. Gripper adjustment in centrifuge position.			Page 6 of 7

Process	Expected Results
8. Select “Finish” to save new values. Picture 12 Check that the values obtained during the previous adjustment fall within the recommended limits: Pos: 1-175 Y Pos: > 200 Z Pos: 26000-28000 X Pos: 12000-16000 9. Press “Start” button from the “Gripper delay” screen. Picture 13 Type ‘100’ in the Centrifuge delay window. Picture 14 10. Remove the tool from the gripper and replace the holders in the centrifuge and the cover of the centrifuge.	--- --- ---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000808	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	20-Jun-2011
Title:	SM DG57-DG517. Gripper adjustment in centrifuge position.			Page 7 of 7

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	Document REGD-0000808_v2 replaces document 57036V0M000.	Uploaded document to GxPharma.
3.0	2011/04/28	Pictures updated.	From versions v3.x, the Diana System installation program does not change the current Windows colour register.
3.0	2011/04/28	Changed the Gripper Adjustment Tool to version 3, code 232245.	There is a new Gripper Adjustment Tool.
3.0	2011/04/28	Added information about recommended values for the adjustment.	In order to check the correct adjustment it is included the values of the recommended limits.
3.0	2011/04/28	- Step 9: Updated the centrifuge value in the Gripper delay configuration from 100 to 110.	It is updated the gripper delay centrifuge to the current value.
3.0	2011/04/28	Updated the Document Change History: - Date format changed to YYYY/MM/DD - Dates changed to document version effective date.	Adoption of ISO 8601 date format to avoid ambiguity.
3.0	2011/04/28	ID-GELSTATION and PROVUE erased from the title.	Service manual intended for Grifols market.
4.0	See document header.	Document REGD-0000808_v4 replaces document REGD-0000808_v3.	Document version changed.
4.0	See document header.	Document title.	Commercial reference substituted by project reference.
4.0	See document header.	Document header.	Increase information to easy document search.