

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000811	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	28-Jun-2011
Title:	SM DG57-DG-517. Reagents camera adjustment.			Page 2 of 8

DG-57 / DG-517 Service Manual

REAGENTS CAMERA ADJUSTMENT

Procedure description
This procedure describes how to adjust Reagents camera.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000811	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	28-Jun-2011
Title:	SM DG57-DG-517. Reagents camera adjustment.			Page 3 of 8

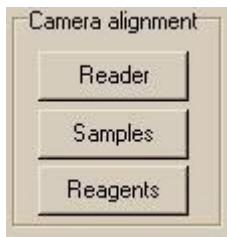
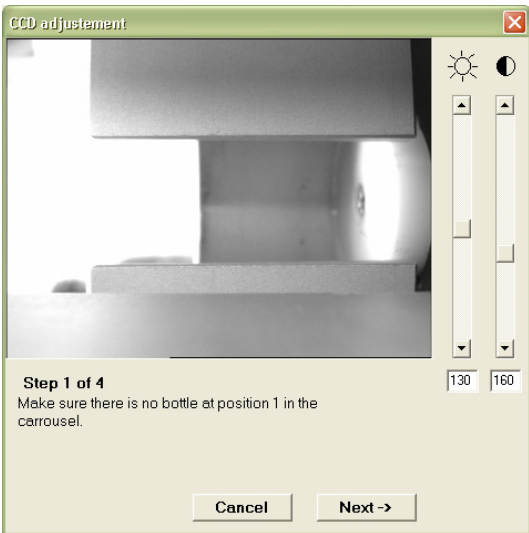
EQUIPMENT

1. An instrument connected with a PC with the software installed.
2. 5ml and 10ml vials.

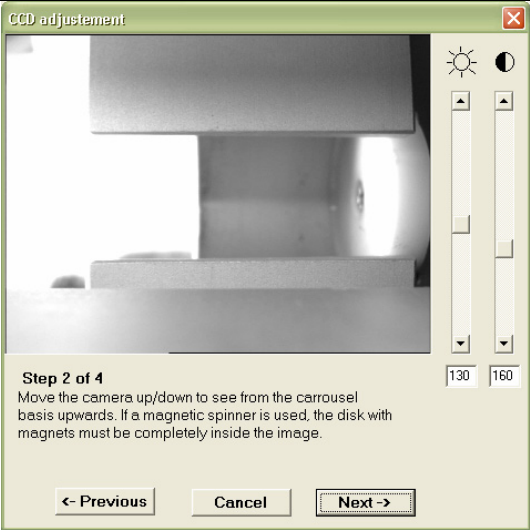
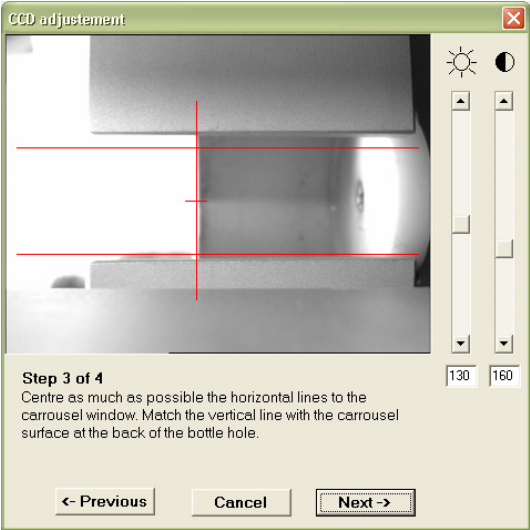
SET UP

1. Turn on the instrument.
2. Start the ST software (Technical Service software).
3. Enter the right Technical Service password (d+g_diana).

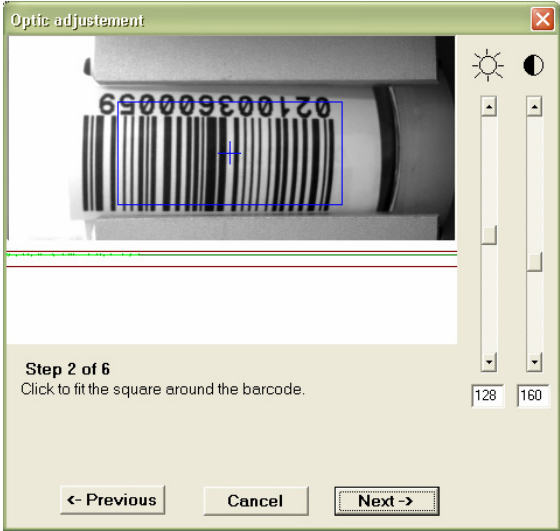
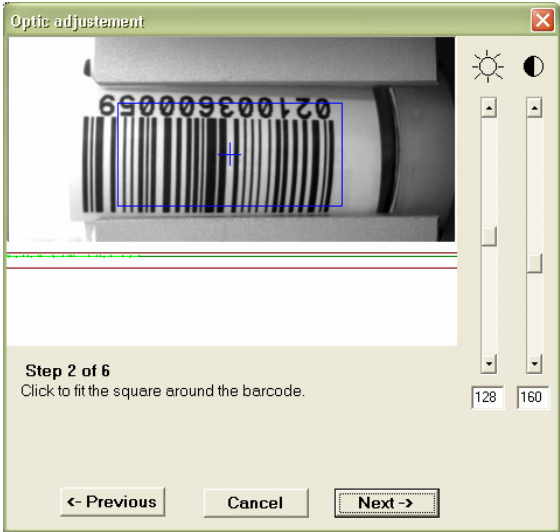
PROCEDURE

Process	Expected Results
Adjusting Reagents camera position 1. In <i>Settings</i> tab, click on <i>Reagents</i> button in <i>Camera alignment</i> section  Picture 1 2. Follow the instructions on the screen. Leave Reagent position 1 empty.  Picture 2 3. Click on <i>Next</i>. 4. Move the camera up/down, right/left to centre the position of the vial in the image. Take as a lower reference the white magnetic disk as shown in the image below. Lock the camera position.	New window appears (Picture 2). --- New window appears (Picture 3). ---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000811	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	28-Jun-2011
Title:	SM DG57-DG-517. Reagents camera adjustment.			Page 4 of 8

Process	Expected Results
 Picture 3 5. Click on <i>Next</i>. 6. Centre the lines to the position shown in the image. The vertical line must fit the top of the carousel at the vial's back wall. The two horizontal lines have to be centered inside the window.  Picture 4 7. Click on <i>Next</i>. A new window appears. Click on <i>Finish</i>. Adjusting Reagents camera optics 8. Disable temperature control: access <i>I/O</i> tab and click on <i>Off</i> checkbox in <i>Heating control</i> section. 9. In <i>Settings</i> tab, click on <i>Reagents</i> button in <i>Camera optics</i> section.	New window appears (Picture 4). --- --- --- New window appears (Picture 8).

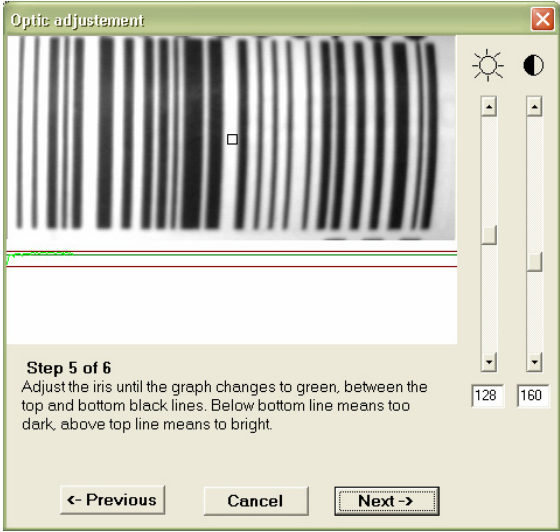
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000811	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	28-Jun-2011
Title:	SM DG57-DG-517. Reagents camera adjustment.			Page 5 of 8

Process	Expected Results
Camera optics Reader Samples Reagents Picture 7 10. Follow the instructions on the screen. Place a 10ml vial in position 1 of Reagents Carousel. Ensure the barcode of the label is facing outwards.  Picture 8 11. Click on <i>Next</i>. 12. Centre the blue square on the barcode. The blue cross has to be in the centre of the barcode.  Picture 9 13. Click on <i>Next</i>.	--- New window appears (Picture 9). --- New window

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000811	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	28-Jun-2011
Title:	SM DG57-DG-517. Reagents camera adjustment.			Page 6 of 8

Process	Expected Results
14. Open the iris until the image gets too bright. Move the focus until the graph gets to its maximum position. Depending on many factors the graph may get higher or lower, but it is possible to find only one peak. Picture 10 15. Click on <i>Next</i>. 16. Place the centre of the square to the brightest part of the label. Avoid the reflection coming from the LED's. Picture 11 17. Click on <i>Next</i>. 18. Move the iris until the graph turns to green. When the graph is red color the image is too bright (on top of upper horizontal line) or too dark (below the lower	appears (Picture 10). --- New window appears (Picture 11). --- New window appears (Picture 12). ---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000811	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	28-Jun-2011
Title:	SM DG57-DG-517. Reagents camera adjustment.			Page 7 of 8

Process	Expected Results
horizontal line)  Picture 12 19. Click on <i>Next</i>. New window appears. Wait until <i>Finish</i> button is available and then click it. 20. Lock the two camera rings. Checking Reagents camera adjustment 21. Place a 10ml vial in any position of the Reagents Carousel. 22. Click on <i>Read Reagent</i> button in <i>Reader</i> tab. 23. Type the desired position in the Carousel and click OK. 24. Click on <i>Read</i> button. The Reagent vial will be read. Ensure that the camera is able to read user's barcodes. 25. Repeat steps 21 to 24 with 5ml vials.	--- --- --- New window appears requesting position in Carousel. New window appears. --- ---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000811	Version:	4.0, CURRENT
	Status:	Effective	Effective Date:	28-Jun-2011
Title:	SM DG57-DG-517. Reagents camera adjustment.			Page 8 of 8

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/06	Document REGD-0000811_v2 replaces document 57039V0M000.	Uploaded document in GxPharma.
3.0	2010/09/21	Updated all the pictures of the document.	The document is updated due to the implementation of new functionalities.
3.0	2010/09/21	Document REGD-0000811_v3 replaces document REGD-0000811_v2.	Document version changed.
4.0	See document header	Updated the Document Change History: -Date format changed to YYYY/MM/DD. -- -Dates changed to document version effective date.	Adoption of ISO 8601 date format to avoid ambiguity.
4.0	See document header.	Added previous version 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
4.0	See document header.	Document REGD-0000811_4 replaces document REGD-0000811_v3.	Document version changed.
4.0	See document header.	Document title.	Commercial reference substituted by project reference.
4.0	See document header.	Document header.	Increase information to easy document search.