

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000814		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Reference image creation			

DG-57 Service Manual

REFERENCE IMAGE CREATION

Procedure description
This procedure describes how to create a reference image. The reference image is required by the imaging software to have a blank reference in order to detect the reactions inside the wells. The image must be taken from a brand new neutral (NaCl) card. Taking a reference image must be done: • After adjusting Reader Camera (corrective action). See related document REGD-0000812. • After the fine adjustment of the Reader camera during a Preventive Maintenance. • During the first installation of an instrument. See related document REGD-0000784. • After a new image acquisition board installation. <i>Note: Document REGD-0000814_v5 replaces document REGD-0000814_v4.</i>

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000814		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Reference image creation			

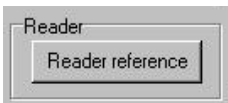
EQUIPMENT

1. An instrument connected with a PC with the software installed.
2. A brand new neutral card.

SET UP

1. Turn on the instrument.
2. Start the ST software (Technical Service software).
3. Enter the Technical Service password (d+g_diana).

PROCEDURE

Process	Expected Results
1. Ensure that the brand new neutral card that will be used is in good condition - It is not expired. - There is no gel/supernatant at the walls of the wells. Previous centrifugation is recommended. 2. Open the instrument door. 3. Insert the card in the first position of the incubator. 4. Access to 'Settings' tab in Technical Service software. 5. Click on 'Reader Reference' button.	--- --- --- --- The Gripper takes the card and reads it. When finished, the card is returned to its initial position.
 Picture 1 There is a time lag of about 15 seconds between reference reading and card's initial repositioning. 6. Open a Windows Explorer window. Access the folder where the software is installed. 7. Look for a file called W1xxxx4.tif, where 'xxxx' is the instrument serial number. 8. Double click on the file to open it. It should look similar to Picture 2.	--- ---
 Picture 2 If the card is in bad condition and there is not enough supernatant liquid in the wells, each well would look like Picture 3. In this case, find an appropriate brand new neutral card and start the procedure again.  Picture 3	---

Confidential Information. It is not allowed the total or partial reproduction of this document without prior authorization.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000814		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Reference image creation			

Process	Expected Results
Note: From Diana System version 4.1 onwards, along with the reference image, the picture taken by the camera is also stored as a TIFF image in the working directory under the name RefCard#, where # will increase from 0 to 9 as new reference images are taken. Once RefCard9 is created, the following images will overwrite the previous ones starting from RefCard0.	

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
3.0	2005/05/06	Document REGD-0000814_v3 replaces document 57041V0M001.	Document reference changed.
3.0	2005/05/06	- Added new reference image name for 8-well card. - Delay warning between reference reading and card positioning.	Document updated.
4.0	July 2011	Added previous version 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
4.0	July 2011	<i>Document REGD-0000814_v4 replaces document REGD-0000814_v3.</i>	Document version changed.
4.0	July 2011	picture 1 and 2 updated.	Pictures updated with last software version and with 8 well cards.
4.0	July 2011	Document title.	Commercial reference substituted by project reference.
4.0	July 2011	Document header.	Increase information to easy document search.
5.0	December 2014	Document header	Added information on the circumstances that require a reference image to be taken
5.0	December 2014	Procedure description	New functionality in Diana System v4.1