

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000820	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	08-Jul-2011
Title:	SM DG57-DG517. Clot detector sensibility adjustment.			Page 2 of 4

DG-57 / DG-517 Service Manual

CLOT DETECTOR SENSITIVITY ADJUSTMENT

Procedure description
This procedure describes how adjust the value of the potentiometer P1 of PCA Clot detector (Spare Part code 231829). <i>Note: Document REGD-0000820_v3 replaces document REGD-0000820_v2.</i>

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000820	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	08-Jul-2011
Title:	SM DG57-DG517. Clot detector sensibility adjustment.			Page 3 of 4

EQUIPMENT


1. Multimeter.



NOTE:
Ground in all measurements is TP2


SET UP

1. Instrument is turned off.
2. Open the upper panel.
3. Turn instrument on.

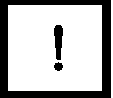

CAUTION!
Short-circuit between any powered part and ground, may damage the instrument.

PROCEDURE

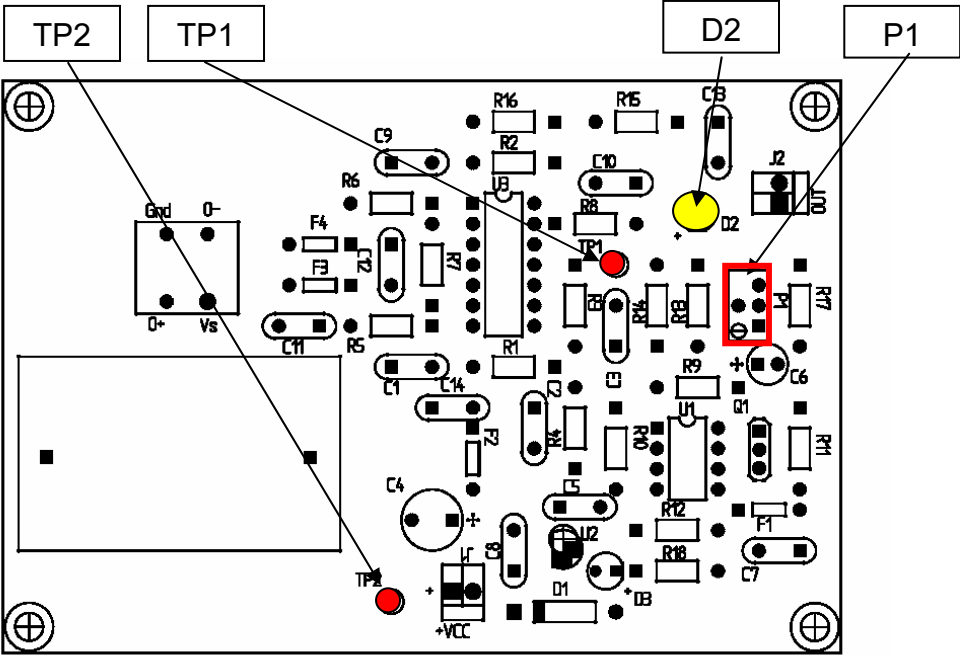Process	Expected Results
Adjusting the clot detector 1. Position the multimeter in TP1 (positive) and TP2 (negative). Adjust P1 until the measured value reaches expected result.  Picture 1 Checking the clot detection 2. Press Syringe file in Technical Service software. 3. Press Init Right and afterwards Execute button.	1,3±0,1 V --- --- Digital pump

