

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000834	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	19-Jul-2011
Title:	SM DG57-DG517. Incubator verification procedure.			Page 2 of 7

DG-57 / DG-517 Service Manual

INCUBATOR VERIFICATION PROCEDURE

Procedure description
This procedure describes how to carry out the verification of the temperature incubator in both temperatures: 24°C and 37°C <i>Note: Document REGD-0000777_v3 replaces document REGD-0000777_v2.</i>

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000834	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	19-Jul-2011
Title:	SM DG57-DG517. Incubator verification procedure.			Page 3 of 7

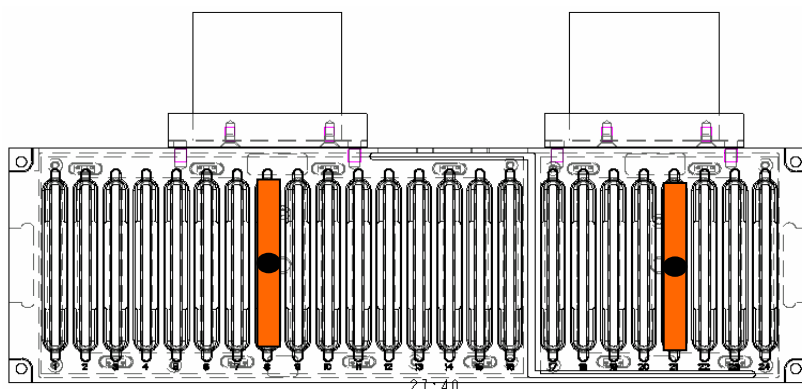
PROCEDURE 1 - RT TEMPERATURE VERIFICATION

EQUIPMENT

1. Card-Probe temperature. Spare part code:
2. Digital temperature Spare part code:

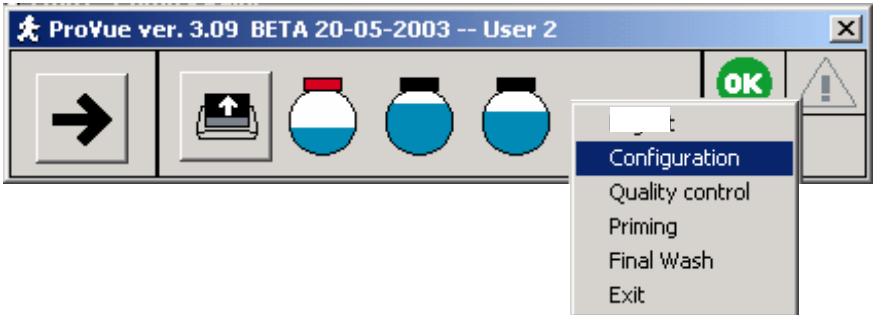
SET UP

1. Turn instrument on
2. Execute WADiana C. Software.
3. Wait until the instrument is initialized properly.
4. Put in the Card-Probe temperature at least 20 minutes before the verification takes places.
5. Place the Card-Probe in position 8 of block 1 and position 21 of block 2.



Before starting the procedure, be sure to close instrument's main door and mount upper cover, so as to make verification test under user conditions.

PROCEDURE

Process	Expected Results
1. Click on secondary mouse button in main windows. 2. Select "Configuration" option. See picture 1.  Picture 1	A pop-up menu appears Incubation process is performed
3. Enable "RT by default at 1st heating block" and "RT by default at 2on heating block". See Picture 2. 4. Click on "OK" button to accept.	--- RT regulation is

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000834	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	19-Jul-2011
Title:	SM DG57-DG517. Incubator verification procedure.			Page 4 of 7

Process	Expected Results
Configuration Working Tests ☒ Prerunning system tests ☐ Automatic disposal of cards ☒ RT by default at 1st heating block ☒ RT by default at 2nd heating block Reagents ☒ Presence verification ☐ BC reading ☐ Volumes verification ☐ Control of expiration dates OK Picture 2 5. Wait 20 minutes. 6. Check than the measure of the digital temperature is inside of Expected Results. 7. If the measure is out of expected result, take a new measure with a maximum time of 40 minutes. 8. If the measure is out of expected result even after 40 minutes, refer to Adjusting and checking the temperature.	activated. --- 21°C to 27°C --- ---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000834	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	19-Jul-2011
Title:	SM DG57-DG517. Incubator verification procedure.			Page 5 of 7

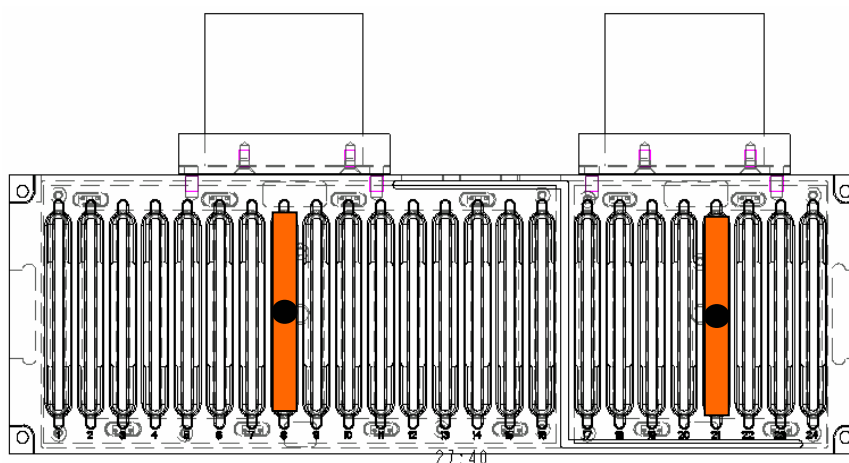
PROCEDURE 2 - 37°C TEMPERATURE VERIFICATION

EQUIPMENT

1. Card-Probe temperature. Spare part code:
2. Digital temperature Spare part code:

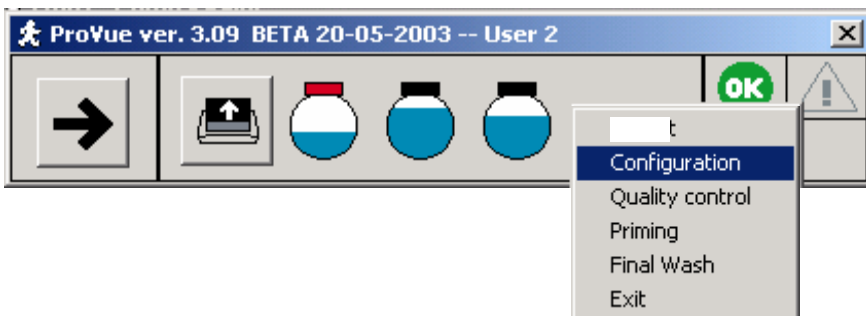
SET UP

1. **Turn** instrument **on**.
2. **Execute** WADiana C. Software.
3. **Wait** until the instrument is initialized properly.
4. **Put** in the Card-Probe temperature at least 20 minutes before the verification takes places.
5. **Place** the Card-Probe in position 8 of block 1 and position 21 of block 2.

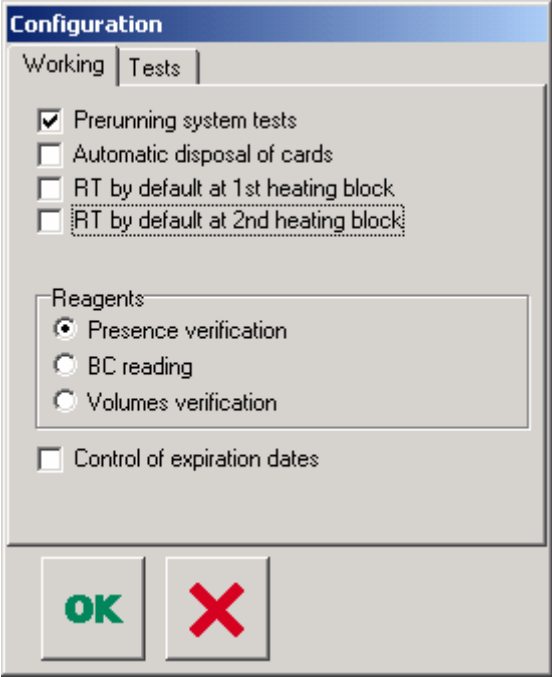


Before starting the procedure, be sure to close instrument's main door and mount upper cover, so as to make verification test under user conditions.

PROCEDURE

Process	Expected Results
1. Click on secondary mouse button in main windows. 2. Select "Configuration" option. See picture 1.  Picture 1. 3. Disable "RT by default at 1st heating block" and "RT by default at 2nd heating block". See Picture 2.	A pop-up menu appears. Incubation process is performed ---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000834	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	19-Jul-2011
Title:	SM DG57-DG517. Incubator verification procedure.			Page 6 of 7

Process	Expected Results
4. Click on “OK” button to accept.  Picture 2.	37°C regulating is activated.
5. Wait 20 minutes. 6. Check than the measure of the digital temperature is inside of Expected Results. 7. If the measure is out of expected result, take a new measure with a maximum time of 40 minutes. 8. If the measure is out of expected result even after 40 minutes, refer to Adjusting and checking the temperature.	--- 36°C to 38°C --- --

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000834	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	19-Jul-2011
Title:	SM DG57-DG517. Incubator verification procedure.			Page 7 of 7

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	Document REGD-0000834_v3 replaces document 57061V0M000.	Document reference changed.
2.0	2005/05/09	Procedure section.	Added user conditions reference in procedures setup section.
3.0	See document header	Added `previous version `data`, `change` and `rationale` inside document change history section.	Update change history for previous document history modification.
3.0	See document header	Document title.	-Commercial reference substituted by project reference. -`user` world deleted.
3.0	See document header	Document REGD-0000777_v3 replaces document REGD-0000777_v2.	Document version changed.
3.0	See document header	Document header.	-Increase information to easy document search. -`user` world deleted.