


GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000837	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	26-May-2011
Title:	SM DG57-DG517. Initialization troubleshooting.			Page 2 of 5

DG-57 / DG-517 Service Manual

INITIALIZATION TROUBLESHOOTING

Procedure description

This procedure describes how to find the cause the software can't initialize the instrument (See Picture 1).

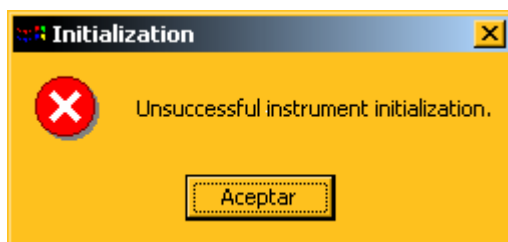

Picture 1

During the instrument initialization different checks are performed and all instrument settings are loaded in order to put the instrument and the software in a controlled state. If an error is detected during this process the general "Unsuccessful instrument initialization" message is shown. To identify the cause of the failure two main methods are used: audit trail file (.log) inspection and compare the performed movements with the expected ones. In order to this methods can be applied this document describes the log file and the expected movements.

Note: Document REGD-0000837_v3 replaces document REGD-0000837_v2.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000837	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	26-May-2011
Title:	SM DG57-DG517. Initialization troubleshooting.			Page 3 of 5


EQUIPMENT

1. An instrument connected with a PC with the software installed

SET UP

1. Turn on the instrument.
2. Execute the ST software (Technical Service software).

PROCEDURE 1 - Log file inspection

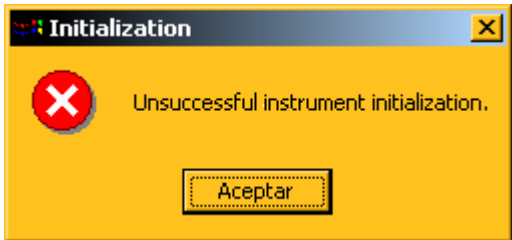Process	Expected Results
1. The message "Unsuccessful instrument initialization" appears, indicating an initialization failure.  Picture 1	The software looks like Picture 1.
2. Press the 'Accept' button. (see Picture 1) 3. Locate the directory of installation of the software and enter in the log subdirectory. 4. Open the Wadiana1.log file with the notepad. 5. Go to the end of the file. 6. Compare the log file with the sample in Annex A 7. The log file must be truncated at a given point with an ErrCod than caused the 'Unsuccessful instrument initialization'. Refer to document 'Error codes list' to get more information about that ErrCod. 8. Once you know which module have the problem you can use the corresponding Checking Procedure but instead of starting Technical Service software, start No Init Technical Service Software. Note than not all procedures can be executed with this option.	The software closes. --- --- --- --- --- ---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000837	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	26-May-2011
Title:	SM DG57-DG517. Initialization troubleshooting.			Page 4 of 5

Annex A- Log file sample.

```
[25/04/2003 11:04:47]-> *** Begin ProVue ver. X.XX ***
[25/04/2003 11:05:00]-> Initialitzing IO
[25/04/2003 11:05:00]-> Initialitzing Cattro plunger
[25/04/2003 11:05:00]-> Initialitzing XYZ
[25/04/2003 11:05:04]-> Initialitzing Centrifuge
[25/04/2003 11:05:04]-> Initialitzing Carrousel
[25/04/2003 11:05:04]-> Reading IO configurations
[25/04/2003 11:05:06]-> *** TAULA Calibration values ***
[25/04/2003 11:05:06]-> Empty Waste : 417
[25/04/2003 11:05:06]-> Full Waste : 2275
[25/04/2003 11:05:06]-> Empty A : 396
[25/04/2003 11:05:06]-> Full A : 2435
[25/04/2003 11:05:06]-> Empty B : 421
[25/04/2003 11:05:06]-> Full B : 2458
[25/04/2003 11:05:06]-> B1 24 : 2222
[25/04/2003 11:05:06]-> B1 37 : 3600
[25/04/2003 11:05:06]-> B2 24 : 2222
[25/04/2003 11:05:06]-> B2 37 : 3600
[25/04/2003 11:05:06]-> *** CARRU Calibration values ***
[25/04/2003 11:05:06]-> Reagents offset : -250
[25/04/2003 11:05:06]-> Reading Centrifuge configurations
[25/04/2003 11:05:08]-> *** CENTRI Calibration values ***
[25/04/2003 11:05:08]-> Base : 58
[25/04/2003 11:05:08]-> Holder 2 offset : 0
[25/04/2003 11:05:08]-> Holder 3 offset : -1
[25/04/2003 11:05:08]-> Holder 4 offset : 0
[25/04/2003 11:05:08]-> Holder 5 offset : 1
[25/04/2003 11:05:08]-> Holder 6 offset : 0
[25/04/2003 11:05:08]-> Holder 7 offset : 1
[25/04/2003 11:05:08]-> Holder 8 offset : 1
[25/04/2003 11:05:08]-> Holder 9 offset : 1
[25/04/2003 11:05:08]-> Holder 10 offset : 2
[25/04/2003 11:05:08]-> Holder 11 offset : 1
[25/04/2003 11:05:08]-> Holder 12 offset : 0
[25/04/2003 11:05:08]-> Reading XYZ configurations
[25/04/2003 11:05:10]-> *** BRAZO Calibration values ***
[25/04/2003 11:05:10]-> X Incubator: 79591
[25/04/2003 11:05:10]-> Y Incubator: 19235
[25/04/2003 11:05:10]-> Z Incubator: 16267
[25/04/2003 11:05:10]-> Z Needle : 1918
[25/04/2003 11:05:10]-> Z Reader : 79
[25/04/2003 11:05:10]-> X Centri : 14563
[25/04/2003 11:05:10]-> Y Centri : 1446
[25/04/2003 11:05:10]-> Z Centri : 26098
[25/04/2003 11:05:10]-> Z Delay Incub : 100
[25/04/2003 11:05:10]-> Z Delay Centri : 160
[25/04/2003 11:05:10]-> Magnetic shaker: Yes
[25/04/2003 11:05:10]-> Reading Vision configurations
[25/04/2003 11:05:13]-> *** VISIO Calibration values ***
[25/04/2003 11:05:13]-> Bright cam0 : 142
[25/04/2003 11:05:13]-> Contrast cam0 : 160
[25/04/2003 11:05:13]-> X cam0 : 382
[25/04/2003 11:05:13]-> Y cam0 : 280
[25/04/2003 11:05:13]-> Bright cam1 : 131
[25/04/2003 11:05:13]-> Contrast cam1 : 160
[25/04/2003 11:05:13]-> X cam1 : 340
[25/04/2003 11:05:13]-> Y cam1 : 286
[25/04/2003 11:05:13]-> Bright cam2 : 203
[25/04/2003 11:05:13]-> Contrast cam2 : 160
[25/04/2003 11:05:13]-> X cam2 : 404
[25/04/2003 11:05:13]-> Y cam2 : 244
[25/04/2003 11:05:13]-> Sample Camera distance : 148
[25/04/2003 11:05:13]-> Serial num. : 100
[25/04/2003 11:05:13]-> Instrument initialitized. Ser. N° 100
[25/04/2003 11:05:15]-> *** TAULA Calibration values ***
[25/04/2003 11:05:15]-> Empty Waste : 417
[25/04/2003 11:05:15]-> Full Waste : 2275
[25/04/2003 11:05:15]-> Empty A : 396
[25/04/2003 11:05:15]-> Full A : 2435
[25/04/2003 11:05:15]-> Empty B : 421
[25/04/2003 11:05:15]-> Full B : 2458
```

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000837	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	26-May-2011
Title:	SM DG57-DG517. Initialization troubleshooting.			Page 5 of 5

```

[25/04/2003 11:05:15] -> B1 24      : 2222
[25/04/2003 11:05:15] -> B1 37      : 3600
[25/04/2003 11:05:15] -> B2 24      : 2222
[25/04/2003 11:05:15] -> B2 37      : 3600
[25/04/2003 11:05:17] -> Door Closed

```

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	<i>Document REGD-0000837_v2 replaces document 57068V0M000.</i>	Document reference changed.
3.0	See document header	Added 'previous version 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
3.0	See document header	<i>Document REGD-0000837_v3 replaces document REGD-0000837_v2.</i>	Document version changed.
3.0	See document header	Removed 'Error codes list' document reference number.	Generalize document explanation to DG-57 and DG-517.
3.0	See document header	Removed software version from log example.	Document updated for current and following versions.
3.0	See document header	Document title.	Commercial reference substituted by project reference.
3.0	See document header	Header title.	Increase information to easy document search.