

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000843	Version:	6.0, CURRENT
	Status:	Effective	Effective Date:	13-Jul-2011
Title:	SM DG57-DG517. Peltier cells replacement.			Page 2 of 5

DG-57 / DG-517 Service Manual

PELTIER CELLS REPLACEMENT

Procedure description
This document describes the procedure to replace the Peltier cells. <i>Note: Document REGD-0000843_v6 replaces document REGD-0000843_v5.</i>

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000843	Version:	6.0, CURRENT
	Status:	Effective	Effective Date:	13-Jul-2011
Title:	SM DG57-DG517. Peltier cells replacement.			Page 3 of 5

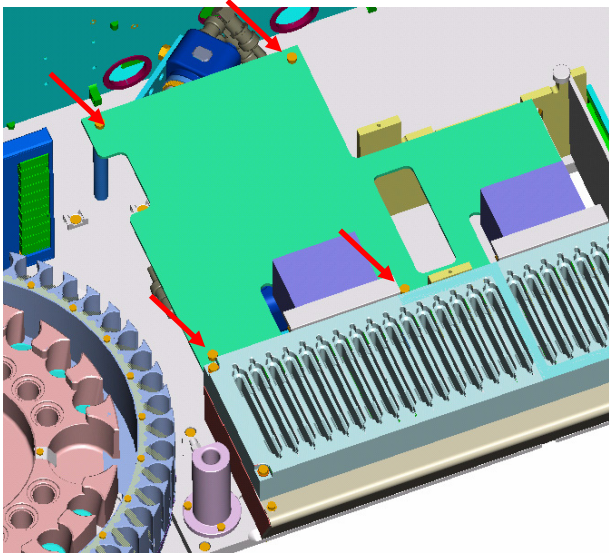
EQUIPMENT

1. 232081 Peltier module (with wire and plug).
2. 230902 Thermal silicon paste (90 ml).
3. Allen Keys

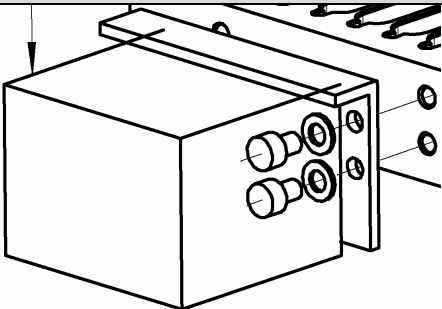
SET UP

1. Instrument switched off

PROCEDURE

Process	Expected Results
1. Remove the front cover incubator by unscrewing the two screws, one at each side.  Picture 1	---
2. Remove the reader cover by unscrewing 4 screws.  Picture 2	---
3. Each Peltier cell is held by 4 screws. Remove them from the cell to be replaced.	---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000843	Version:	6.0, CURRENT
	Status:	Effective	Effective Date:	13-Jul-2011
Title:	SM DG57-DG517. Peltier cells replacement.			Page 4 of 5

Process	Expected Results
 Picture 3 4. Disconnect the cable that connects the cell to the incubator. 5. Replace new peltier element (232081). 6. Place Thermal Silicon (230902) in the new peltier surface to be in contact with the incubator. Be aware not to leave areas without it, but also not to spread too much silicon. 7. Connect the connector to the incubator block. 8. Screw the four screws that hold the Peltier cell. Clean any rests of Thermal silicon that may overflow by its sides. 9. Place again the reader cover and the incubator front cover. 10. Use document REGD-0000799 to verify this replacement.	--- --- --- --- --- ---

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000843	Version:	6.0, CURRENT
	Status:	Effective	Effective Date:	13-Jul-2011
Title:	SM DG57-DG517. Peltier cells replacement.			Page 5 of 5

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
4.0	2005-05-09	Document REGD-0000843_v4 replaces document 57073V0M000.	Document reference changed
5.0	2011/04/28	Added previous version 'data', 'change' and 'rationale' inside Document change history section.	Update Change History for previous document version modification.
5.0	2011/04/28	Eliminated the reference to ID-GelStation and ProVue distributions.	Service manual intended for Grifols market'.
5.0	2011/04/28	Updated the Document Change History: - Date format changed to YYYY-MM-DD. - Dates changed to document version effective date.	Adoption of ISO 8601 date format to avoid ambiguity.
5.0	2011/04/28	Modified spare part names and added spare part code number inside 'equipment' section.	Document change to easy the way to locate part numbers.
5.0	2011/04/28	Inside procedure section: - One step added between previous processes steps 4 and 5. - Added spare part code numbers. - Modified previous step 5 description.	Document modified to easy procedure.
6.0	See document header.	<i>Document REGD-0000843_v6 replaces document REGD-0000843_v5.</i>	Document version changed.
6.0	See document header.	Document title.	Commercial reference substituted by project reference.
6.0	See document header.	Document header.	Increase information to easy document search.