

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000845		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Closed OS setup			

DG-57 Service Manual

CLOSED OS SETUP

Procedure description
This produce describes how to setup the Closed Operating System for Windows XP OS. This configuration does not apply to Windows 7 OS. Note: Document REGD-0000845_v5 replaces document REGD-0000845_v4.

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000845		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Closed OS setup			

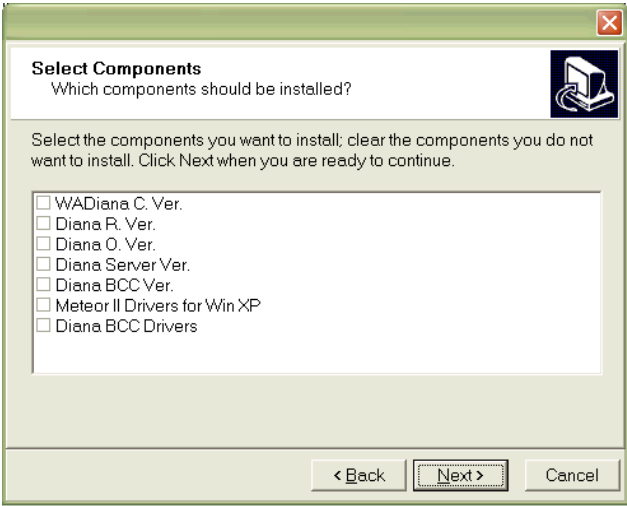
EQUIPMENT

1.

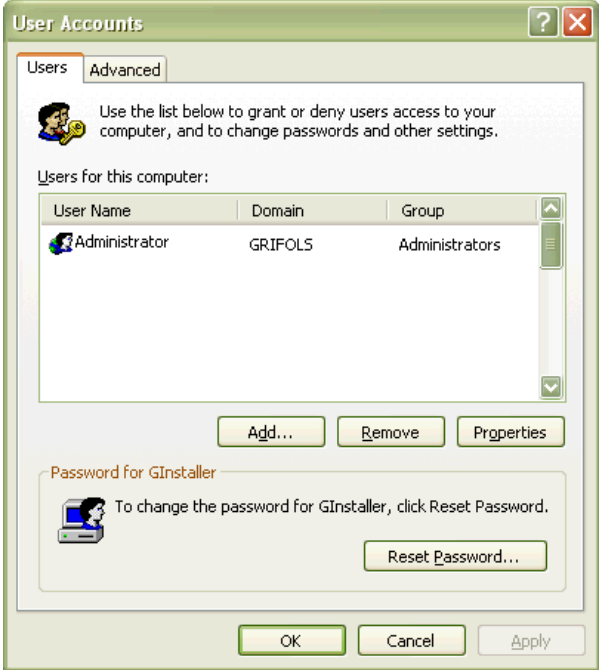
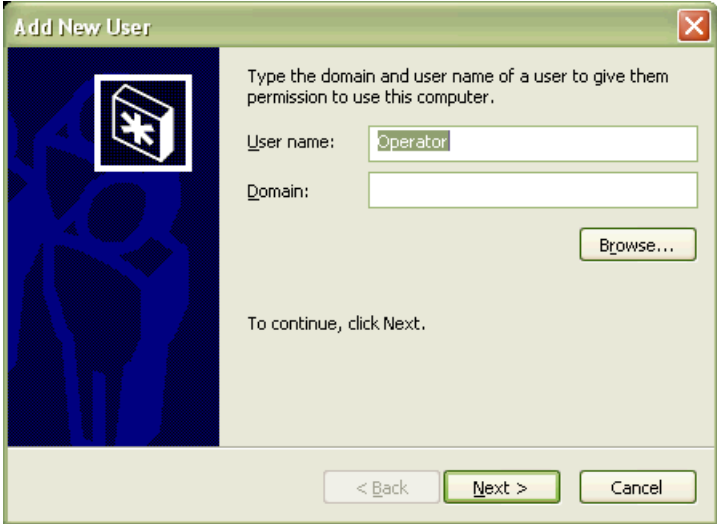
SET UP

1. PC turned on.

PROCEDURE

Process	Expected Results
1. Install Microsoft Windows XP with the latest Service Pack (3 or superior) at a NTFS partition. 2. Login to the system as a local Administrator. 3. Install the latest software in the predefined directory.	--- --- ---
 Picture 1	

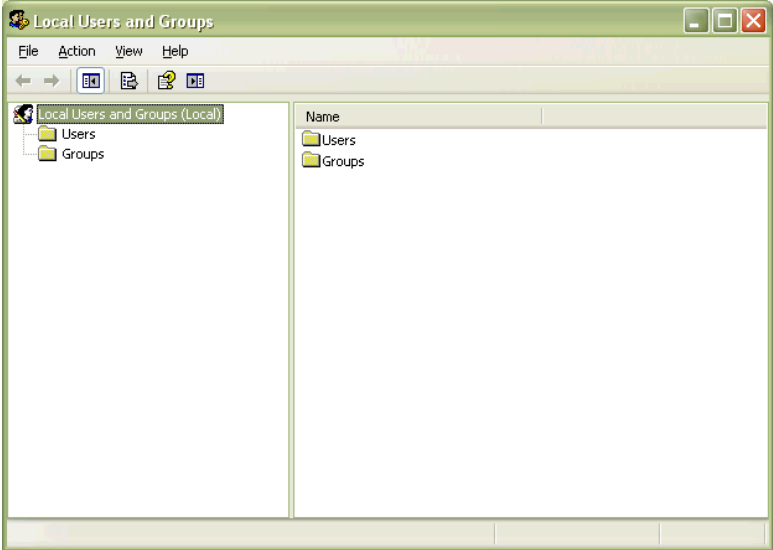
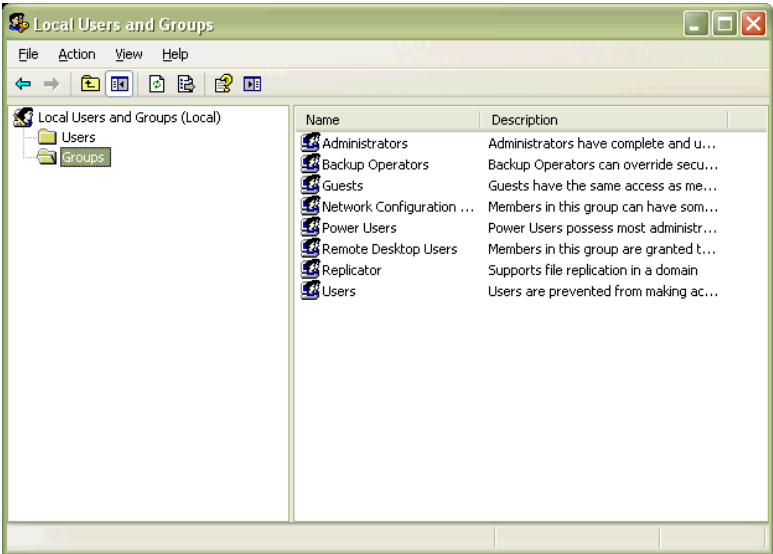
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000845		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Closed OS setup			

Process	Expected Results
4. Now to create the user or users to be in closed system click Start->->Control Panel ->User Accounts  Picture 2 5. Click on Add to create a new user like picture 3  Picture 3 6. Click Next and enter the password for this closed system user.	

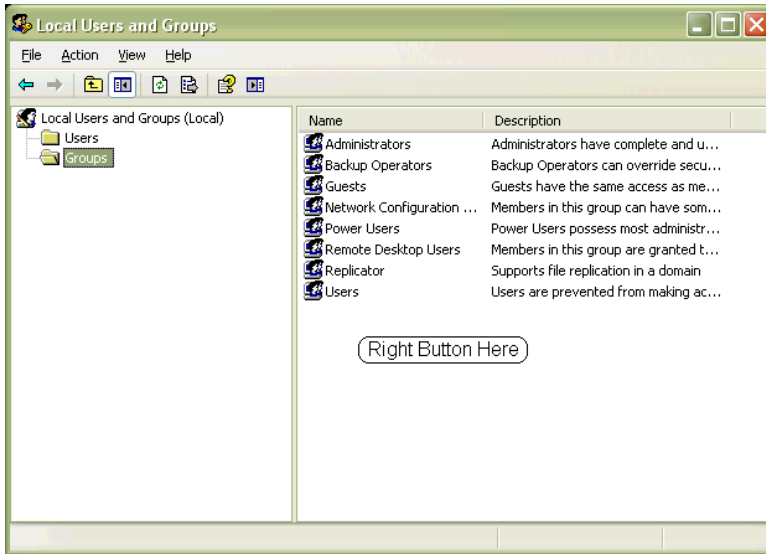
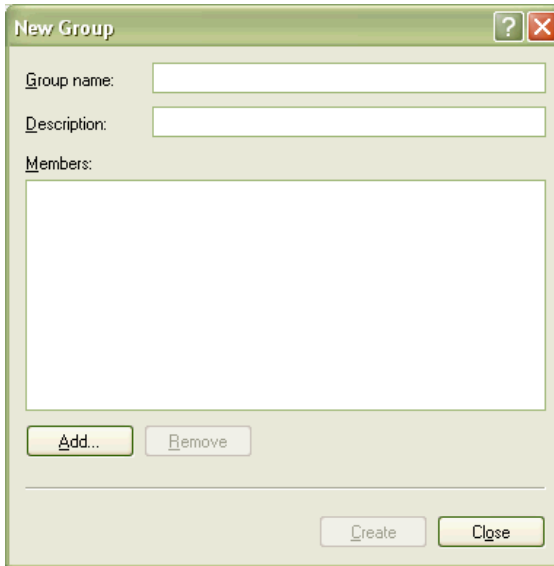
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000845		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Closed OS setup			

Process	Expected Results
7. Select type of user (level of access) other, select Administrators :  Picture 4 8. Now click Finish and click the Advanced button and over the Advanced tab at Advanced User Management part click Advanced button.  Picture 5	

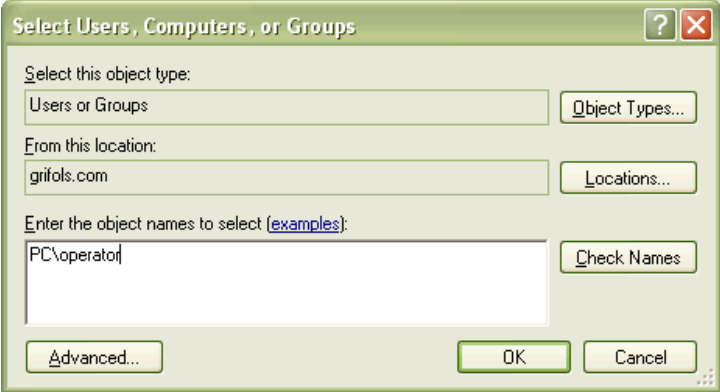
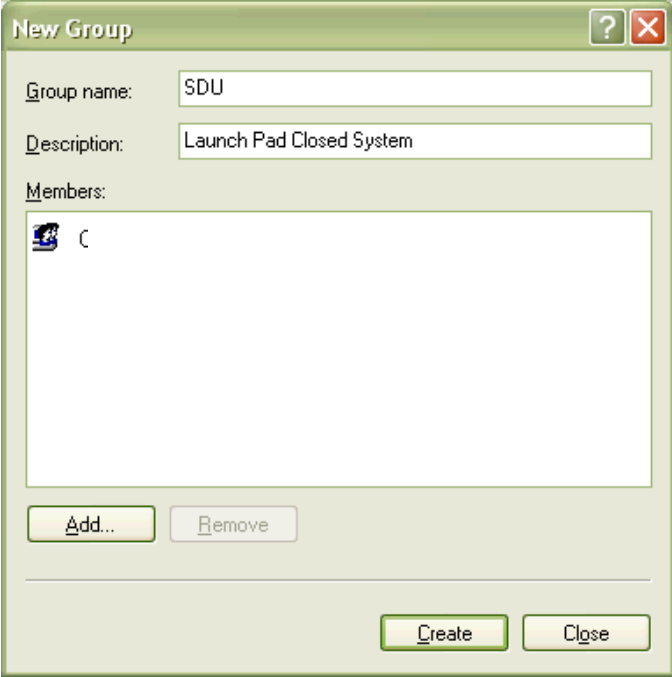
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000845		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Closed OS setup			

Process	Expected Results
9. Click over groups.  Picture 6  Picture 7	

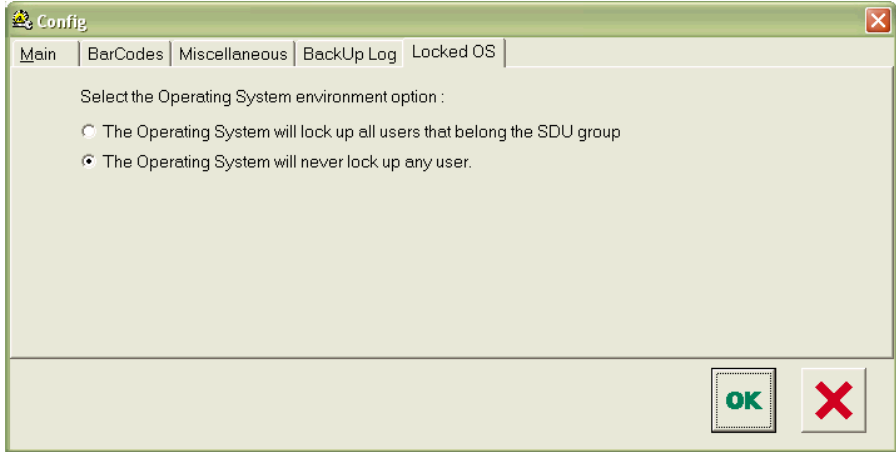
GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000845		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Closed OS setup			

Process	Expected Results
10. Right button over the white part of the application (where the image says).  Picture 8	
11. Select New Group and fill the Group Name with 'SDU' as image 9.  Picture 9	

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000845		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Closed OS setup			

Process	Expected Results
12. Click Add and select all users with Closed Operating System. All users belonging the SDU group will be OS closed.  Picture 10 13. Select the previous created user 'operator' (as example) and Add clicking on Add... and click OK button.  Picture 11 14. Click Create, Close, Close, OK.	

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000845		
	Version:	5.0; Most-Recent; Effective; CURRENT		
	Status:	Effective	Purpose:	Uncontrolled Copy
Class:	Control del Diseño		Effective date:	19-Dec-2014
Title:	SM DG57. Closed OS setup			

Process	Expected Results
15. Click Start->Programs->Sistema Diana->Config. Select the Locked OS tab and click on the radio button named 'The Operating System will lock up all users that belong the SDU group. Click OK button.  Picture12 16. Restart the computer.	

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	Document REGD-0000845_v2 replaces document 57075V0M000	Uploaded document to GxPharma.
3.0	2011/04/28	Updated procedure to Windows XP.	The current OS is the Windows XP.
3.0	2011/04/28	Updated the Document Change History: - Date format changed to YYYY/MM/DD - Dates changed to document version effective date.	Adoption of ISO 8601 date format to avoid ambiguity.
3.0	2011/04/28	ID-GELSTATION and PROVUE erased from the title.	Service manual intended for Grifols market.
4.0	July 2011	Document REGD-0000845_v4 replaces document REGD-0000845_v3.	Document version changed.
4.0	July 2011	Document title.	Commercial reference substituted by project reference.
4.0	July 2011	Document header.	Increase information to easy document search.
5.0	December 2014	Procedure description	Configuration is only available for Windows XP OS

Confidential Information. It is not allowed the total or partial reproduction of this document without prior authorization.

Print date:	19-Dec-2014	Page	9 of 9
--------------------	-------------	-------------	--------