


AmiCORE Parts Lookup Tool




Applicable for AmiCORE software: 2.1 (and later), China 2.0C (and later)

Special Tools

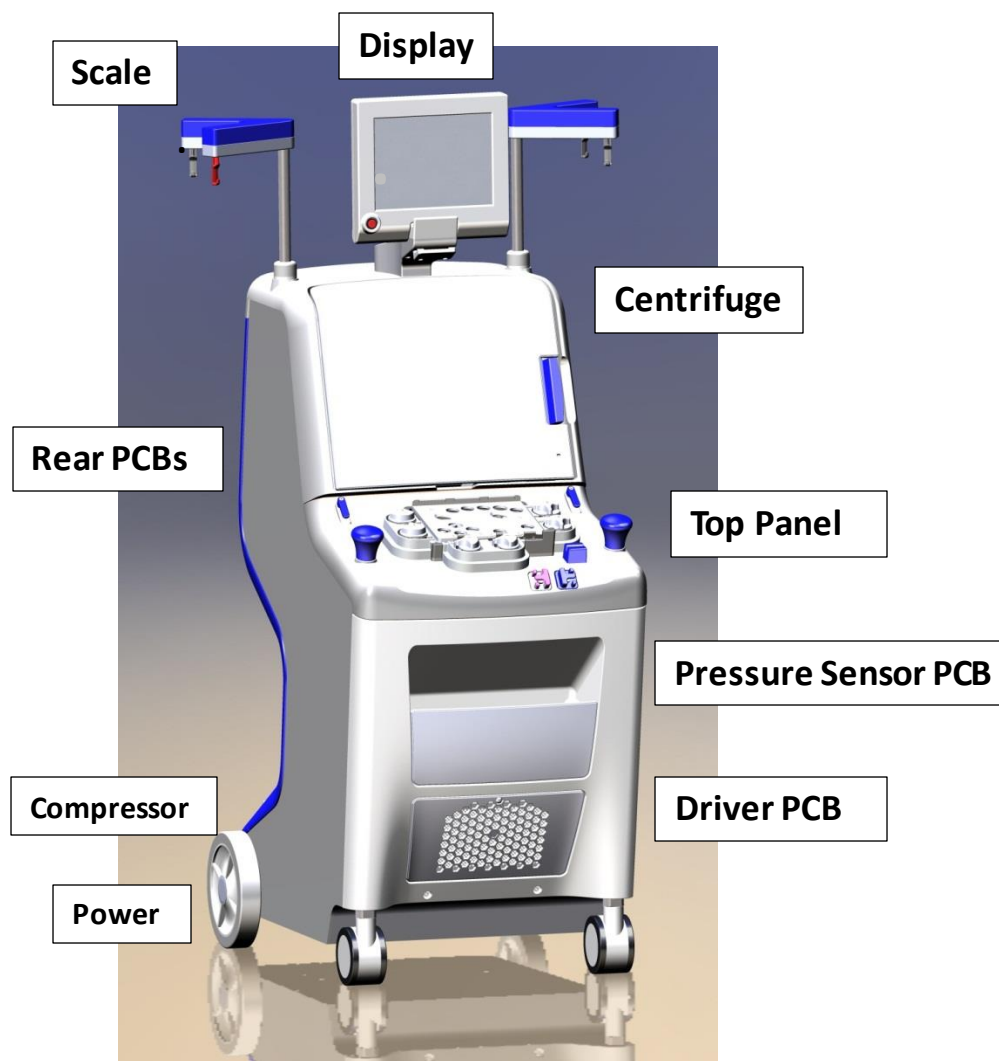


Rear PCB

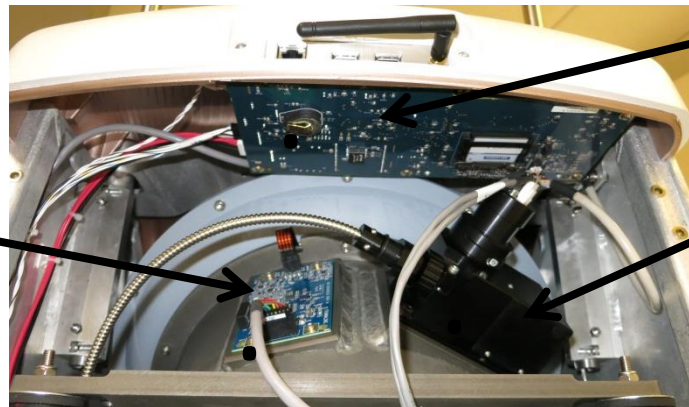

Display/
Comm

ID Assembly

Leak & Temp

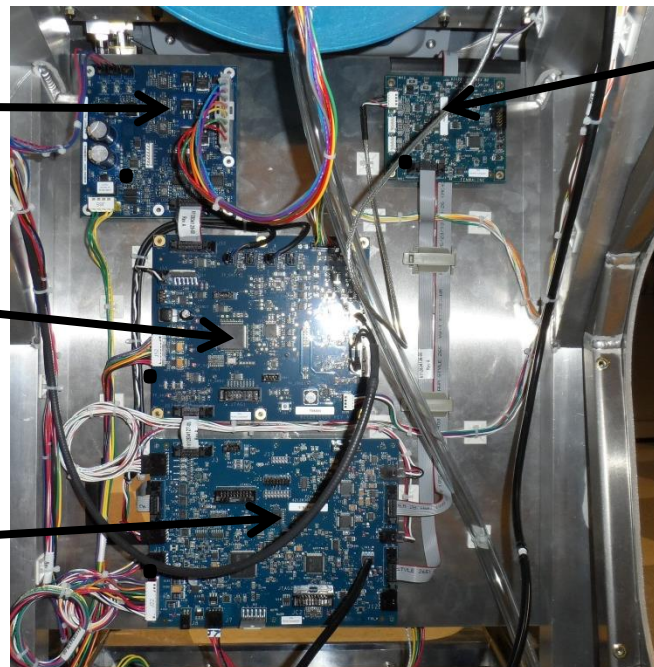


Weight

Centrifuge


Interface

Main

Power Distribution
Board

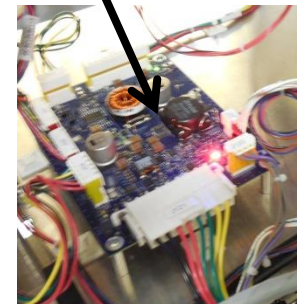


Driver PCB (located in front)

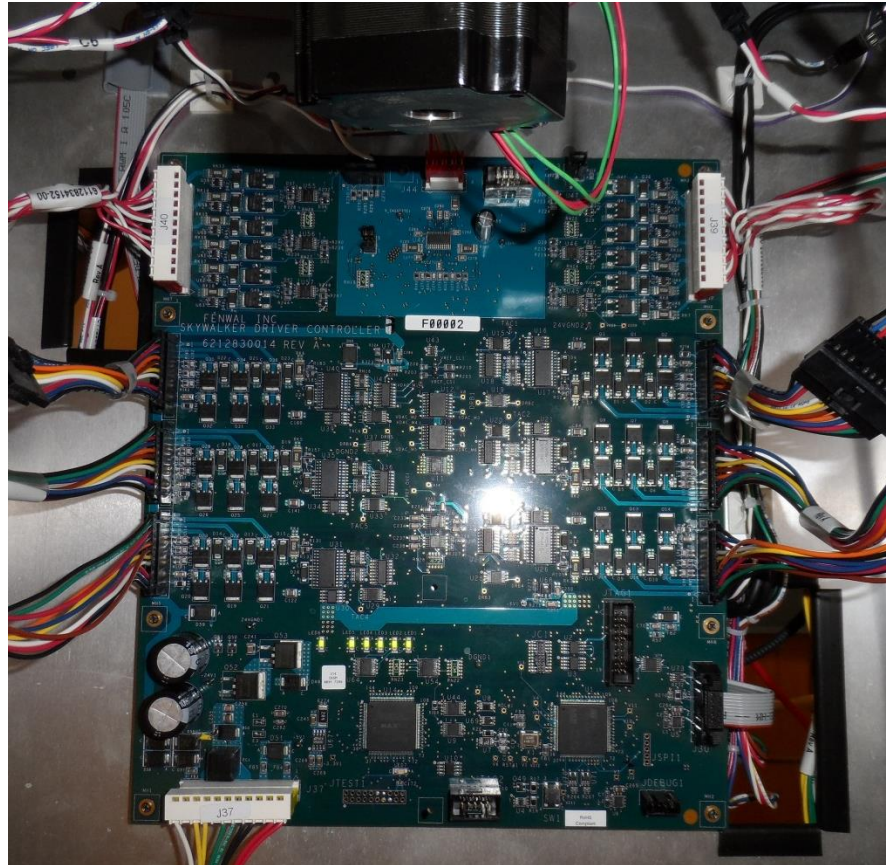

Part	Part #
Driver PCB	6412830018S



Top Panel

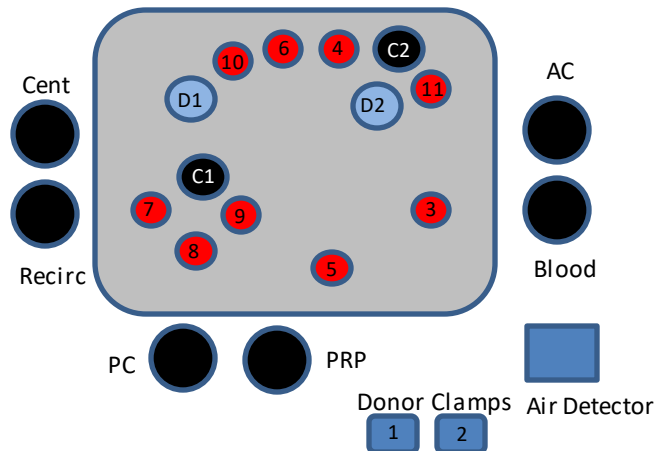

Part	Part #
Valve (Solenoid)	6112835126S
Donor Draw Clamp (purple)	6112836088S
Donor Return Clamp (blue)	6112836086S
Centrifuge Pressure Sensor (Pos Only Pressure)	6212833764S
Donor Pressure Sensor insert (Vent cap blue shell qty=3)	033214658S
Donor Pressure Sensor O ring	6212835710S
Centrifuge, PC, PRP, Blood Pump (26:1 gear ratio)	6212837093S
Recirc, AC Pump (46.1:1 gear ratio)	6212837094S
Top Plate (toolless top plate)	6212837975S
Air Detector	6212680066S
Pump Head	6112833254S



Display





Part	Part #
Display Assembly (Complete assembly with cables)	6112838067S




Main Controller Board:

Part	Part #
Main Controller Board	6412830004S

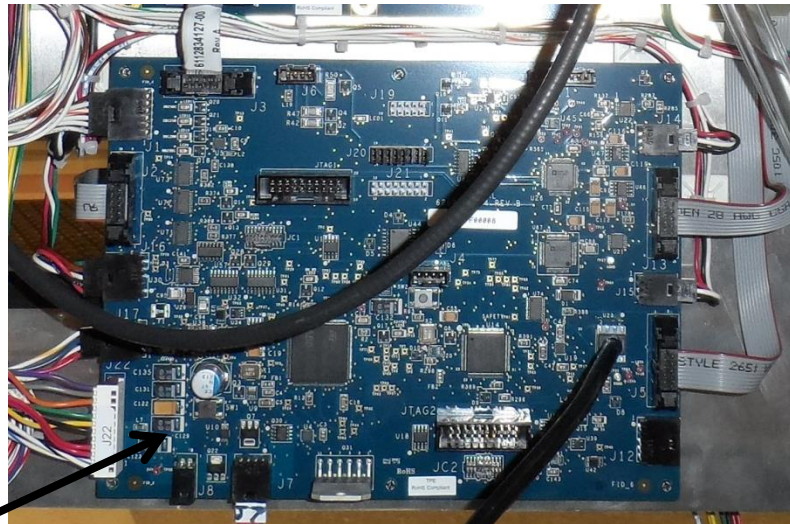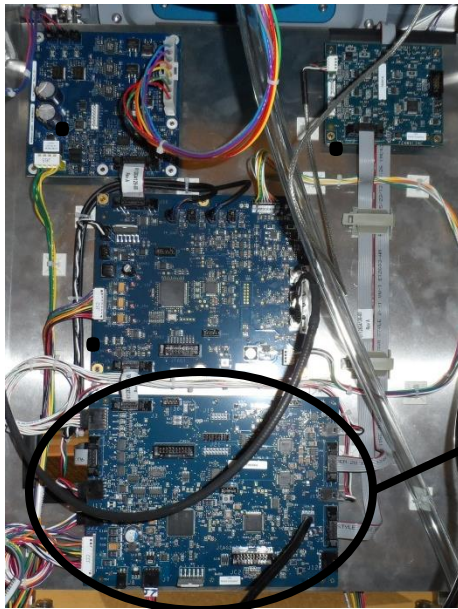


Interface Controller Board:

Part	Part #
Interface Controller Board (low gain)	64128300365

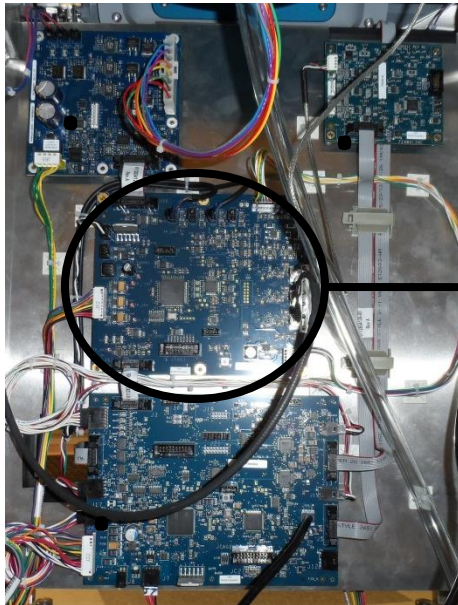



Leak Detection & Interface Detector:

Part	Part #
Leak/Temp PCB	6412830010S
Interface Detector Assembly (includes fiber optic cables)	6212833856S
Fan assembly	6112836158S
Speaker	6212836495S

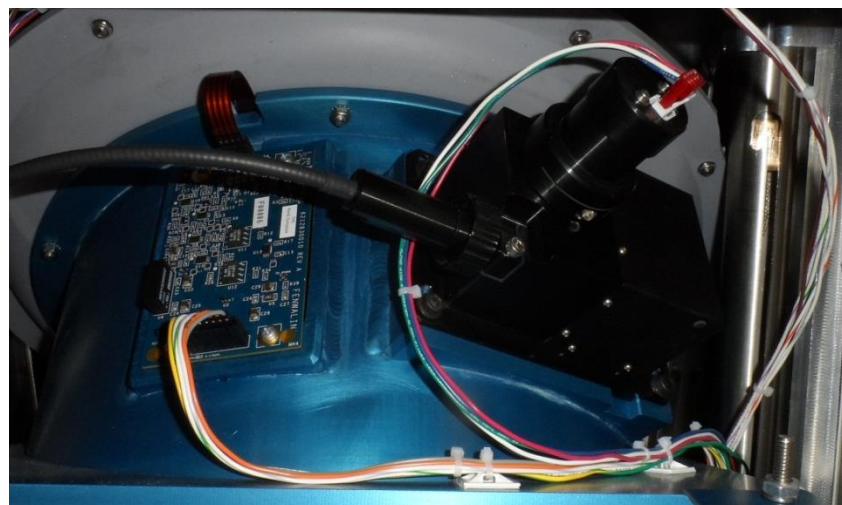



Centrifuge Controller Board:

Part	Part #
Centrifuge Controller Board	6412830007S

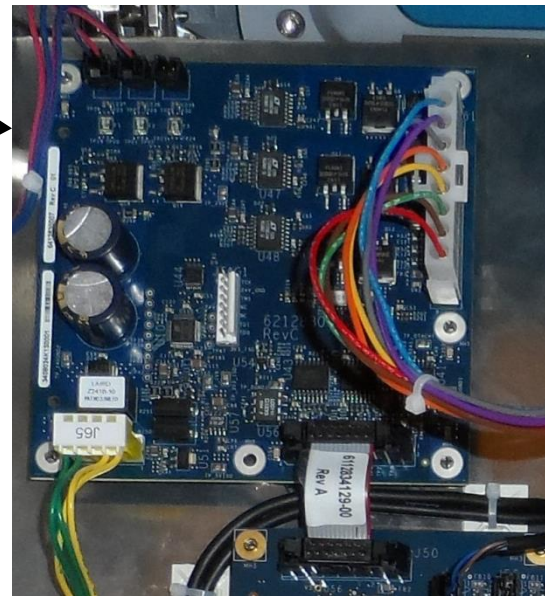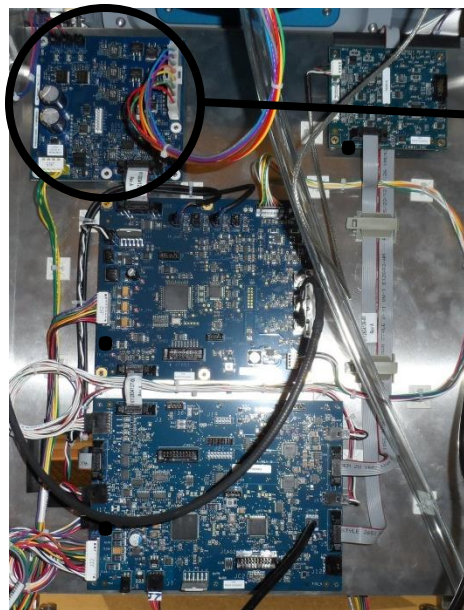

Weight Scale Board:

Part	Part #
Weight Scale Board	6412830011S

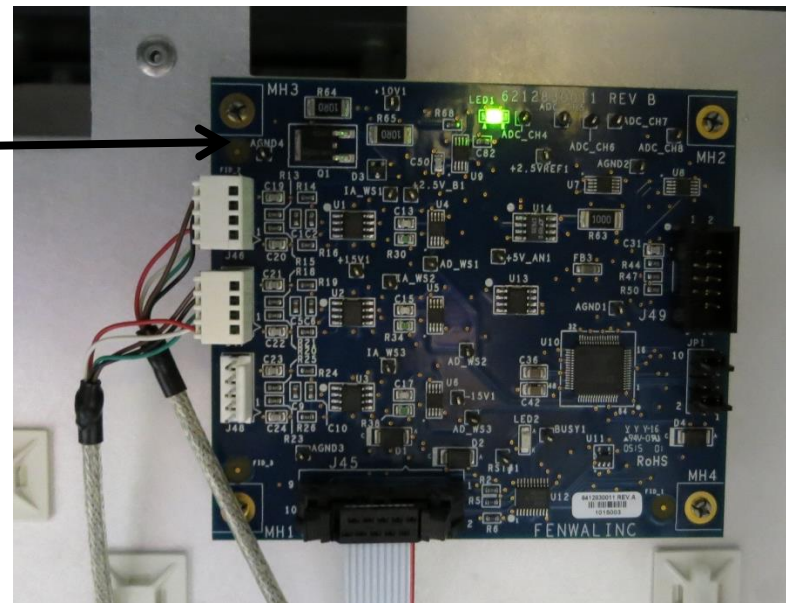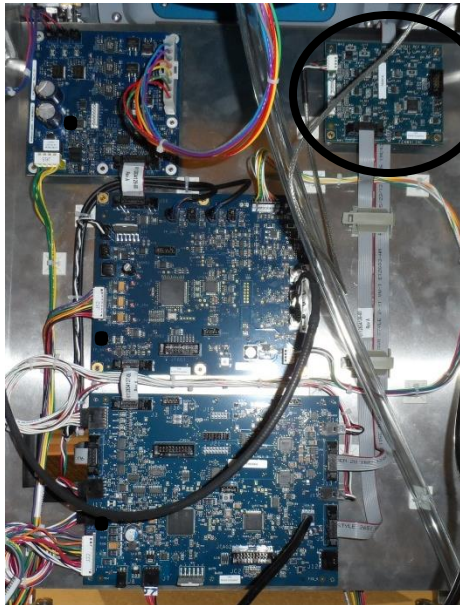


Weight Scale Assemblies




Part	Part #
Weight Scale (load cell 10lb capacity)	6212834023S
Plasma Hook (yellow)	6212836198S
ACD Hook (red)	6212683334S
White Hook	6212836199S
Gray Hook	6212683333S
Kit with 2 scales and pole	6212835842S



Power Supply Distribution

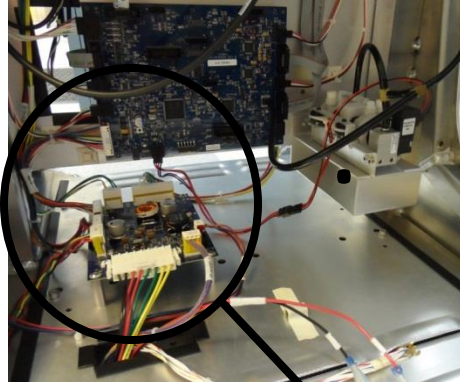


Part	Part #
Power Supply Distribution PCB (Battery Charger)	6412830012S

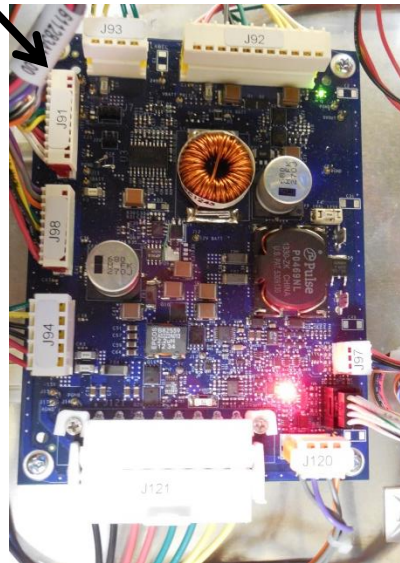


Power Supply



Power Supply

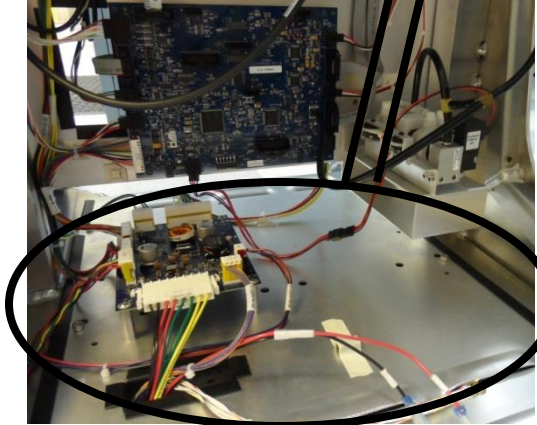


Remove panel w power supply distribution board

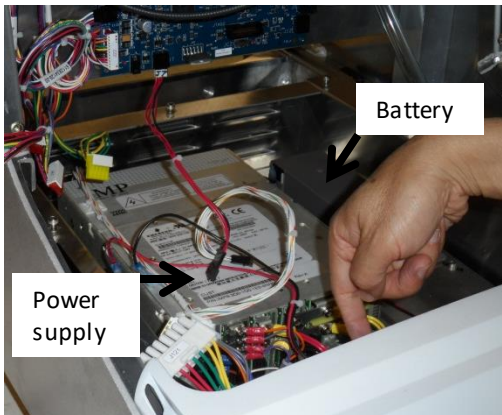


Part	Part #
Power Supply (Modular 1200w)	6212834204S
Battery (12v)	851743001S





Compressor

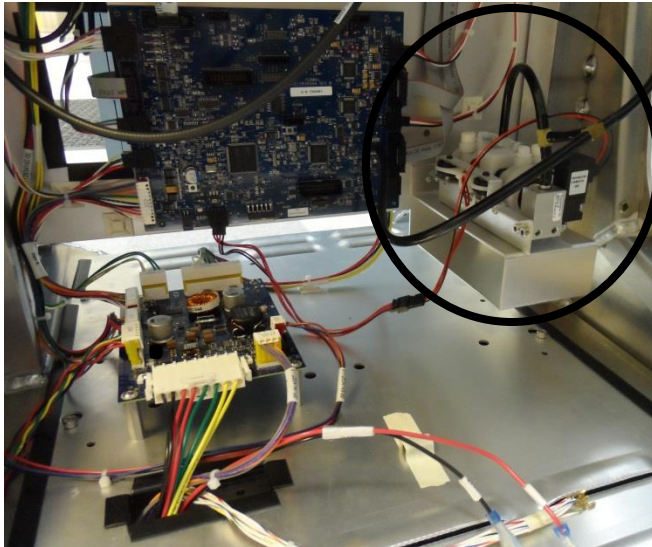

Part	Part #
Compressor (Cuff Compressor Assembly)	6112833313S

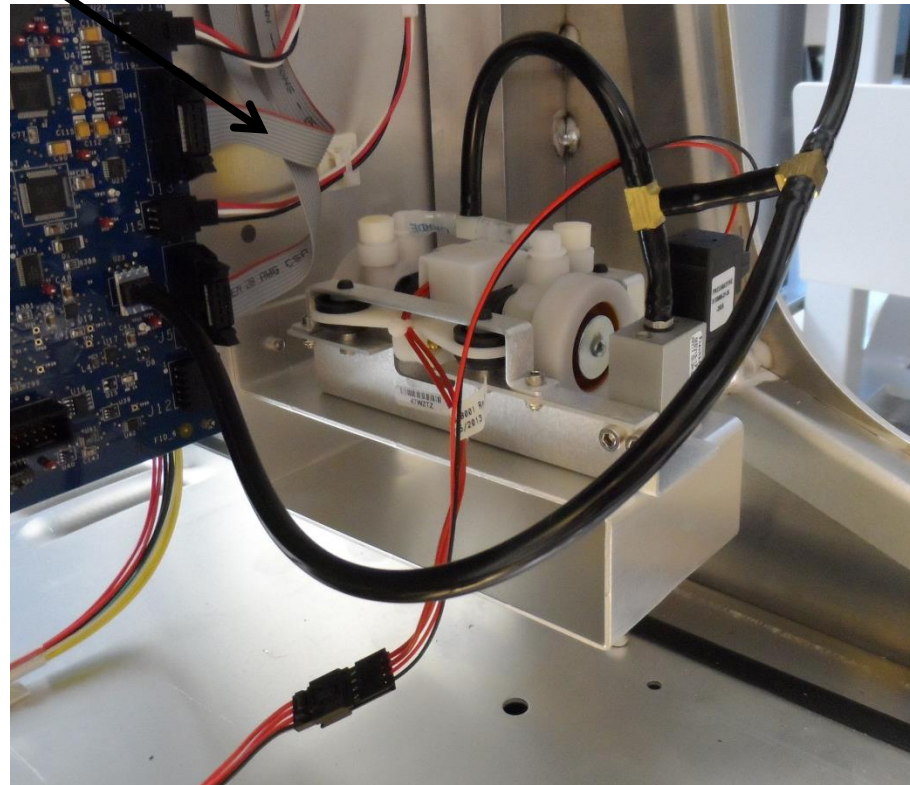

Centrifuge

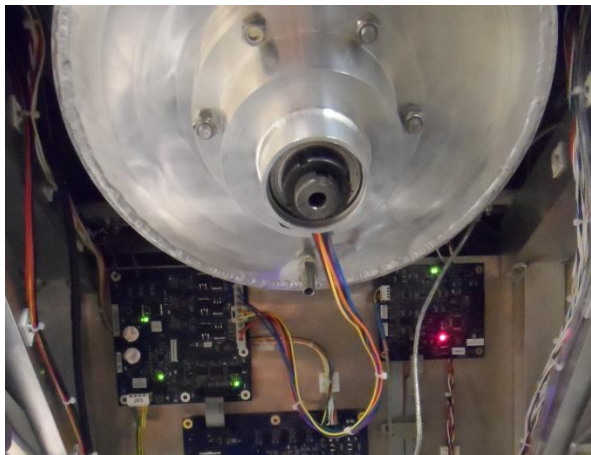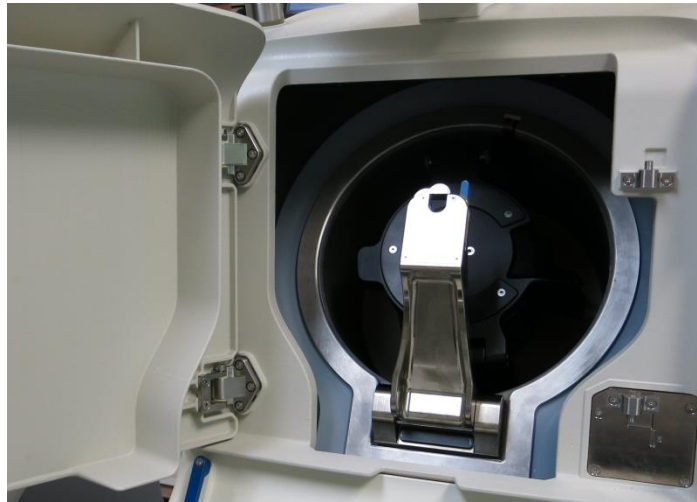

Door Lock



Isolators



Part	Part #
0 Omega Arm (Plate and Handle)	6112838050S
Bowl aka Spool Holder (Bowl 2 Omega)	A10416001S
Spool	6112685306S
Leak Detector (Flex circuit leak sensor)	6412830009S
Motor (Centrifuge Motor Assembly)	6112833801S
Rotor Assembly	6112833736S
Reflector (Shared with Amicus)	6212688547S



Door Lock




Part	Part #
Door Lock Assembly	6112835441S
Door Handle Assembly	6112835435S

Isolators



Part	Part #
Vibration Isolator	6212833858S

Pressure Sensor Board

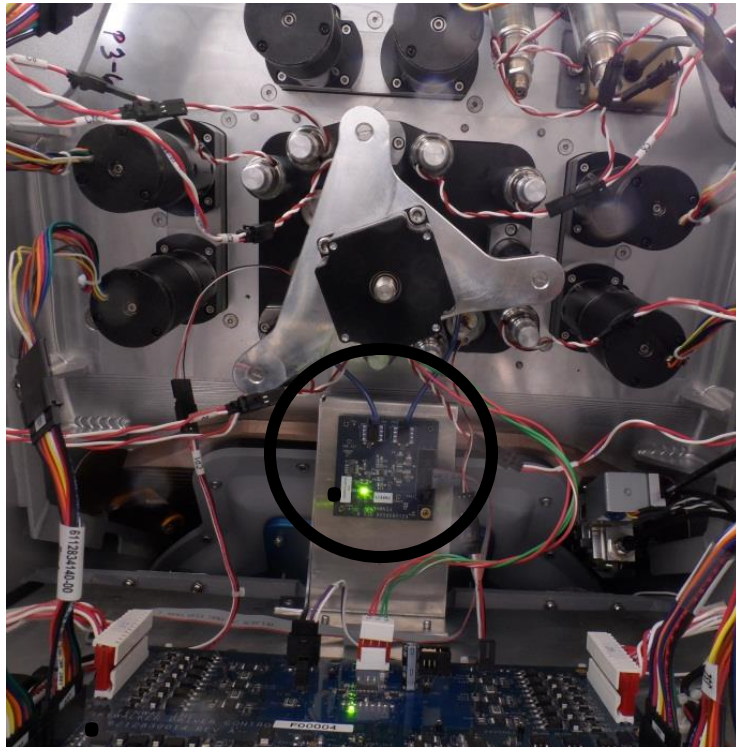


Part	Part #
Pressure Sensor Board	64128300085



Display & Communications Board

Part	Part #
Display SBC Assembly	6412830026S

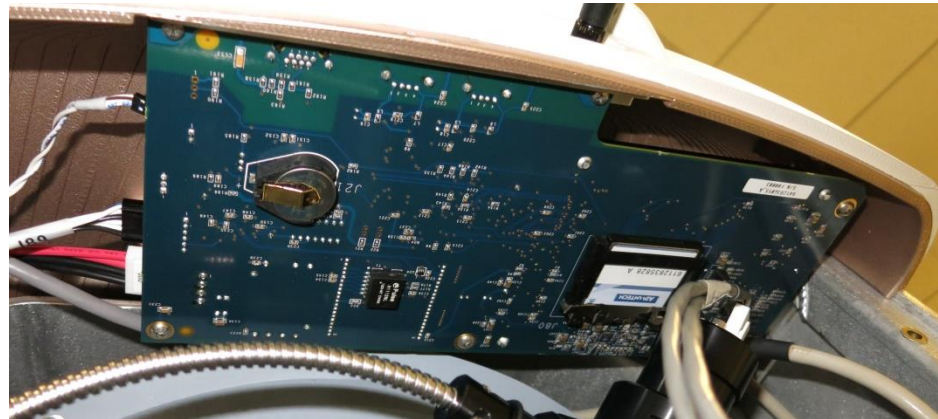

AmiCORE Specialty Tools



Clamp
Removal



Donor Pressure Sensor
Cover Replacement



Rotor Removal



Bowl Tower
Removal



Interface Alignment



O Omega Arm Alignment



Weight Scale
Installation



Donor P.S.
Tubing Removal



Donor P.S.
Housing Removal



Cent PS
Removal








Optical Fiber Continuity Tester

6212835817S

Interface Detector Alignment Fixture

6112835595S

Use: Alignment of Interface Detector Assembly

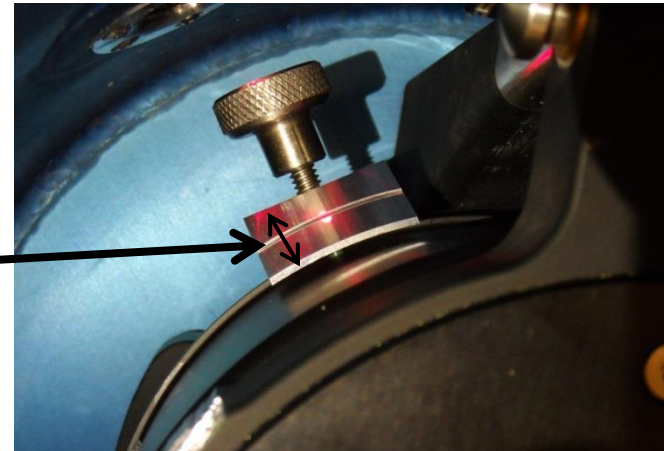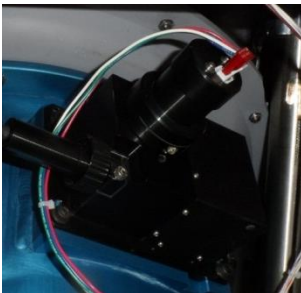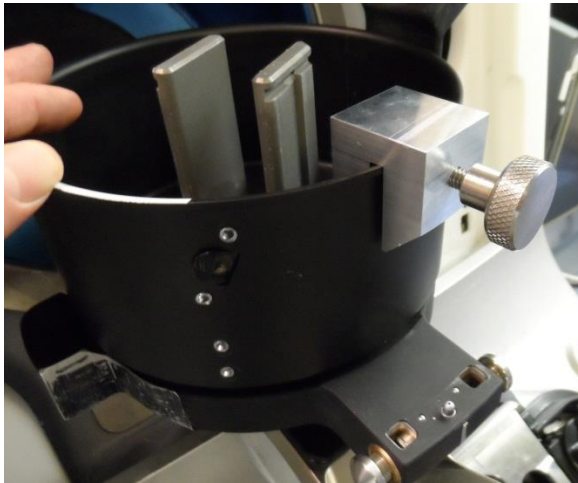



Fixture Pry Bar

6212836633S

Fixture Pressure Sensor Adaptor

6212836632S

Use: Replacement donor pressure sensor cover







Tool Spanner Large

6212834829S


Use: Removal of rotor assembly

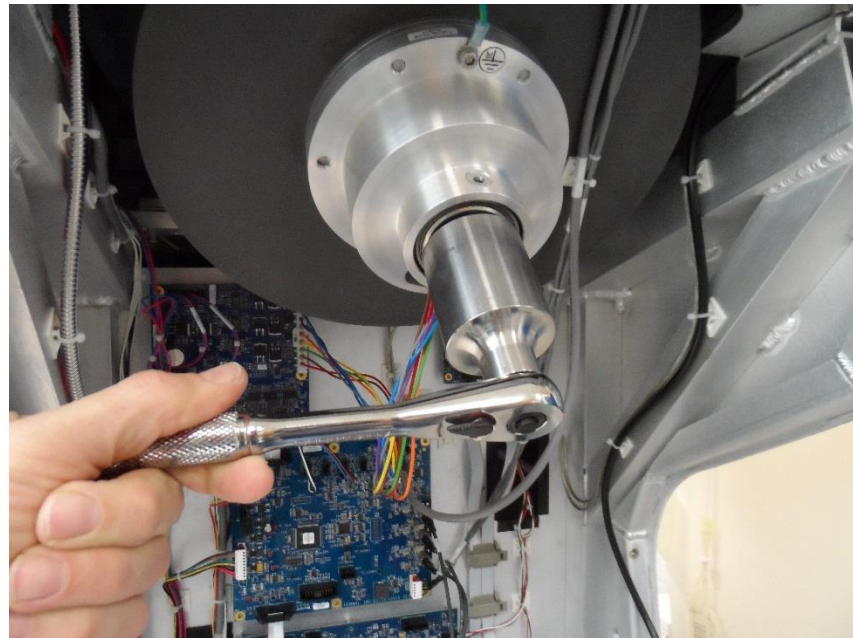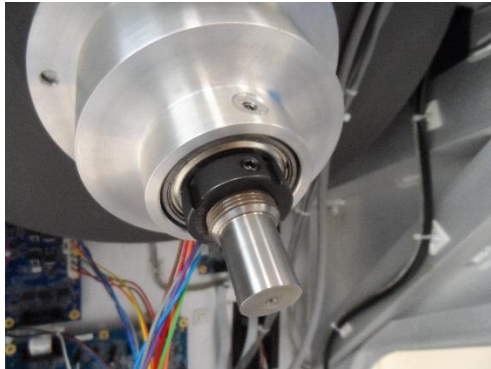



Tool Spanner Small

6212834828S

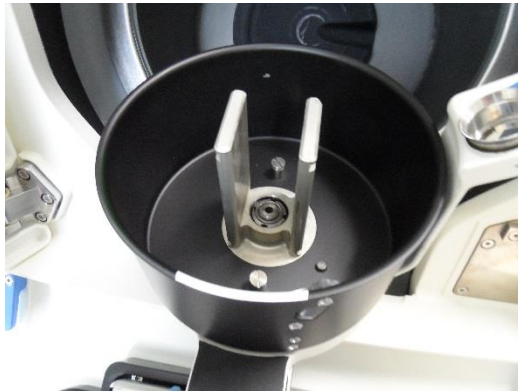



O Omega Arm Alignment Fixture


6112835748S

Use: Alignment of O Omega Arm



- 1) Remove rotor
- 2) Install alignment fixture
- 3) Close O Omega arm and ensure it aligns with top post of tool. If not, loosen four O Omega screws and align with post then retighten



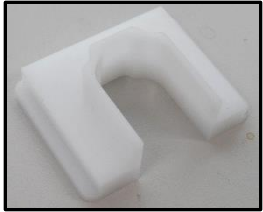

Weight Scale Alignment

6212836720S

Use: Alignment of weight scale assembly







1/4 In Socket Custom

6212836705S

Use: Removal of donor pressure sensor tubing and tubing nut

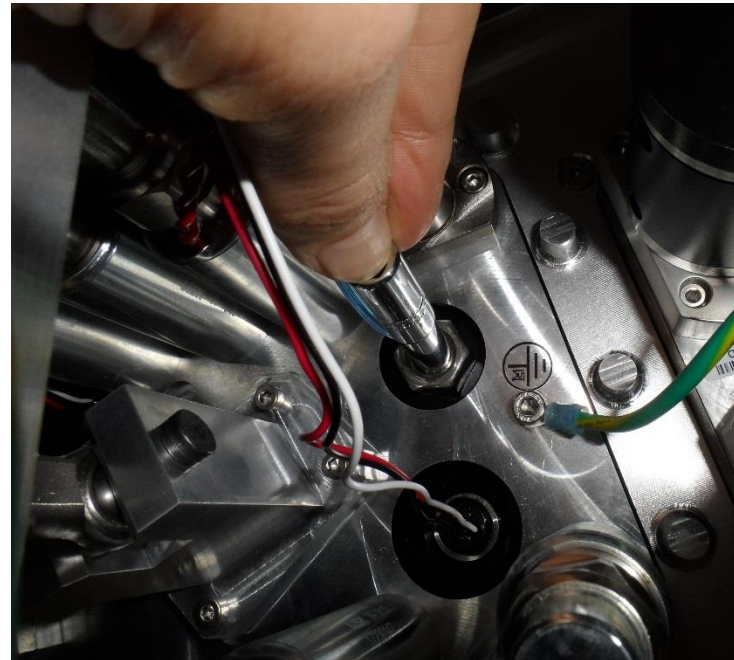
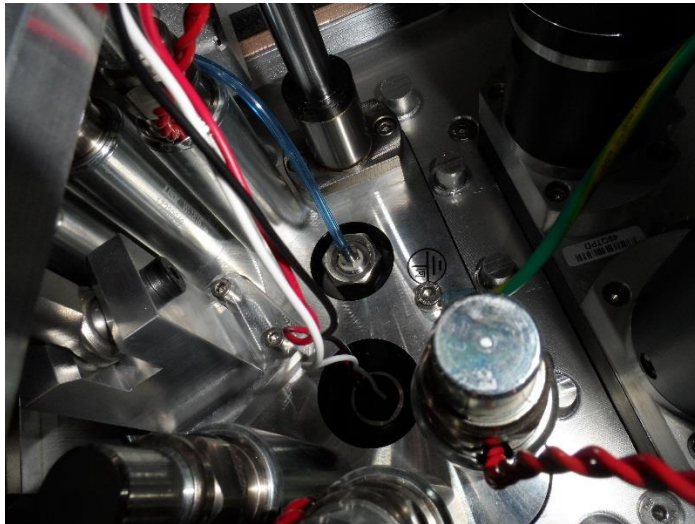




19 mm Socket Custom

6212836707S

Use: Removal of complete donor pressure sensor housing with tubing and tubing nut

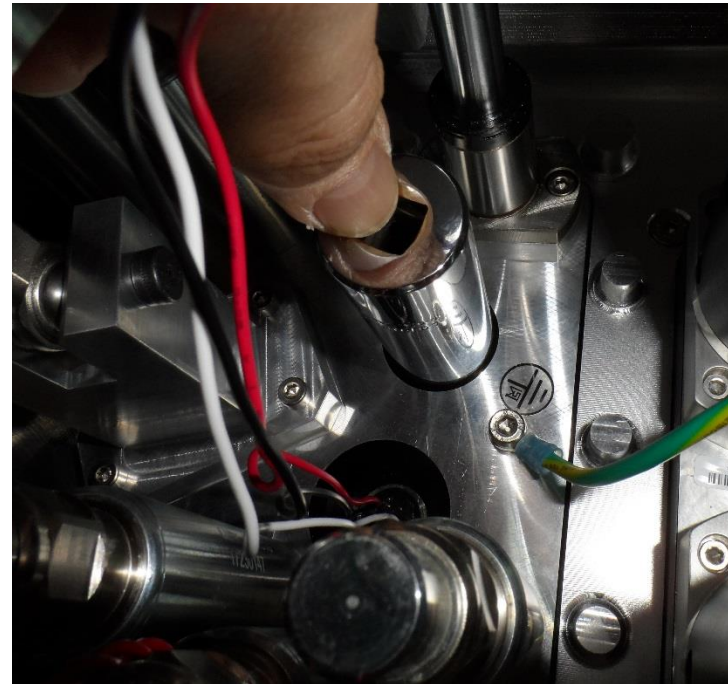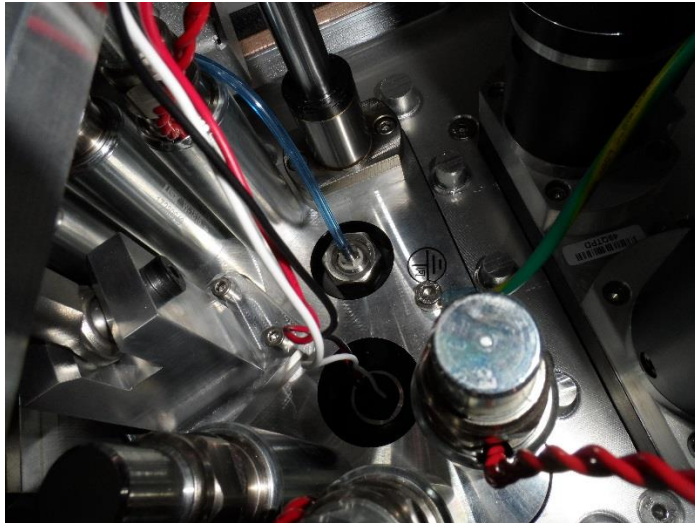




11/16 in Socket Custom

6212836706S


Use: Removal of centrifuge pressure assembly

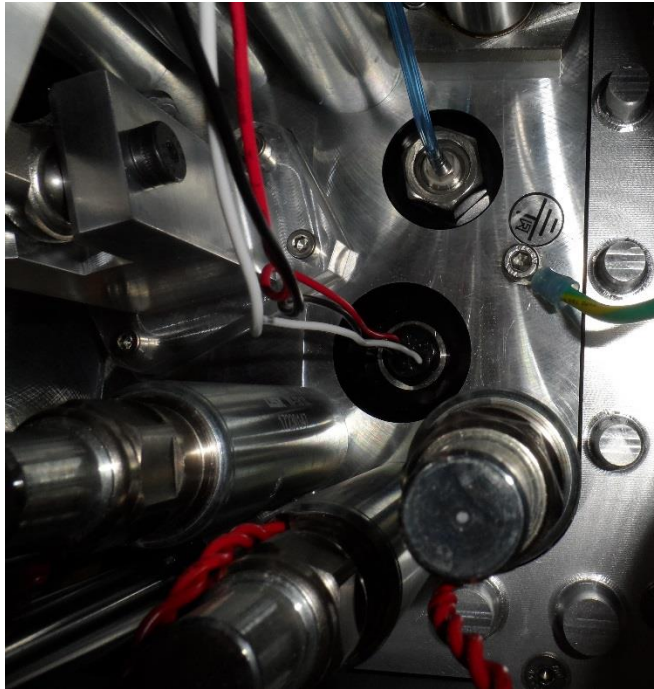



Clamp Removal Tool (shared Amicus) A40025001

Use: Removal of nut on donor clamps

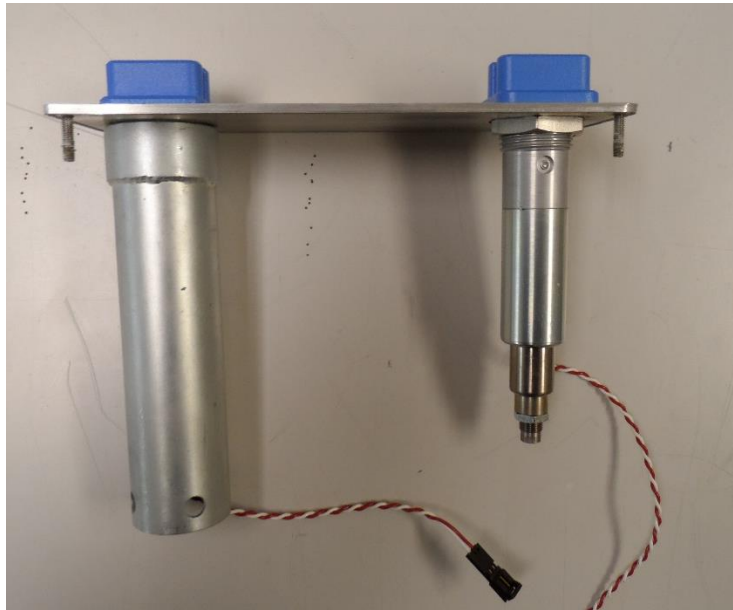

REVISION HISTORY

Revision	Date	Notes
A	9/3/20	Initial Release

