

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000775	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	19-May-2011
Title:	SM DG57-DG517. Unpacking and location remarks.			Page 2 of 5

DG-57 / DG517 Service Manual

UNPACKING AND LOCATION REMARKS

Procedure description
This document describes how to unpack the instrument and locate it to its final position. <i>Note: Document REGD-0000775_v3 replaces document REGD-0000775_v2.</i>

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000775	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	19-May-2011
Title:	SM DG57-DG517. Unpacking and location remarks.			Page 3 of 5

EQUIPMENT

1. Screwdriver (electrical is recommended)
2. Cutter



SET UP

1. ---

PROCEDURE 1: Location remarks


Process	Expected Results
1. Do not place the equipment so that it is difficult to access for maintenance or servicing. 10-20 cm space for service on the table beside the right and left sides of the instrument is recommended.	---
2. A wheel table with brake capability is recommended.	---
3. The inclination of the equipment can be adjusted by means of its feet.	---
4. Ambient Light: the instrument should not be located in an area of direct sunlight.	---

PROCEDURE 2: Unpacking the instrument

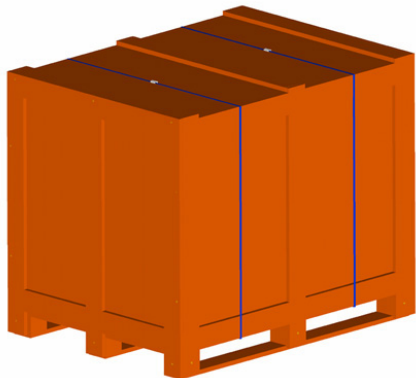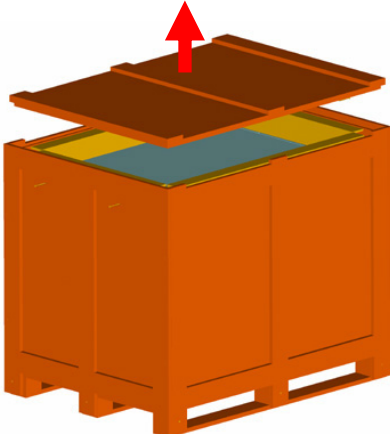
Process	Expected Results
5. Remove bands and top cover. To do that, 9 screws must be removed.   Picture 1 Picture 2  DANGER! Screws may be sharp. The cover must be lifted by two people.	---
6. Remove all items placed on top of the instrument.	---
7. Remove laterals and protective foams.	---

Process

Expected Results

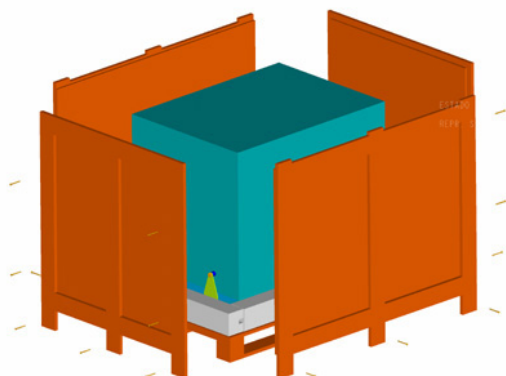

Picture 3

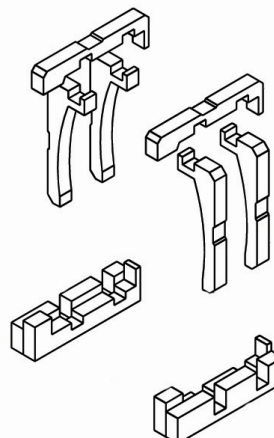

Picture 4

8. Open the insulating bag that covers the instrument.




DANGER!

The use of the cutter may be dangerous.

The instrument is ready to be placed in its final position

9. Lift the instrument to its final position using the handles at both sides of the instrument. 4 people are necessary to lift the instrument.

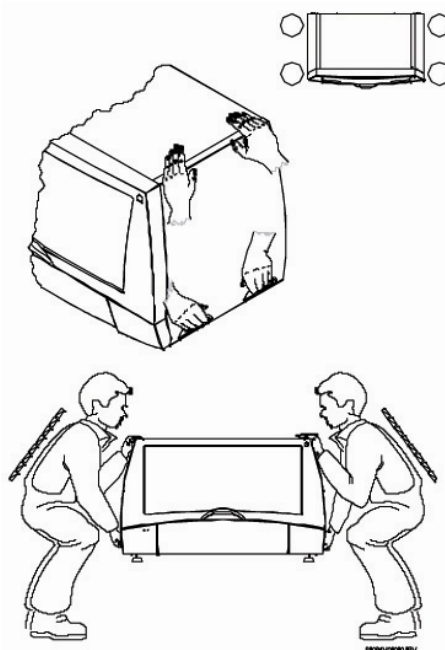

Picture 5

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000775	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	19-May-2011
Title:	SM DG57-DG517. Unpacking and location remarks.			Page 5 of 5

Process	Expected Results
 DANGER! The weight of the equipment is considerable, so it is advisable to strictly follow these instructions to place it into position. It must be lifted by four people, holding it in the manner shown (hands in the handles and with their backs erect and flexing the legs).	

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	Document REGD-0000775_v2 replaces document 57003V0M000.	Document reference changed.
3.0	See document header.	Added `previous version `data`, `change` and `rationale` inside document change history section.	Update change history for previous document history modification.
3.0	See document header.	Document REGD-0000775_v3 replaces document REGD-0000775_v2.	Document version changed.
3.0	See document header.	Replaced commercial name by project name.	Delete commercial name.
3.0	See document header.	Header title has changed.	Previous header title has been substituted by project and document title.