

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000798	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	16-Jun-2011
Title:	SM DG57-DG517. PCA Samples-reagents-2 and PCA Centrifuge-2 fuse status ver			Page 2 of 4

DG-57 / DG-517 Service Manual

PCA SAMPLES-REAGENTS-2 AND PCA CENTRIFUGE-2 FUSE STATUS VERIFICATION

Procedure description
This procedure describes how check the fuses in Samples-Reagents and Centrifuge PCA's. <i>Note: Document REGD-0000798_v3 replaces document REGD-0000798_v2.</i>


GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000798	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	16-Jun-2011
Title:	SM DG57-DG517. PCA Samples-reagents-2 and PCA Centrifuge-2 fuse status ver			Page 3 of 4

EQUIPMENT

1. Multimeter.

SET UP


1. Turn Instrument on.
2. Open the rear panel.
3. Turn Instrument off.

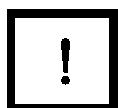

CAUTION!

Short-circuit between any powered part and ground, may damage the instrument.

PROCEDURE

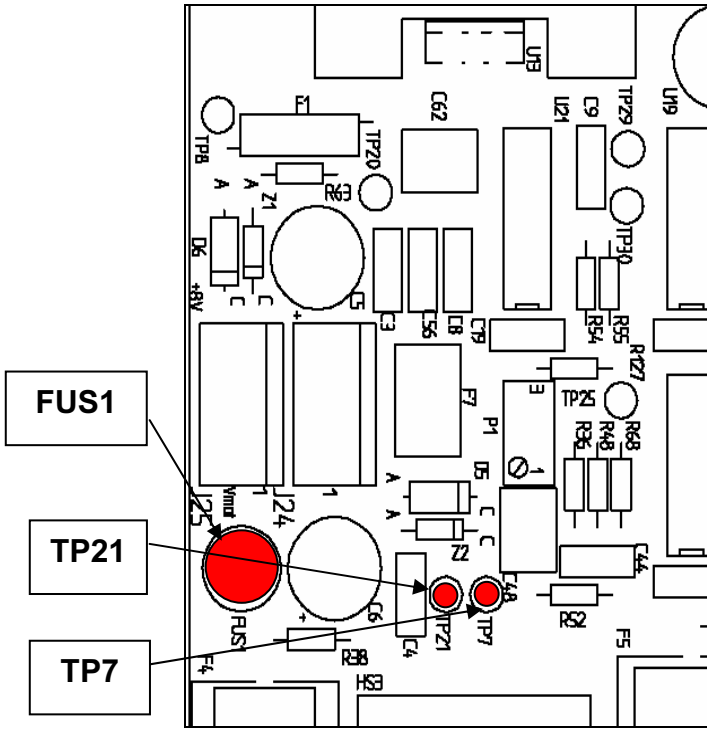Process	Expected Results
Checking the 24V fuse 1. To check the fuse status, measure the voltage between TP21 (negative) and TP7 (positive) with the multimeter (see picture 1)	23.95 ± 0,05 V
 Picture 1 2. If the measured value is not the expected, then the FUS1 may be broken.	

GRIFOLS Diagnostic Grifols, S.A.	Reference:	REGD-0000798	Version:	3.0, CURRENT
	Status:	Effective	Effective Date:	16-Jun-2011
Title:	SM DG57-DG517. PCA Samples-reagents-2 and PCA Centrifuge-2 fuse status ver			Page 4 of 4

DOCUMENT CHANGE HISTORY

Document change history			
Version	Date	Change	Rationale
2.0	2005/05/09	<i>Document REGD-0000798_v2 replaces document 57026V0M000.</i>	Document reference changed.
3.0	See document header.	Added previous version 'data', 'change' and 'rationale' inside document change history section.	Update change history for previous document history modification.
3.0	See document header.	<i>Document REGD-0000798_v3 replaces document REGD-0000798_v2.</i>	Document version changed.
3.0	See document header.	Document title.	Commercial reference substituted by project reference.
3.0	See document header.	Document header.	Increase information to easy document search.